DATA EXPORT INTERFACE PROGRAMMING GUIDE

IntelliVue Patient Monitor X2, MP Series, MX Series

Patient Monitoring

4535 643 22831

PHILIPS

1	About this Guide	9
	Who Should Use this Guide?	9
	About the Data Export Interface	10
	Data Export Interface Features	10
	Note About Changes in Rev. G.0	10
	Manufacturer's Information	11
	Trademark Information	11
2	Connecting to the Network	13
	Connecting to the Network via a LAN Interface	13
	Connection via Hub/Switch	14
	Connection with Cross-over Cable	14
	Avoiding Current Leakage	15
	Using the Monitor with an Installed, Wired Network	15
	Configuring the LAN Interface	16
	Configuring the Network Address	16
	Configuring the LAN Data Export Setting	16
	Configuring the Network Setting	17
	Connecting to the IntelliVue MP20-90 or MX Series MIB/RS232 Interface	18
	LEDs on the MIB/RS232 Board (MP20-90 only)	19
	Connecting to the IntelliVue MP5 Monitor MIB/RS232 Interface	21
	Configuring the IntelliVue Monitor MIB/RS232 Interface	22
3	Protocol Concept	23
	Supported Transport Protocols	23
	UDP/IP Protocol	23
	Fixed Baudrate Protocol	24
	Auto Speed Protocol	24
	Protocol Model	24
	Protocol Dialog	25
	Connect Indication	25
	Association Request	26
	Association Result	26
	MDS Create Event Report	26
	MDS Create Event Result	26
	Poll Data Request	26
	Poll Result	26
	Association Release Request	27
	Association Release Result	27
	Association Abort	27
	More Information	27
	Connection Time-out Mechanism	27
	Network Load Consideration	28
	Input Data	28

Output Data	28
4 Definition of the Transport Protocols	29
Transport Protocols for the LAN Interface	29
UDP/IP	29
IP Address	29
UDP Port Number	29
Transport Protocols for the MIB/RS232 Interface	30
The Fixed Baudrate Protocol, RS232 Port Settings	30
Framing	30
Header Information	30
Frame Check Sequence Field	31
Transparency	31
Frame Abort	31
The AutoSpeed Protocol	32
Establishing a Connection	32
5 Definition of the Data Export Protocol	35
Definitions Shared by Protocols	35
Byte Order	35
Byte Alignment	35
Bit Order	35
Common Data Types	36
Basic Data Types	36
Absolute Time	36
Relative Time	36
OID Type	37
Private OID	37
TYPE	37
Handle	38
Global Handle	38
Managed Object Identifier	38
Attribute Value Assertion	38
Attribute List	39
String	39
Variable Label	40
TextId	40
FLOAT-Type	40
Protocol Command Structure	41
Session/Presentation Header	42
Remote Operation Header	42
Remote Operation Invoke	43
Remote Operation Result	43
Remote Operation Linked Result	44
Remote Operation Error	45
Command Header	47

Event Report	47
Event Report Result	48
Action	49
Action Result	49
Get	50
Get Result	50
Set	50
Set Result	51
Command Structure Summary	51
Protocol Commands	52
Notation	52
Device Discovery Messages	53
CONNECT INDICATION EVENT	53
Connection Startup	53
MDS CREATE EVENT	54
MDS CREATE EVENT RESULT	55
Specific Data Access Commands	55
SINGLE POLL DATA REQUEST	55
SINGLE POLL DATA RESULT	56
EXTENDED POLL DATA REQUEST	59
Attribute: Time Periodic Data Poll	60
Limiting the Number of Objects in the Poll Result	61
Attribute: Number of Prioritized Objects	61
EXTENDED POLL DATA RESULT	61
Keep Alive Message	63
Specify Objects in the Poll Result	63
GET PRIORITY LIST REQUEST	63
GET PRIORITY LIST RESULT	63
SET PRIORITY LIST REQUEST	64
SET PRIORITY LIST RESULT	64
6 Definition of the Association Control Protocol	65
<hr/>	
Protocol Command Structure	65
Protocol Commands	65
Session Headers	67
Message Encoding	67
Association Request Message	67
Attribute: Poll Profile Support	69
Attribute: Poll Profile Extensions	71
Release Request Message	72
Abort Message	72
Message Parsing	72
Association Response Message	73
Refuse	73
Release Response	73

7 Attribute Data Types and Constants Used

Numeric Objects	75
Numeric Object Attributes	75
Attribute: Handle	75
Attribute: Type	76
Attribute: Numeric Observed Value	76
Attribute: Compound Numeric Observed Value	77
Attribute: Absolute Time Stamp	78
Attribute: Relative Time Stamp	78
Attribute: Label	78
Attribute: Label String	78
Attribute: Display Resolution	78
Attribute: Color	79
Attribute: Metric Specification	79
Attribute MetricModality	81
Attribute Groups	81
Dynamic Context Changes	81
Wave Objects	82
Wave Object Attributes	82
Attribute: Handle	82
Attribute: Type	82
Attribute: Metric Specification	82
Attribute: Sample Array Specification	82
Attribute: Sample Array Fixed Value Specification	83
Attribute: Sample Period	84
Attribute: Label	84
Attribute: Label String	84
Attribute: Metric State	84
Attribute: Unit Code	85
Attribute: Color	85
Attribute: Measure Mode	85
Attribute: Metric Info Label	85
Attribute: Metric Info Label String	85
Attribute: Scale and Range Specification	86
Attribute: Sample Array Physiological Range	86
Attribute: Visual Grid	86
Attribute: Sample Array Calibration Specification	87
Attribute: Sample Array Observed Value	87
Attribute: Compound Sample Array Observed Value	87
Attributes Groups	88
Enumeration Objects	88
Enumeration Object Attributes	88
Attribute: Handle	88
Attribute: Type	89
Attribute: Metric Specification	89

Attribute: Label	89
Attribute: Label String	89
Attribute: Color	89
Attribute: Enum-Observed-Value	90
Attribute Groups	90
System Objects	91
System Objects Attributes	91
Attribute: Handle	91
Attribute: System Type	91
Attribute: System Model	91
Attribute: System ID	91
Attribute: Nomenclature Version	92
Attribute: System Localization	92
Attribute: System Specification	93
Attribute: Mds General System Info	94
Attribute: Production Specification	95
Attribute: MDS Status	96
Attribute: Bed Label	96
Attribute: Operating Mode	96
Attribute: Application Area	97
Attribute: Date and Time	97
Attribute: Relative Time	97
Attribute: Altitude	98
Attribute: Line Frequency	98
Attribute: Association Invoke ID	98
Attribute Groups	98
Alert Monitor Object	99
Attributes of the Alert Monitor Object	99
Attribute: Handle	99
Attribute: Type	99
Attribute: Device Alert Condition	99
Attribute: Device T-Alarm List	101
Attribute: Device P-Alarm List	102
Attribute Groups	103
Patient Demographics Object	103
Attributes of the Patient Demographic Object	103
Attribute: Handle	103
Attribute: Pat Demo State	103
Attribute: Patient Type	104
Attribute: Patient Paced Mode	104
Attribute: Given Name	104
Attribute: Middle Name	104
Attribute: Family Name	104
Attribute: Patient ID	105
Attribute: Encounter ID	105
Attribute: Patient Sex	105

Attribute: Date of Birth	105
Attribute: Patient Height	106
Attribute: Patient Weight	106
Attribute: Patient Age	106
Attribute: Patient BSA	106
Attribute: Patient BSA Formula	106
Attribute: Notes1	107
Attribute: Notes2	107
Attribute Groups	107
Patient Conflict Handling	107
Connect Indication Attributes	107
Attribute: System Type	108
Attribute: Protocol Support	108
Attribute: System Localization	108
Attribute: IP Address Information	109
Partition IDs	110
Object Classes	111
Physiological Identifier	115
Numerics	115
Enumerations	180
Waves	181
Attribute IDs	191
Component IDs	194
Unit Codes	195
Alert Codes	204
ECG/HR/Arrhy	205
ST	206
QT Analysis	206
Resp	207
Derived Measurements	207
C.O./CCO	208
EEG	209
BIS	209
NMT	210
Temp	211
Invasive Pressure	212
SpO ₂	212
SvO ₂	214
CO ₂	214
AGM	216
System	220
AlarmMgr	221
NBP	222
TcGas	222
VueLink	224
Battery	225

Telemetry	225
Spirometry	227
Predictive Temp	228
Protocol Watch	228
Metabolics	230
Intellibridge	231
Short Range Radio	231
Private Unicode Characters	252
List of Constants Used Within the Protocol Definition	252
RO Types	252
ROLS Identifier	253
ROSE Commands	253
ROER Error Values	253
Action and Event Types	253
Protocol Identification	253
Association Control	254
Label Mapping Table	255
8 Building a Computer Client	279
Interfacing the LAN interface with UDP/IP	279
Setting Up the BootP Server	279
Parsing the Connect Indication Message	279
Interfacing the MIB/RS232 Interface with the Fixed Baudrate Protocol	280
Interfacing the MIB/RS232 Interface with the AutoSpeed Protocol	281
Establishing an Association	282
Accessing Data	283
Message Frequencies	283
Single and Extended Polling	283
Availability of Data	284
Parsing the Poll Result	284
Parsing AttributeLists	285
Interpreting Data from Numerics	286
Interpreting Data from the Alert Monitor	286
Interpreting Wave Data	286
9 Troubleshooting	289
Further Troubleshooting	290
10 Protocol Examples	291
Data Export Protocol Examples	291
CONNECT INDICATION EVENT	291
MDS CREATE EVENT	291
MDS CREATE EVENT RESULT	292
SINGLE POLL DATA REQUEST	292
SINGLE POLL DATA RESULT	292
SINGLE POLL DATA RESULT (LINKED)	293

EXTENDED POLL DATA REQUEST	294
EXTENDED POLL DATA RESULT	294
GET PRIORITY LIST REQUEST	294
GET PRIORITY LIST RESULT	295
SET PRIORITY LIST REQUEST	295
SET PRIORITY LIST RESULT	296
AttributeList	296
Association Control Protocol Examples	298
ASSOCIATION REQUEST	298
ASSOCIATION RESPONSE	299
REFUSE	300
RELEASE REQUEST	301
RELEASE RESPONSE	302
ASSOCIATION ABORT	303
User Data	304

About this Guide

This Programming Guide is for use with the Philips IntelliVue X2, MP Series and MX Series patient monitors, hereafter referred to as the IntelliVue monitor. It describes the functionality in the monitor software version \leq J.xx.xx .

The information in this Programming Guide describes the capability of the Data Export Interface. It is the responsibility of the user to create applications using the capability provided.

This device is not intended for home use.

In this guide

- A warning alerts you to a potential serious outcome, adverse event or safety hazard. Failure to observe a warning may result in death or serious injury to the user or patient.
- A caution alerts you where special care is necessary for the safe and effective use of the product. Failure to observe a caution may result in minor or moderate personal injury or damage to the product or other property, and possibly in a remote risk of more serious injury.

Who Should Use this Guide?

This programming guide is intended to be used by software professionals and biomedical engineers at medical research clinics or industrial institutions.

To successfully create an application, users should have a **good** working knowledge of:

- Advanced software application design.
- C and/or C++ Programming Language.
- General digital communications theory.
- Local Area Network configuration guidelines and communication protocols.
- RS232 communication protocols and the IrDA protocol.

Given this background knowledge, this Programming Guide provides the information necessary to create your own applications.

Philips cannot provide any technical assistance for individual programming efforts.

About the Data Export Interface

This document describes the IntelliVue Data Export Interface. Using a communication interface protocol, data from the Philips IntelliVue Patient Monitor can be transferred via the Local Area Network (LAN) Interface or Medical Information Bus (MIB/RS232) Interface to an external Computer.

By creating basic applications using the IntelliVue Data Export Interface, the following data can be accessed from the IntelliVue monitor:

- All measurement numerics (**excluding** second level numerics such as those obtained in the Calculations windows (Hemodynamic Calculations, Oxy Calculations, Ventilation Calculations) or those exclusively shown in the VueLink/IntelliBridge device data windows).
- Alarm data (real-time update rates up to 1024 ms).
- Wave data (see “Interpreting Wave Data” on page 286 for details)
- IntelliVue monitor system data.
- Patient demographic data entered by the user in the IntelliVue monitor.

The IntelliVue Data Export Interface cannot be accessed via the Local Area Network when the IntelliVue monitor is connected to the Philips LAN, e.g. to an Information Center (central station). Communication via the MIB/RS232 Interface is always possible (except with MP2/X2).

-
- CAUTION**
- Although alarm data can be accessed using the protocol, it must not be used as a real-time alarming system due to the delays in message transfer and the possibility of data loss.
 - The computer client (the interfacing system) and/or the user of the communication system must comply with applicable data privacy regulations.
-

Data Export Interface Features

- The IntelliVue Data Export Interface uses the Local Area Network (LAN) and MIB/RS232 interfaces.
- The LAN interface uses the standard UDP/IP transport protocol.
- The MIB/RS232 interface can be configured to use either a fixed or a variable baudrate protocol.
- The Data Export Protocol is a connection-oriented, message-based request/ response protocol on top of the transport protocol. The UDP and fixed baudrate transport protocols are connection-less, whereas the variable baudrate protocol is connection-oriented.
- The LAN interface supports automatic configuration of the network IP address with the standard BootP protocol.

Note About Changes in Rev. G.0

IntelliVue Rev. G.0 and higher differs from Rev F.0 and lower in the nomenclature of some numeric and wave labels. The labels that previously resided in the namespace NOM_EMFC are now merged into the NOM_SCADA namespace and the new defined NOM_SETTING namespace.

For details on identifying the software revision of the client interface protocol, please refer to “Building a Computer Client” on page 279

Manufacturer's Information

The information contained in this document is subject to change without notice. Philips makes no warranty of any kind with regard to this material, including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose.

Philips shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance or use of this material.

Philips assumes no responsibility for the use or reliability of its software on equipment that is not furnished by Philips. Philips reserves the right to change the protocol described in this document without prior notice. Philips makes no warranty of any kind with regard to software applications that are created by the user. Philips assumes no responsibility for adverse interaction of an external Computer Interface communicant with the Philips IntelliVue Series Patient Monitor.

This document contains proprietary information that is protected by copyright. All rights are reserved. No part of this document may be photocopied, reproduced or translated to another language without prior written consent of Philips.

Philips Medizin Systeme Boeblingen Gmbh
Hewlett-Packard-Str. 2
71034 Böblingen
Germany

© Copyright 2011. Koninklijke Philips Electronics N.V. All Rights Reserved.

Trademark Information

Microsoft ® is a U.S. registered trademark of Microsoft Corp. Windows 2000® is a U.S. registered trademark of Microsoft Corp.

INFRARED DATA ASSOCIATION (IrDA) is a trademark of the Infrared Data Association in the USA and other countries.

Connecting to the Network

The Philips IntelliVue Series monitor uses a standard IEEE802.3 10BaseT (10MBit/s) Local Area Network interface for the Data Export Capability.

The Data Export Interface via LAN is not available when the IntelliVue monitor is connected to the Philips LAN (e.g. to the Philips Information Center central station). Only devices approved for use with the Philips network may be connected to the Philips LAN.

Connecting to the Network via a LAN Interface

The IntelliVue monitor connects to the network using a standard unshielded LAN cable with an RJ45 connector. The network cable must be plugged to the orange-framed LAN connector of the IntelliVue monitor. Note that for IntelliVue MP2/X2 the LAN connector is located on the external power supply.

WARNING In order to maintain the galvanic isolation of the IntelliVue monitor, it is essential that UTP (Unshielded Twisted Pair) LAN cables must be used to connect the IntelliVue monitor to other devices.

The following LAN cables supplied by Philips can be used to connect the IntelliVue monitor:

- M3199AI #J10 - 3ft (0.91m), Part No. M3199-60103 (12NC: 453563337391)
- M3199AI #J11 - 7ft (2.1m), Part No. M3199-60104 (12NC: 453563337401)
- M3199AI #J12 - 12ft (3.6m), Part No. M3199-60105 (12NC: 453563337411)

The maximum cable length between the IntelliVue monitor and the Computer Client should never exceed 330ft (100m) in total.

Connection via Hub/Switch


The IntelliVue monitor and the Computer Client are connected to a standard Ethernet switch or hub using UTP LAN cables.

NOTE In order to avoid high latency and data loss and to ensure data privacy, a dedicated network that is exclusively used for patient data collection by IntelliVue monitor devices and Computer Clients must be used.

Connection with Cross-over Cable

You can connect the IntelliVue monitor directly to the Computer Client, without a network hub or network switch, by using a UTP network crossover cable. In this case, the connection is a point-to-point connection only (one IntelliVue monitor connects to one Computer Client).

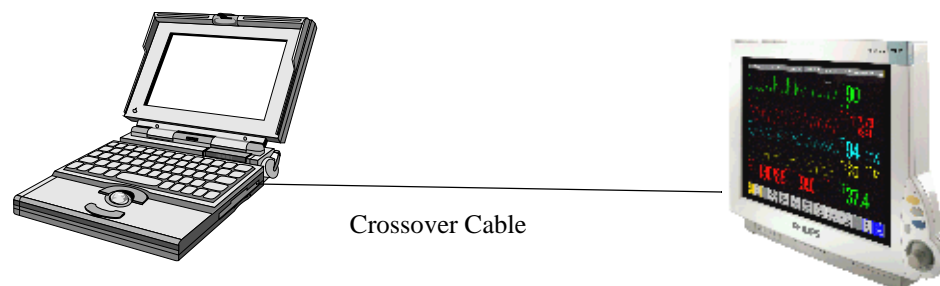The following cross-over LAN cables supplied by Philips can be used to connect an IntelliVue monitor:

- M3199-60101 (453563337371) 3Ft UTP Crossover cbl Orange, 0,9m
- M3199-60102 (453563337381) 12Ft UTP Crossover cbl Orange, 3,6m

Avoiding Current Leakage

You must use Unshielded Twisted Pair (UTP) LAN cables to connect the IntelliVue monitor to other devices.

The Computer Client and network infrastructure devices typically are not classified as medical devices and must be located outside the patient vicinity. The patient vicinity is defined as an area within 6ft (1.85m) of the perimeter of the patient's bed or within 7.5ft (2.3m) of the floor.


- If the Computer Client is installed in the patient vicinity and connected to the monitoring device, it must be correctly isolated from the mains power supply by an isolation transformer.
- If the Computer Client is installed in the patient vicinity and a network switch or hub is used to connect it to a monitoring device, it must be correctly isolated from the mains power supply by an isolation transformer.

WARNING All external devices in the patient vicinity must comply with IEC 60601-1:1988/A1:1991A2:1995 or EN 60601-1:1990/A1:1993/A2:1995. This applies also to all signal connections, entering the patient vicinity. Additional safety equipment, e.g. isolation transformers might be used.

The installation procedures e.g. for electrical connections as documented in the Instructions for Use must be strictly followed.

Using the Monitor with an Installed, Wired Network

The following diagram shows an overview of a possible LAN installation which provides galvanic isolation of the IntelliVue monitor:

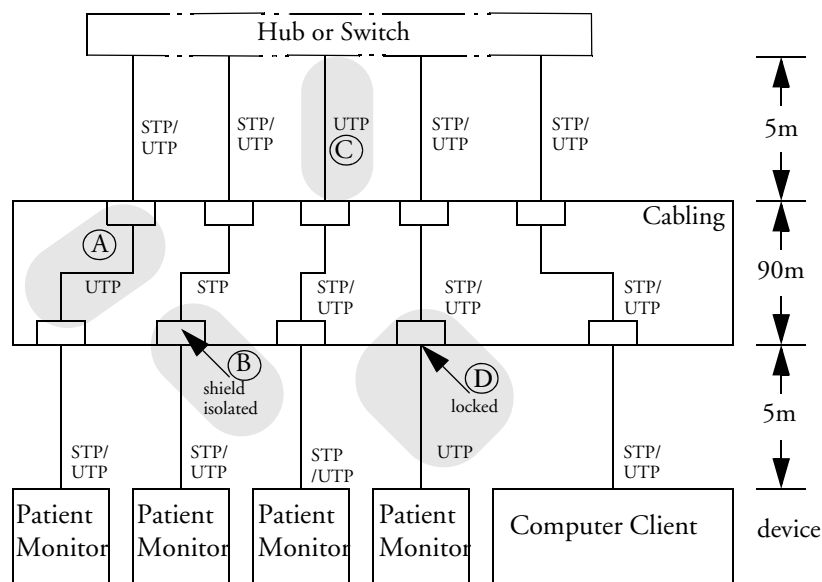If required by regulations valid in your hospital, the installation must comply to EN60601-1-1:1993/A1:1996 or IEC 60601-1-1:1992/A1:1995.

The maximum cable length between the IntelliVue monitor and the Computer Client should never exceed 330ft (100m) in total.

Note regarding MP2 and X2:

The MP2 and X2 allow Data Export via the LAN interface only because they do not have an RS232 port. The LAN interface is only available if the MP2 or X2 is used in combination with the M8023A External Power Supply. If the X2 is connected to a host monitor, the data export can be performed via the host monitor.

IntelliVue Rev. G.0 or higher allows the combination of a LAN interface and one MIB/RS232 port for Data Export. Only one connection is able to request wave data at a time, the other connection responds with a notification that wave polling is not possible. The first connection to request a successful wave poll receives the wave data.

WARNING In order to maintain the galvanic isolation of the IntelliVue monitor, it is essential that the shield is not connected from the IntelliVue monitor through to the hub or switch. At least one of the following precautions must be taken:

- UTP (Unshielded Twisted Pair) LAN cables are used in the wall.
 - If STP (Shielded Twisted Pair) LAN cables are used in the wall, do not connect the shield of the cable from the IntelliVue monitor to the wall socket. Ensure that the shield of the STP cable in the wall is isolated from the other contacts. For a reference voltage of 250V, a clearance of at least 2.5 mm and a creepage distance of at least 4.0 mm is required. Cutting the shield back and covering it with a nonconducting shroud will fulfill this requirement.
 - Ensure that only UTP cables are used in the wiring closet for connections to the hub or switch.
 - Use only UTP cables such as M3199AI #J10/J11/J12 to connect the IntelliVue monitor to the wall socket. To avoid these cables being replaced by non-UTP cables, the connector which goes into the wall socket must be modified so that it cannot be removed without using tools. This can be done by cutting off the part of the plug lock which normally extends beyond the socket.
-

Configuring the LAN Interface

Configuring the Network Address

No explicit configuration of the network addresses (IP addresses etc.) is required. The IntelliVue monitor uses the standard BootP protocol to acquire an IP address and subnet mask from a BootP server in the network. If you are using a DHCP server, make sure the server supports BootP clients.

Without a working BootP/DHCP server in the network, the IntelliVue monitor will show a technical alarm (INOP) “Unsupported LAN”, indicating that no (valid) IP address has been received.

- NOTE**
- For IntelliVue Software Revision E and later it is possible to manually enter the IP address used by the patient monitor by entering service mode and accessing the bed information window from the main setup menu.
 - With IntelliVue Software Revision H and later, use of the DHCP protocol is also supported. This requires configuration of the network interface. See the IntelliVue Configuration Guide for details.

Configuring the LAN Data Export Setting

The data that can be exported via the LAN interface is configurable. You can choose between the following options: all, anonymous data, off. In case of anonymous data the patient name and given

name are not included in the data stream.

To change the CentralMon configuration switch, first switch to configuration mode

To configure the LAN Data Export Setting, in Configuration Mode,

Select **Main Setup** to enter the Main Setup menu.

Select **Global Settings**

Select **LAN Data Export** and toggle the appropriate setting.

Configuring the Network Setting

The Central Monitoring setting on the IntelliVue monitor determines whether the monitor requires a connection to the Philips Information Center (central station). If **Central Monitoring** is set to **Mandatory**, the monitor issues a technical alarm (INOP) if a network is detected without an Information Center (central station). If you are connecting the IntelliVue monitor to a Computer Client, **Central Monitoring** should be set to **Optional**.

To do this, in Configuration Mode,

- 1 Select **Main Setup** to enter the Main Setup menu.
- 2 Select **Network**
- 3 Select **Central Monitoring** and toggle to the appropriate setting:

Mandatory The IntelliVue monitor should be connected to an Information Center.
An INOP is displayed if no connection is available.


Optional The IntelliVue monitor can be connected to an Information Center.
An INOP is only displayed if the connection to the Information Center is lost.
No INOP is displayed if no connection is found at power on.

- 4 After the configuration, make sure you have stored all the active settings and leave Configuration Mode. You do not need a password to return to Monitoring Mode.

For further details on configuration, please refer to the IntelliVue configuration guide (M8000-9306X).

Connecting to the IntelliVue MP20-90 or MX Series MIB/RS232 Interface

The IntelliVue monitor MIB/RS232 interface provides an eight-pin RJ-45 modular jack. (MP70 and MX Series RS-232 connector shown as an example below).

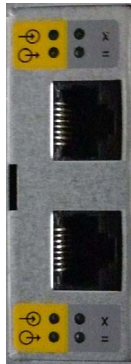MP20-90

MX Series


For the cable connection an eight conductor #24 American Wire Gauge (AWG) unshielded twisted-pair (UTP) cable must be used. The cable must follow ANSI/TIA/EIA-568-A-1995 Category 5 (CAT-5). The cable length must not exceed 65ft (20m). Straight-through pinning must be used.

The physical specification of the MIB/RS232 Interface follows the standard IEEE 1073.3.2. Refer to the standard for more information on cables and pin assignment.

The MIB/RS232 interface provides a RS232 port with the following pin assignment. This table is valid when the MIB/RS232 Interface is in DCC (Device Communication Controller) mode (DCC LED on the MIB/RS232 board is on - see below for details).

Computer Client	Pin and Signal Direction	IntelliVue monitor
	1 <=	dDPWR
GND	4 <=>	GND
RxD	5 <=	TxD
TxD	7 =>	RxD

The pins of the RJ45 are counted from 1 for the lowest pin to 8 for the highest pin when looking at the RS232/MIB interface board.

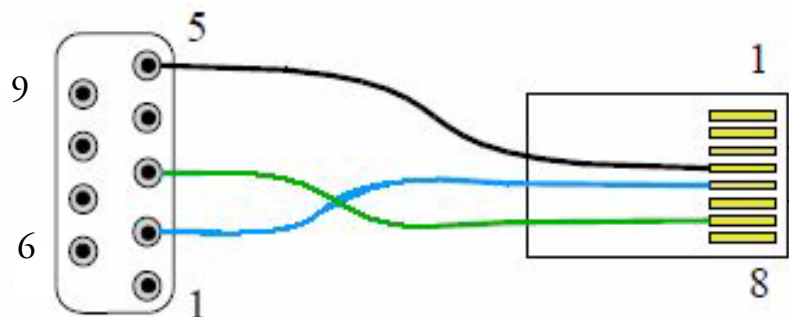LEDs on the MIB/RS232 Board (MP20-90 only)

There are four LEDs per port on the MIB/RS232 board which provide information on the configuration of the respective board. The MIB functionality is indicated by the LEDs in the yellow fields, other functionality (e.g. use for AGM or touch) is indicated by the LEDs in the grey fields. Only one LED is lit at a time.

LED	Meaning
yellow, arrow in	MIB BCC (Bedside Communication Controller) Mode
yellow, arrow out	MIB DCC (Device Communication Controller) Mode
grey, =	RS232 Mode, RX/TX lines straight
grey, X	RS232 Mode, RX/TX lines crossed

NOTE The drawings and descriptions of the RS232/MIB board above apply to the IntelliVue MP60/70 monitors. Location and orientation of the board may vary, depending on the monitor purchased.

Please note that Data Export will only function with the MIB/RS232 interface in DCC mode.

The TxD and RxD lines are the RS232 receive and transmit lines. The signals are referenced to the round (GND). The dDPWR can be used to power an external device with low power consumption. Refer to the Power Output specification in the table below.

Other applications in the IntelliVue monitor may be configured to use the MIB/RS232 Interface. These applications may use pins which are not used by the Data Export interface. Unused pins should not be connected. The IntelliVue monitor provides multiple RJ-45 connectors. Make sure, to use the correct connector with a port configured for Data Export.

The configuration of a specific MIB/RS232 port can be viewed in config mode and altered in **service mode**. To alter the configuration of an MIB port select **Main Setup** then **Hardware** then **Interfaces**. This brings up the MIB/RS232 card configuration. The port that you are using must be set to **DtOut1** for the "Data Out" function. If the MIB/RS232 port is configured for data export the yellow arrow out LED will be lit.

IntelliVue Rev. G.0 or higher allows the configuration of either a second MIB/RS232 port or the combination LAN interface and MIB/RS232 port for Data Export. The **DtOut2** driver is used to connect a second port to Data Export. Only one connection is able to request wave data at a time, the other connection responds with a notification that wave polling is not possible. The first connection to request a successful wave poll receives the wave data.

Parameter	Limit
Driver (TxD)	
Driver load output voltage (3 kOhm to 7 kOhm load)	$5\text{ V} \leq V_{out} \leq 15\text{ V}$
Driver open-circuit voltage	$ V_{out} \leq 25\text{ V}$
Driver short-circuit current (to +/- 15 V)	$ I_{osv} \leq 100\text{ mA}$
Receiver (RxD)	
Receiver input resistance	3 kOhm to 7 kOhm
Maximum receiver input voltage	+/- 25 V
Receiver threshold	+/- 3V
Power output (dDPWR)	
Minimum output voltage	4.75 V
Maximum output voltage	5.25 V
Minimum guaranteed output current	100 mA
Maximum typical output current	150 mA

REPEATED INFORMATION: If the Computer Client is not classified as a medical device, it must be located outside the patient vicinity. The patient vicinity is defined as an area within 6ft (1.85m) of the perimeter of the patient's bed or within 7.5ft (2.3m) of the floor.

WARNING All external devices in the patient vicinity must comply with IEC 60601-1:1988/A1:1991A2:1995 or EN 60601-1:1990/A1:1993/A2:1995. This applies also to all signal connections, entering the patient vicinity. Additional safety equipment, e.g. isolation transformers might be used.

The installation procedures e.g. for electrical connections as documented in the User's Guide must be strictly followed.

If it is installed in patient vicinity, the Computer Client, connected to the instrument, must be correctly isolated from the mains power supply by an isolation transformer. The MIB/RS232 interface provides galvanic isolation of the IntelliVue monitor from a connected device.

Connecting to the Intellivue MP5 Monitor MIB/RS232 Interface

The physical specification of the MP5 RS232 Interface follows the standard IEEE 1073.3.2. Refer to the standard for more information on cables and pin assignment. Note that the MP5 monitor's RS232 interface is always configured as a BCC device.


The MP5 RS232 interface provides an RS232 port with the following pin assignment.

Computer Client	Pin and Signal Direction	IntelliVue monitor
GND	4 <=>	GND
TxD	5 =>	RxD
RxD	7 <=	TxD

The TxD and RxD lines are the RS232 receive and transmit lines. The signals are referenced to the ground (GND).

To connect a PC via RS232 to the MP5 monitor, use a cable configured as shown below.

The pins of the RJ45 connector are counted from 1 for the highest pin to 8 for the lowest pin looking directly at the pins with the cable leaving the connector to the left.

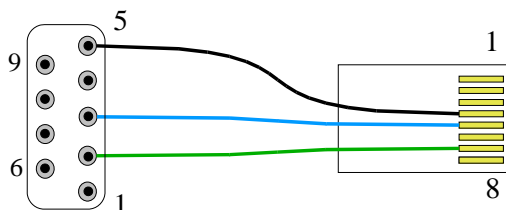The MIB/RS232 port must be set to **DtOut1** for the "Data Out" function. See page 19 for details on how to change the configuration of the MIB/RS232 port.

REPEATED INFORMATION: If the Computer Client is not classified as a medical device, it must be located outside the patient vicinity. The patient vicinity is defined as an area within 6ft (1.85m) of the perimeter of the patient's bed or within 7.5ft (2.3m) of the floor.

WARNING All external devices in the patient vicinity must comply with IEC 60601-1:1988/A1:1991A2:1995 or EN 60601-1:1990/A1:1993/A2:1995. This applies also to all signal connections, entering the patient vicinity. Additional safety equipment, e.g. isolation transformers might be used.

The installation procedures e.g. for electrical connections as documented in the User's Guide must be strictly followed.

If it is installed in patient vicinity, the Computer Client, connected to the instrument, must be correctly isolated from the mains power supply by an isolation transformer. The MIB/RS232 interface provides galvanic isolation of the IntelliVue monitor from a connected device.

Configuring the IntelliVue Monitor MIB/RS232 Interface

The MIB/RS232 interface supports different transport protocols. To change the MIB/RS232 interface configuration, in Configuration Mode,

- 1 Select **Main Setup**
- 2 Select **Hardware**
- 3 Select **Data Export** and select the required setting:

AutoSpeed Transport protocol with baudrate negotiation, based on the IrDA protocol.

Fix 19200 Transport protocol with a fixed baudrate of 19200 baud.

Fix 115200 Transport protocol with a fixed baudrate of 115200 baud.

- 4 Exit Configuration Mode. You do not need a password to return to Monitoring Mode.

IntelliVue Rev. G.0 and higher devices that allow two MIB/RS232 interfaces have two "Data Export" options. You can select **DtOut1** and/or **DtOut2**. These can be configured in service mode only. Please refer to the respective service guide for information on how to access service mode and the required password.

For further details on configuration, please refer to the IntelliVue configuration guide (M8000-9306X).

Protocol Concept

The Protocol is based on a Client/Server Model. The Personal Computer (*Client*) maintains a logical connection with the Philips IntelliVue Series Patient Monitor (*Server*). Communication occurs by sending and receiving Command messages.

Supported Transport Protocols

The Data Export functionality in the IntelliVue monitor can be accessed via the LAN interface or via the MIB/RS232 interface. While the Association Control and Data Export Protocol is the same for both interfaces, the underlying transport protocol varies.

- For the LAN interface the transport protocol is the standard UDP/IP protocol.
- For the MIB/RS232 interface, two transport protocols are supported:
 - a fixed baudrate protocol at 19200 or 115200 baud and
 - a protocol with baudrate negotiation (Auto Speed) based on the IrDA protocol with a baudrate from 9600 baud to 115200 baud.

Association Control and Data Export Protocol		
UDP/IP	RS232	RS232
	Fixed Baudrate	Auto Speed
LAN Interface	MIB/RS232 Interface	

UDP/IP Protocol

The transport protocol uses the Universal Datagram Protocol/ Internet Protocol (UDP/IP). The protocol is based on the Request For Comment (RFC) internet standard. UDP is defined in RFC 768; IP is defined in RFC 760.

The UDP/IP transport protocol is part of the internet protocol suite. Drivers and necessary hardware are available for all relevant computing platforms. It provides for a simple exchange of messages (Datagrams) across a Local Area Network. The maximum size of user data in a protocol message can be negotiated at connection time between the IntelliVue monitor and the Computer Client.

Fixed Baudrate Protocol

The Fixed Baudrate Protocol provides a transport protocol with minimal overhead and complexity. It is intended for Computer Clients which cannot use the Auto Speed Protocol. The protocol operates at a fixed baudrate and can be used with standard RS232 concentrators. It provides packet-oriented data exchange and checksum protection on top of the RS232 protocol. For the specification of the Fixed Baudrate Protocol see “Transport Protocols for the MIB/RS232 Interface” on page 30.

Auto Speed Protocol

The Auto Speed Protocol is based on the IrDA protocol. It offers a reliable transport layer with checksum protection and a retry mechanism in the case of transmission problems. The baudrate can be negotiated in a range from 9600 baud to 115200 baud. For the specification of the AutoSpeed Protocol see “Transport Protocols for the MIB/RS232 Interface” on page 30.

Protocol Model

The protocol is based on an object-oriented modelling concept. All information available through the Data Export Protocol is modelled as attribute values of information objects.

The following information object classes are supported by the IntelliVue monitor:

- Medical Device System (MDS)
The MDS object contains attributes representing dynamic state information (e.g. current operating mode) and static device specific identification information (e.g. Serial Numbers).
- Alert Monitor
The Alert Monitor object contains attributes representing the current technical and patient alarms, as e.g. displayed on the IntelliVue monitor.
- Numeric
Numeric objects contain attributes representing the state and value of numerical measurements (e.g. Heart Rate).
- Waves
Realtime sample array objects contain attributes representing the state and value of wave data (e.g. ECG).
- Patient Demographics
The Patient Demographics object contains attributes representing patient information stored in the IntelliVue monitor (e.g. Patient Name).

The object attributes can be accessed by a poll of the MDS object, which allows a query of the sets of attribute values from all objects of a specified class.

The method can be called by sending a command message from a Computer Client to the IntelliVue monitor.

Protocol Dialog

The following diagram shows the protocol dialog between the IntelliVue monitor Data Export server and a Computer Client:

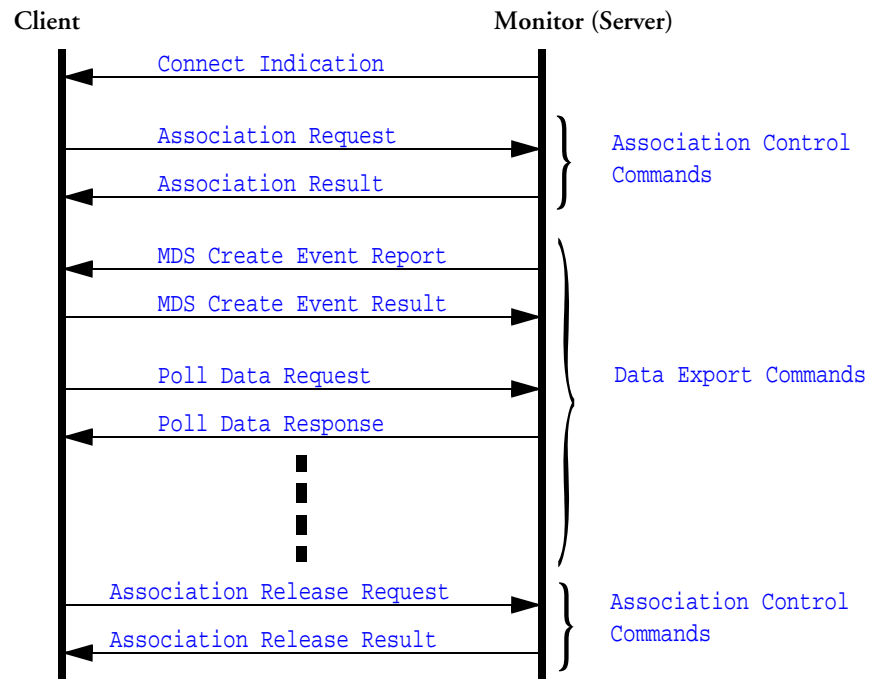
Figure 1 Protocol Dialog

Transport protocol-specific messages are not shown in the diagram. The Connect Indication message is only available on the LAN interface.

The Philips IntelliVue Series Patient Monitor processes global commands and sends response messages to the requests made by the client of the personal computer.

The messages shown in the diagram are explained in the following sections.

Connect Indication

The Connect Indication message is only sent on the LAN interface.

As soon as the IntelliVue monitor has received a valid IP address from the BootP server in the network, it sends out the Connect Indication message on its LAN interface. The message is a periodic subnet broadcast message that allows Computer Clients to find the IntelliVue monitor on the network. The message contains a set of device-related information, e.g., serial numbers, network addresses, internal states.

The IntelliVue monitor resends the Connect Indication message until a logical connection to a central station has been established. The IntelliVue monitor uses the retransmit strategy described in RFC 951. The resend period starts with 4 seconds and is doubled with each resend. The maximum resend period is about 64 seconds. The actual resend period contains a random component to avoid network congestion e.g. after a power failure.

Association Request

To establish a logical connection, the Computer Client sends the Association Request message to the IntelliVue monitor.

The Association Request can be used to set optional features of the logical connection between Computer Client and IntelliVue monitor.

Association Result

The IntelliVue monitor processes the Association Request and sends an Association Result. The result can be either a refuse message or an accept message.

The Computer Client must parse the Association Result to find out which protocol features can be used for this association.

MDS Create Event Report

If the IntelliVue monitor accepts the association, it sends a MDS Create Event Report after the positive Association Result message.

The MDS Create Event Report contains information about the system and its configuration.

MDS Create Event Result

The Computer Client must confirm the reception of the MDS Create Event Report. If the IntelliVue monitor does not receive a MDS Create Event Result message, the association is aborted.

Poll Data Request

After establishing an association, the Computer Client can send Poll Data Requests to access the data within the IntelliVue monitor.

The Poll Data Request contains a data-type parameter, which defines the specific type of requested data. The following data types are supported:

- Numeric Measurements
- Wave data
- Alerts (patient alarms and technical alarms)
- Patient Demographics
- System Attributes (e.g. dynamic state information, serial numbers, versions, etc.)

Only one type of data can be accessed per Poll Data Request.

Poll Result

Depending on the status of the IntelliVue monitor and the options set during the establishment of the logical connection (association phase), a Poll Data Request message can return:

- a single Poll Data Reply
- multiple, linked Poll Data Replies, if the size of the requested data exceeds the maximum size of a transport layer message
- a continuous number of periodic Poll Data Replies for a time period defined by the Computer Client (supported for Numeric Measurements, Waves and Alerts only).

Association Release Request

When the Computer Client wants to close an association, it can send a Release Request.

Association Release Result

The IntelliVue monitor parses the Release Request. If the Release Request is syntactically correct, the IntelliVue monitor sends an Association Release Result, indicating that the Association has been released.

Association Abort

In the case of communication problems, such as time-out, the IntelliVue monitor can send an Association Abort message. This message indicates that the association has been closed. A Computer Client should use the Association Release Request which provides a confirmation.

More Information

- For more details on the association control commands, such as Association Request, Association Result, Association Abort etc., please refer to the section “Definition of the Association Control Protocol” on page 65.
- For more detail on the data export commands, such as Poll Data Request, MDS Create Event Report, MDS Create Event Result, etc., please refer to the section “Definition of the Data Export Protocol” on page 35.

Connection Time-out Mechanism

The IntelliVue monitor automatically closes the connection if it detects a connection time-out condition. The connection time-out value is derived from the minimum poll period that is negotiated during the connection establishment phase.

A connection time-out period is 3 times the negotiated minimum poll period time. However, the minimum connection time-out is 10s, the maximum connection time-out period is 130s.

If the IntelliVue monitor does not receive a protocol message within the connection time-out period, the device closes the connection to the Computer Client by sending an Association Abort message. After that, a new connection can be established from the Computer Client to the IntelliVue monitor.

Network Load Consideration

Input Data

The IntelliVue monitor accepts a specific amount of input data per association. If the Computer Client sends more than the specified number of messages, the IntelliVue monitor will discard messages to avoid an unreasonably high system load. A Computer Client should be able to handle the loss of messages.

Message Type	Messages per Second
Association Control	1
Poll Request - Numerics (observed values)	1
Poll Request - Numerics (other attributes)	1
Poll Request - Enumerations	1
Poll Request - Waves	1
Poll Request - Alert Monitor	1
Poll Request - Patient Demographics	1
Poll Request - Medical Device System	1

The IntelliVue monitor will send a Remote Operation Error message if it receives a poll request for an object while it is still processing another poll request for the same object.

Output Data

The IntelliVue monitor processes the received message and sends the corresponding results. In rare cases, it can take up to several seconds until the response message is returned, and Poll Requests may be lost.

To avoid poll requests or poll responses getting lost, it is strongly recommended that the Computer Client uses the extended poll method to poll real-time numerics.

Definition of the Transport Protocols

Transport Protocols for the LAN Interface

UDP/IP

The Protocol uses the Universal Datagram Protocol/ Internet Protocol (UDP/IP) as the transport protocol. The protocol is based on the following internet standards (Request For Comment, RFC):

UDP is defined in RFC 768.

IP is defined in RFC 760.

The UDP/IP transport protocol is part of the internet protocol suite. Drivers and necessary hardware are available for all relevant computing platforms.

It provides for a simple exchange of messages (Datagrams) across a Local Area Network.

The maximum size of user data in a protocol message can be negotiated at connection time between the IntelliVue monitor and the Computer Client.

The upper limit for the negotiated user data size (MTU, Maximum Transport Unit) is 1364 bytes, the lower limit for the negotiated MTU is 500 Bytes. The maximum size of a UDP message sent by the IntelliVue monitor is 1380 bytes.

IP Address

The IP Address and the subnet mask necessary for communicating with the IP Protocol is set using the BootP protocol defined in the Internet RFC 951.

In order to communicate with the Philips IntelliVue Series Patient Monitor, a BootP server must exist in the network. The BootP server must be configured so that it answers BootP Request messages from the IntelliVue monitor.

UDP Port Number

The UDP Port Number used by the IntelliVue monitor for the Protocol can be extracted from the Connect Indication broadcast message used for Device Discovery (see “CONNECT INDICATION EVENT” on page 53). The current Protocol version uses the fixed UDP port 24105.

All messages sent from the Computer Client to the IntelliVue monitor must use this port number as the destination port number.

The Computer Client can choose any available source port for the communication. Once the Computer Client has chosen a source port, it must not use any other port. Protocol messages from another source port will be regarded as messages from a different Computer Client).

Any messages sent from the IntelliVue monitor back to the Computer Client use the source port number set by the Computer Client in first message (the Association Request message, see “Association Request Message” on page 67).

Transport Protocols for the MIB/RS232 Interface

The Fixed Baudrate Protocol, RS232 Port Settings

Each transmitted byte consists of one start bit, 8 data bits (no parity) and one stop bit. The baudrate can be set to 115kBit/s or 19.2kBit/s.

Flow control is not supported (same behavior as UDP). The monitor limits the number of Frames which will be processed in a given time. The monitor will process up to 4 frames within 128ms. If a client sends more frames, additional frames are ignored. (Implementation Note: the monitor allows 5 frames within 128ms, the additional frame is required because of possible jitter.)

A client system must be able to handle the loss of messages, because the Fixed Baudrate Protocol does not guarantee the reliable transmission of messages.

Framing

BOF	Hdr	User Data	FCS	EOF
-----	-----	-----------	-----	-----

The framing structure is the same as for AutoSpeed protocol. A frame starts with a single BOF.

BOF	Beginning Of Frame (0xC0)
Hdr	Header Information
User Data	Association Control or Data Export Command message
FCS	16 bit Frame Check Sequence using CRC-CCITT algorithm
EOF	End Of Frame (0xC1)

Header Information

The *Hdr* field is defined as follows:

```
typedef struct
{
    u_8      protocol_id;
    u_8      msg_type;
    u_16     length;
} FrameHdr;
```

The *protocol_id* field contains ID and version information. It can be used to define different service access points. Data Export uses the ID 0x11.

The *msg_type* field defines the type of message which is being sent. The value 0x01 indicates an Association Control or Data Export Command message, future message types could be used for flow control, lifetick, message confirmation etc.

The *length* field contains the length of the appended user data in bytes (without transparency characters).

If a client receives messages with an unknown *protocol_id* or *msg_type*, it should ignore the message.

Frame Check Sequence Field

The Frame Check Sequence Field can be used to detect transmission errors. The field contains a 16 bit CRC-CCITT cyclic redundancy check (not the popular XMODEM variation of CRC-CCITT). The CRC is computed from the *Hdr* and *User Data* field. Refer to "Serial Infrared Link Access Protocol (IrLAP)" Version 1.1 for the actual computation method of the CRC. A code snippet for the FCS algorithm can be found in the Network Working Group Request for Comment: 1171 (PPP protocol). The one's complement of the CRC is transmitted, rather than the CRC itself. The CRC is transmitted LSB first.

If the CRC is not correct, a client system should ignore the message.

Transparency

The contents of the *Hdr* and *User Data* fields is unrestricted. This can lead to problems if a BOF or EOF character appear in the *Hdr*, *User Data*, or FCS field. A Control Escape byte is defined as 0x7D. The sender must examine each byte in the *User Data* and FCS fields; for each byte with the value 0xC0, 0xC1, 0x7D it does the following:

- insert a 0x7D byte proceeding the byte
- complement bit 5 of the byte (XOR with 0x20).

Frame Abort

The sending station may abort the transmission of a frame by sending a control escape character followed by a EOF character (0x7D 0xC1) without sending the FCS field.

Examples The examples below do not include the *Hdr* field. For a correct message, the framing algorithm must be applied to the *Hdr* and *UserData* field of the message.

1 If a Computer Client wants to send the data:

"0x3a 0x71"

The CRC for this data would be:

"0xd9 0x64"

after building the one's-complement and byte-swapping, this results in:

"0x9b 0x26"

The whole frame would be:

"0xc0 0x3a 0x71 0x9b 0x26 0xc1"

2 If a Computer Client wants to send the data:

"0x3a 0x91"

The CRC for this data would be:

“0x3e 0x6a”

after building the one's-complement and byte-swapping, this results in:

“0x95 0xc1”

The whole frame would be:

“0xc0 0x3a 0x91 0x95 0x7d 0xe1 0xc1”

Note that byte “0xc1” in the CRC is a reserved character and must be escape. This results in “0x7d 0xe1”.

The AutoSpeed Protocol

The AutoSpeed Protocol follows the definition of the Transport Protocol defined in the standard IEEE 1073.3.2: IEEE Standard for Medical Device Communications - Transport Profile - IrDA Based Cable Connection.

For a description of the IrDA Protocol refer to the specifications of the Infrared Data Association (www.irda.org):

- IrDA, Serial Infrared Link Access Protocol (IrLAP), Version 1.1, June 16, 1996
- IrDA, Link Management Protocol (IrLMP), Version 1.1, Oct. 20, 1996
- IrDA, Tiny TP: A Flow-Control Mechanism for use with IrLMP, Version 1.1, Oct. 20, 1996

Commercial IrDA stacks are available for most operating systems. Some operating systems, like Microsoft® Windows 2000® and Linux, come with an off-the-shelf IrDA stack.

The Data Export protocol resides as a packet oriented client on top of the IrDA TinyTP layer.

Establishing a Connection

A connection is created using the following steps:

- Discovery

The Computer Client sends an IrLAP discovery request to find out if a device is physically connected. The IntelliVue monitor answers with an discovery response message. The discovery procedure is done at a fixed baudrate of 9600 baud.

- Open an IrLAP connection

When the Computer Client finds a connected system, it can send an IrLAP Set Normal Response Mode message to establish a logical IrLAP connection. The IntelliVue monitor sends an response message. During this procedure parameters of the IrLAP connection, like baudrate, data size, etc. are negotiated.

- Open an IAS port

The Information Access Service (IAS) is provided by the IrLMP layer. It provides a database with device information which can be queried by the client. Before accessing the service, the client must connect to the special IrLMP service access point (SAP) 0.

- Perform an IAS query

The IrLMP layer does not specify a well-known SAP for the Data Export Protocol, hence the client should query the IAS database to find the SAP for the Data Export Protocol. The database contains the attribute “IrDA:TinyTP:LsapSel” under the object class “IEEE:1073:3:2:MDDL”. The attribute specifies the SAP for the Data Export Protocol on the IrDA TinyTP layer as an integer value.

- Close the IAS port

After performing the IAS query, the Computer Client should close the IAS port again with an IrLMP disconnect message.

- Open a Tiny TP connection

After retrieving the number for the TinyTP SAP, the client system can open a connection on this SAP. This is done with an IrLMP connect request message which contains a TinyTP connect in its user data.

- Send an Association Request

After the transport layer connection has been established, the Computer Client can send an Association Request message to start a Data Export session.

- Send a Release Request

When the client has no need for further communication, it can send a Release Request message to terminate the Data Export session.

- Close the IrLAP connection

After the Data Export session has been closed, the Computer Client should also close the TinyTP SAP. This can be done by sending an IrLMP disconnect message or by closing the whole IrLAP connection.

Definition of the Data Export Protocol

Definitions Shared by Protocols

Byte Order

The protocol data structures use the Network Byte Order. This means that bytes of a multi-byte data structure are transmitted on the network with the most significant byte first (as in big-endian data storage). This may or may not match the order in which numbers are normally stored in memory for a particular processor.

If the Computer Client is not using big-endian storage internally (many common Personal Computer Platforms use little-endian storage), protocol data structures (message structures) must be transformed before they are sent to an IntelliVue monitor or after they have been received from an IntelliVue monitor.

Byte Alignment

The Association Control and Data Export protocols assume no data alignment. However, most data types used in this guide have an even length for performance reasons. Many compilers use different alignment modes by default. Make sure that the compiler uses the right alignment when parsing and formatting protocol messages.

Bit Order

The index for bits starts with zero for the most significant bit.

MSB										LSB					
0	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15

Common Data Types

Basic Data Types

The C data types defined here make use of the following basic types:

```
u_8      unsigned 8 bit wide integer
u_16     unsigned 16 bit wide integer
u_32     unsigned 32 bit wide integer
i_8      signed 8 bit wide integer
i_16     signed 16 bit wide integer
i_32     signed 32 bit wide integer
```

The mapping of these types to data types used in a Computer Client application is machine specific and compiler dependent.

Absolute Time

The Absolute Time data type is used whenever data is time stamped and a resolution of 1s is sufficient.

```
typedef struct {
    u_8      century;
    u_8      year;
    u_8      month;
    u_8      day;
    u_8      hour;
    u_8      minute;
    u_8      second;
    u_8      sec_fractions;
} AbsoluteTime;
```

The individual u_8 fields are BCD encoded, they are not encoded as regular integer values. E.g. the year 99 (decimal) is coded as 0x99. An invalid time is marked with 0xff in all positions.

Note that the time resolution in IntelliVue monitor with this format is 1 second. The *sec_fractions* element in the structure is not used.

Relative Time

The Relative Time is a high resolution time marker which defines a time relative to an event (e.g. power-on). It is used to position events (a particular event message) relative to each other with a higher resolution. It is defined as follows:

```
typedef u_32    RelativeTime;
```

The resolution of the *RelativeTime* is 1/8ms (125us). The IntelliVue monitor sets the Relative Time with a precision of 2 ms. The Computer Client can calculate the absolute time (wall clock) from a known relation between Absolute Time and Relative Time with a precision of about 1s. For more information on the time mapping refer to “MDS CREATE EVENT” on page 54.

OID Type

For the identification of all protocol elements (e.g. physiological meaning, alert codes, units of measure), the *OIDType* (Object Identifier Type) is used.

```
typedef u_16          OIDType;
```

Values for the *OIDType* (the nomenclature) are listed at the end of the section “Attribute Data Types and Constants Used” on page 75. Independent value ranges (partitions) exist, e.g. for physiological identifiers, alert condition identifiers, units of measurement etc.

Private OID

For the identification of private or manufacturer specific elements, a special type is used.

```
typedef u_16          PrivateOID;
```

Values for the *PrivateOIDs* are listed whenever a *PrivateOID* is used. Refer to the section “Attribute Data Types and Constants Used” on page 75 for a complete list of identifiers.

TYPE

Whenever it is not clear from the context, from which nomenclature value range the *OIDType* comes, the TYPE data type is used. Here, the nomenclature value range (the partition) is explicitly identified.

```
typedef u_16          NomPartition;
#define NOM_PART_OBJ          1
#define NOM_PART_SCADA        2
#define NOM_PART_EVT          3
#define NOM_PART_DIM          4
#define NOM_PART_PGRP         6
#define NOM_PART_INFRASTRUCT  8

typedef struct {
    NomPartition  partition;
    OIDType       code;
} TYPE;
```

The *code* values are grouped in the following partitions:

- NOM_PART_OBJ: Object oriented element, device nomenclature
- NOM_PART_SCADA: Types of measurement and place of the measurement
- NOM_PART_EVT: Codes for alerts
- NOM_PART_DIM: Units of measurement
- NOM_PART_PGRP: Identification of parameter groups
- NOM_PART_INFRASTRUCT: Infrastructure for Data Export applications

The *code* is only unique in a given partition. The values for the *OIDType* are defined in the section “Attribute Data Types and Constants Used” on page 75.

Handle

Object instances, e.g. Numeric object instances, are identified with a 16bit wide ID, the object Handle:

```
typedef u_16          Handle;
```

Global Handle

Handles are unique within the context of a particular system. The Protocol supports multiple measurement servers, where each measurement server assigns object handles independently. To assure handle uniqueness across system boundaries, the Global Handle contains an additional identifier for the source system, e.g., each measurement server has a distinct context id. The context id is assigned dynamically when a measurement server is connected.

```
typedef u_16          MdsContext;

typedef struct {
    MdsContext      context_id;
    Handle          handle;
} GlbHandle;
```

Managed Object Identifier

The Managed Object Identifier is a fully qualified object identifier which contains an identifier for the object class (e.g. Numeric object) together with a Global Handle.

```
typedef struct {
    OIDType          m_obj_class;
    GlbHandle         m_obj_inst;
} ManagedObjectId;
```

Attribute Value Assertion

Object attributes are represented in the form of data record structures which contain an identifier for the attribute, a length field for parsing and the actual value of the attribute.

The structure of such an attribute record is the Attribute Value Assertion, which is defined as follows:

```
typedef struct {
    OIDType          attribute_id;
    u_16             length;
    u_16             attribute_val;
} AVAType;
```

The *attribute_id* identifies the type of the attribute. The length field contains the size of the *attribute_val* field in bytes. The *attribute_val* field itself is only a placeholder in this structure. The parsing algorithm must assign the attribute value to the correct data structure based on the value of the *attribute_id*.

Attribute List

Typically, object instances have multiple attributes which are captured in a list with the following data type:

```
typedef struct{
    u_16          count;
    u_16          length;
    AVAType       value[1];
} AttributeList;
```

The count field contains the number of Attribute Value Assertion elements in the list.

The length field contains the size of the list (the value array) in bytes.

The value field itself again is only a placeholder data structure. A parser must be used to interpret the data structure. Refer to “Protocol Examples” on page 291 for an example of an *AttributeList*.

String

The text string is preceded by a length field, followed by the value. The *length* field denotes the number of octets in *value*. If the length is zero, no octets are appended. The *String* data type is used for UNICODE encoded texts.

```
typedef struct {
    u_16          length;
    u_16          value[1];
} String;
```

Where possible, the real string lengths have been included in this document. However, these string lengths may change in future releases, producing discrepancies between the actual string lengths and this document.

The *String* uses the same language as the IntelliVue monitor. The IntelliVue monitor uses UNICODE for the *String* data type (see “Connect Indication Attributes” on page 107). The *String* may contain code values from the UNICODE private use area (0xE000 to 0xF8FF). The Computer Client most likely will not support these characters. The following codes are frequently used:

```
#define SUBSCRIPT_CAPITAL_E_CHAR      0xE145
/* SUBSCRIPT CAPITAL E                */
#define SUBSCRIPT_CAPITAL_L_CHAR      0xE14C
/* SUBSCRIPT CAPITAL L                */
#define LITER_PER_CHAR                0xE400
/* LITER PER - used in 4 char unit "l/min" */
#define HYDROGEN_CHAR                 0xE401
/* HYDROGEN - Used in 4 char unit "cmH2O"   */
#define ALARM_STAR_CHAR                0xE40D
/* ALARM STAR ""                       */
#define CAPITAL_V_WITH_DOT_ABOVE_CHAR 0xE425
/* CAPITAL_V_WITH_DOT_ABOVE (V with dot) */
#define ZERO_WIDTH_NO_BREAK_SPACE_CHAR 0xFFEF
/* The character 0xFFEF is used as FILL character.
For each wide asian character, a FILL character is
appended for size calculations. */
```

Variable Label

The string is preceded by a length field, followed by the value. If the length is zero, no octets are appended. The *VariableLabel* data type uses 8 bit ASCII encoding for the text. The *length* of a *VariableLabel* is always even.

```
typedef struct {
    u_16    length;
    u_8     value[1];
} VariableLabel
```

Where possible, the real string lengths have been included in this document. However, these string lengths may change in future releases, producing discrepancies between the actual string lengths and this document.

TextId

The *TextId* type is a 32bit wide private ID.


```
typedef u_32    TextId;
```

FLOAT-Type

For floating point numbers, a special 32bit wide format is used. For message parsing and for the definition of the message structures a 32bit wide placeholder structure is defined here.

```
typedef u_32    FLOATType;
```

The FLOAT-Type must be interpreted as follows:

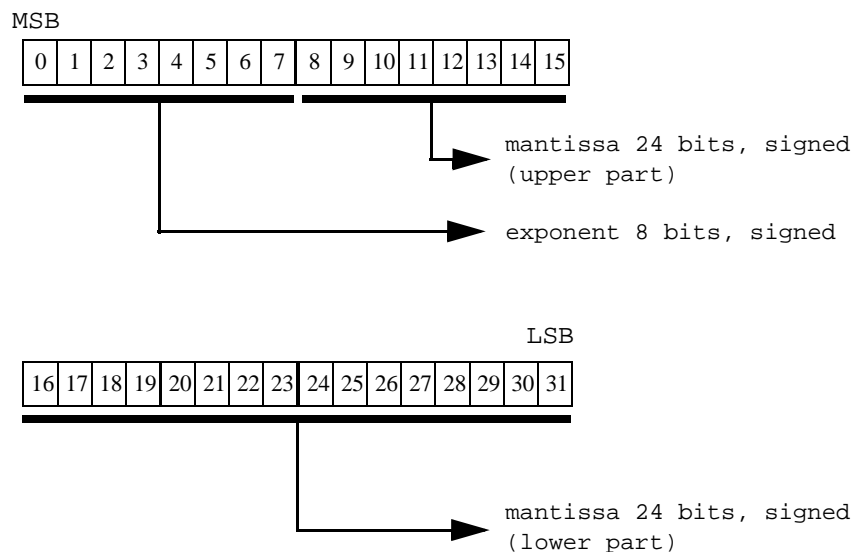The number represented is $(\text{mantissa}) \cdot (10^{\text{exponent}})$. Both the exponent and mantissa are in 2's complement form. The mantissa is not necessarily normalized.

There are four special values of the mantissa that can be represented:

NaN (Not a Number), which has a mantissa of $+(2^{23}-1)$ (0x7fffff)

NRes (Not at this resolution), which has a mantissa of $-(2^{23})$ (0x800000)

+/- INFINITY, which have mantissa of +/- $(2^{23}-2)$ (0x7ffffe, 0x800002).

The exponent is not important in these cases. This leaves the following ranges for normal number representation:

$-128 \leq \text{exponent} \leq 127$

$-(2^{23}-3) \leq \text{mantissa} \leq +(2^{23}-3)$

Definition on the number of the valid digits for the presentation on the IntelliVue monitor's display:

1.) If the exponent < 0 , then the integer value of the exponent shows the number of valid digits after the point:

Examples:

value = 0xfd007d00: exponent = -3, mantissa = 32000 ➡ 32.000

value = 0xff000140: exponent = -1, mantissa = 320 ➡ 32.0

2.) If the exponent ≥ 0 , then the number of valid digits after the point is zero.

Examples:

value = 0x01000140: exponent = 1, mantissa = 320 ➡ 3200

value = 0x02000020: exponent = 2, mantissa = 32 ➡ 3200

Protocol Command Structure

Protocol Command messages, as defined in this section, are the data structures that are transported within the transport layer message (UDP datagram, IrDA message or Fixed Baudrate Protocol message). The generic structure is common for messages sent from the Computer Client to the IntelliVue monitor (e.g. Poll Request messages) and messages sent from the IntelliVue monitor to the Computer Client (e.g. Poll Result messages).

The Protocol Command messages represent the ISO/ OSI layers 5 - 7 (session layer, presentation layer, application layer). The message that transports a Protocol Command contains a checksum. Computer Clients should validate this checksum to detect corrupted messages.

The Protocol command messages used to establish the logical connection (association) between the IntelliVue monitor and a Computer Client follow the definitions of the ACSE Standard (ISO/IEC 8649 and ISO/IEC 8650).

For the Protocol Commands during the logical connection, the message structure is layered and has the following basic format:

Session/Presentation Header
Remote Operation Header
Command Header
Command- Specific Parameter Data

The Session Header and Presentation Header are small fields only which contain fixed values for the life time of the logical connection between the IntelliVue monitor and the Computer Client.

The Remote Operation Header allows to distinguish between the different types of command messages, command response messages and error messages.

The Command Header contains the common part of the Command data structure identified in the Remote Operation Header.

Command-specific parameters or data are appended to the generic message structure.

Session/Presentation Header

Each protocol message starts with a common data structure representing the session and presentation protocol, defined as follows:

```
typedef struct {
    u_16    session_id;    /* contains a fixed value 0xE100 */
    u_16    p_context_id; /* negotiated in association phase */
} SPpdu;
```

- session_id

This field identifies a Protocol message. The field contains a fixed value 0xE100. Conceptually, this field represents the session header.

- p_context_id

The presentation context identifier is negotiated during the exchange of the association messages.

The Computer Client can use the first byte of the *session_id* to distinguish between Data Export protocol commands and Association Control protocol commands.

If a Computer Client encodes the Association Control protocol commands as suggested in “Definition of the Association Control Protocol” on page 65, the *context_id* for the Data Export protocol commands is 2.

Remote Operation Header

A protocol message is considered a remote operation. There are different types of operations as defined below. The different operations are described by a common operation header data structure:

```
typedef struct {
    u_16    ro_type;          /* ID for operation */
    #        define          ROIV_APDU 1
    #        define          RORS_APDU 2
    #        define          ROER_APDU 3
    #        define          ROLRS_APDU 5
    u_16    length;          /* bytes to follow */
} ROapdus;
```

- ro_type

This field defines which type of remote operation is appended.

The following remote operation types exist:

Remote Operation Invoke (ROIV_APDU) invokes (calls) a remote operation.

Remote Operation Result (RORS_APDU) returns the result of a remote operation

Remote Operation Error (ROER_APDU) returns an error for a remote operation.

Remote Operation Linked Result (ROLRS_APDU) returns parts of the result of a remote operation. It is used when the size of the complete result exceeds the maximum size of one message.

- length

This field defines the remaining number of bytes in the message.

Remote Operation Invoke

A Remote Operation Invoke message is defined as follows:

```
typedef struct {
    u_16    invoke_id;      /* identifies the transaction */
    CMDType command_type;   /* identifies type of command */
    u_16    length;         /* no. of bytes in rest of message */
} ROIVapdu;
```

- invoke_id

The invoke identifier is used to reference the specific operation while it is being processed. Result messages or error messages will use this identifier as a reference. Therefore, the invoke identifier should be unique while the operation transaction is in process.

- command_type

The command type identifier defines what command data type is appended to this structure.

- length

This field defines the remaining number of bytes in the message.

Remote Operation Result

A Remote Operation Result message is a response to an Operation Invoke message requiring confirmation.

The message is defined as follows:

```
typedef struct {
    u_16    invoke_id;      /* mirrored back from op. invoke */
    CMDType command_type;   /* identifies type of command */
    u_16    length;         /* no of bytes in rest of message */
} RORSapdu;
```

- invoke_id

The invoke identifier is mirrored back from the related Remote Operation Invoke message that triggered this result. This field allows to relate the response message to the original request.

- command_type

The command type identifier defines what command data type is appended to this structure.

- length

This field defines the remaining number of bytes in the message. This length is not larger than the negotiated Maximum Transport Unit (MTU). For larger messages, the Remote Operation Linked Result mechanism will be used.

Remote Operation Linked Result

In some cases, the total data that must be returned as a result of a command may exceed the maximum message size. In these cases, multiple Remote Operation Linked Result messages are used.

These messages are defined as follows:

```
typedef struct {
    RorlsId linked_id;    /* see below */
    u_16   invoke_id;    /* see below */
    CMDType command_type; /* identifies type of command */
    u_16   length;       /* no of bytes in rest of message */
} ROLRSapdu;
```

- linked_id

The linked identifier identifies each Remote Operation Linked Result message in a sequence of linked messages (see below).

- invoke_id

The invoke identifier is mirrored back from the related Remote Operation Invoke message that triggered this result. This field allows to relate the response message to the original request.

- command_type

The command type identifier defines what command data type is appended to this structure.

- length

This field defines the remaining number of bytes in the message.

If the size of the result data exceeds the maximum message size, a combination of Remote Operation Linked Result Messages and Remote Operation Result messages is used, with the following rules:

- For all response messages except the very last one:
 - the ROLRS_APDU message type is used
 - the linked identifier is set by the responder to the *RorlsId* data type
 - the invoke identifier is the value of the invoke identifier of the associated Operation Invoke
- For the very last message:
 - The RORS_APDU message type is used

The invoke identifier in this response is the value of the invoke identifier of the associated Operation Invoke.

The following data type is used for the linked identifier:

```
typedef struct {
    u_8   state;
    #     define RORLS_FIRST 1    /* set in the first message */
    #     define RORLS_NOT_FIRST_NOT_LAST 2
    #     define RORLS_LAST 3     /* last RORLSapdu, one RORSapdu
                                   to follow */
    u_8   count;                 /* counter starts with 1 */
} RorlsId;
```

The first Remote Operation Linked Result message sets the state RORLS_FIRST.

The last Remote Operation Linked Result message sets the state RORLS_LAST. Note that there is one more Remote Operation Result message to follow.

All other Remote Operation Linked Result messages set the state RORLS_NOT_FIRST_NOT_LAST.

Examples:

- If a total of 3 messages are needed, the first message is a Remote Operation Linked Result with state RORLS_FIRST and count field 1. The second message is a Remote Operation Linked Result with state RORLS_LAST and count field 2. The third message is a Remote Operation Result message.
- If a total of 2 messages are needed, the first message is a Remote Operation Linked Result with state RORLS_LAST and count field 1. The second message is a Remote Operation Result message.

The *count* field starts with 1 for the first of the linked messages and is increased with each following message.

When a message is split, each message contains a full command data structure (see “Command Header” on page 47).

If the messages contain data from several objects, the Computer Client can not assume that all data belonging to one object is sent within one message. In some cases it can happen that the data belonging to one attribute of a given object must be sent in multiple messages (see the description of the available data in the section “Attribute Data Types and Constants Used” on page 75). This may only occur for attributes which are encoded in the form of a list (e.g Device T-Alarm List).

Object data which did not fit in one message is guaranteed to continue in the next linked message.

Remote Operation Error

If an error is detected at the Remote Operation level, an error message is returned:

```
typedef struct {
    u_16    invoke_id;
    u_16    error_value;
    #       define      NO_SUCH_OBJECT_CLASS          0
    #       define      NO_SUCH_OBJECT_INSTANCE       1
    #       define      ACCESS_DENIED                 2
    #       define      GET_LIST_ERROR                7
    #       define      SET_LIST_ERROR                8
    #       define      NO_SUCH_ACTION                 9
    #       define      PROCESSING_FAILURE            10
    #       define      INVALID_ARGUMENT_VALUE        15
    #       define      INVALID_SCOPE                 16
    #       define      INVALID_OBJECT_INSTANCE       17
    u_16    length;
} ROERapdu;
```

- *invoke_id*

The invoke identifier is mirrored back from the related Remote Operation Invoke message that triggered this result. This field allows to relate the response message to the original request.

- error_value

The error values have the following meaning:

GET_LIST_ERROR: Get operation failed. A *GetListError* is appended to the message.

SET_LIST_ERROR: Set operation failed. A *SetListError* is appended to the message.

NO_SUCH_ACTION: Unknown action type. The object class ID and action type are appended to the message.

NO_SUCH_OBJECT_CLASS: There is no such object class in the system. An *OIDType* with the class ID is appended to the message.

NO_SUCH_OBJECT_INSTANCE: The object instance does not exist. The *ManagedObjectId* of the instance is appended.

ACCESS_DENIED: Computer Client has not required privileges to perform the operation. No data is appended.

PROCESSING_FAILURE: Generic error indicating an invalid request. A *ProcessingFailure* is appended to the message.

INVALID_ARGUMENT_VALUE: The argument of the ROSE message was not valid. An Action result is appended.

INVALID_SCOPE: The scope is not valid for the operation. The value of the scope is appended.

INVALID_OBJECT_INSTANCE: Wrong object instance. The *ManagedObjectId* of the instance is appended.

- length

This field defines the remaining number of bytes in the message.

The *GetListError* and *SetListError* structures are defined as follows:

```
typedef struct {
    ManagedObjectId    managed_object;
    struct {
        u_16           count;
        u_16           length;
        GetError        value[1];
    } getInfoList;
} GetListError;

typedef struct {
    ErrorStatus        errorStatus;
    OIDType            attributeId;
} GetError;

typedef struct {
    ManagedObjectId    managed_object;
    struct {
        u_16           count;
        u_16           length;
        SetError        value[1];
    } setInfoList;
} SetListError;

typedef struct {
    ErrorStatus        errorStatus;
    ModifyOperator      modifyOperator;
    OIDType            attributeId;
} SetError;

typedef u_16          ErrorStatus;
#define ATTR_ACCESS_DENIED      2
#define ATTR_NO_SUCH_ATTRIBUTE  5
#define ATTR_INVALID_ATTRIBUTE_VALUE 6
```

```
#define ATTR_INVALID_OPERATION      24
#define ATTR_INVALID_OPERATOR      25
```

The *ProcessingFailure* is defined as follows:

```
typedef struct {
    OidType      error_id;
    u_16         length;
} ProcessingFailure;
```

Additional data with error information can be appended to the *ProcessingFailure*. The default *error_id* is 0 with no appended data.

Command Header

In each protocol message, a Command data structure is appended. The specific Command is identified by the value of the *CMDType* field in the Remote Operation Invoke/ Result/ Linked Result data structures.

The following Command types are used in the Protocol:

```
typedef u_16      CMDType;
#define           CMD_EVENT_REPORT      0
#define           CMD_CONFIRMED_EVENT_REPORT 1
#define           CMD_GET               3
#define           CMD_SET               4
#define           CMD_CONFIRMED_SET     5
#define           CMD_CONFIRMED_ACTION  7
```

The following command types are used:

CMD_EVENT_REPORT: An Event Report is used for an unsolicited event message.

CMD_CONFIRMED_EVENT_REPORT: The Confirmed Event Report is an unsolicited event message for which the receiver must send an Event Report Result message.

CMD_GET: The Get operation is used to request attribute values of managed objects. The receiver responds with a Get Result message.

CMD_SET: The Set operation is used to set values of managed objects.

CMD_CONFIRMED_SET: The Confirmed Set operation is used to set attribute values of managed objects. The receiver responds with a Set Result message.

CMD_CONFIRMED_ACTION: The Confirmed Action is a message to invoke an activity on the receiver side. The receiver must send an Action Result message.

For confirmed messages, the receiver must send the appropriate result message. For both the confirmed and unconfirmed Event Report, an *EventReportArgument* is appended.

If the result message is not received within 3 seconds, the IntelliVue monitor resends the message. If the message has not been confirmed after sending it 3 times (2 resend tries), the association is aborted by the IntelliVue monitor.

Event Report

The Event Report command (CMD_EVENT_REPORT) is used for unsolicited messages from the sending device to the receiving device. It is appended to the Remote Operation Invoke message. In the Protocol the Event Report may require a response from the receiver (if a response is required, the CMD_CONFIRMED_EVENT_REPORT Command identifier is used).

The Event Report message uses the following data structure:

```
typedef struct {
    ManagedObjectId managed_object; /* ident. of sender */
    RelativeTime    event_time;    /* event time stamp */
    OIDType         event_type;    /* identification of event */
    u_16            length;        /* size of appended data */
} EventReportArgument;
```

- **managed_object**
Identifies the object that generates the unsolicited Event Report command.
- **event_time**
The relative time (in 1/8ms time ticks) of the event.
- **event_type**
Identifies the event type and thus the data structure that is appended.
- **length**
This field defines the remaining number of bytes in the message (which is the size of the event specific data appended to this data structure).

Event-specific data is appended to the data type.

Event Report Result

The Event Report Result command is used as a response message to the Event Report message. It is appended to the Operation Result message with the *command_type* CMD_CONFIRMED_EVENT_REPORT.

The Event Report Result uses the following data structure:

```
typedef struct {
    ManagedObjectId managed_object; /* mirrored from EvRep */
    RelativeTime    current_time;  /* result time stamp */
    OIDType         event_type;    /* identification of event */
    u_16            length;        /* size of appended data */
} EventReportResult;
```

- **managed_object**
Identifies the object to which the response is sent back. This field must be mirrored back from the Event Report message.
- **event_time**
The relative time (in 1/8ms time ticks) of the event result.
- **event_type**
Identifies the event type and thus the data structure that is appended. This field must contain the same value as the Event Report.
- **length**
This field defines the remaining number of bytes in the message (which is the size of the event specific result data appended to this data structure).

Event-specific data is appended to the data type.

Action

The ACTION command (CMD_CONFIRMED_ACTION) is used to call a Protocol specific method in the receiver. The Protocol uses this command to call the *Data Poll* method which returns device data. The ACTION command is appended to the Operation Invoke message.

The Action command uses the following data structure:

```
typedef struct {
    ManagedObjectId managed_object; /* addressed object */
    u_32             scope;         /* fixed value 0 */
    OIDType          action_type;   /* identification of method */
#define NOM_ACT_POLL_MDIB_DATA      3094
#define NOM_ACT_POLL_MDIB_DATA_EXT 61755
    u_16             length;        /* size of appended data */
} ActionArgument;
```

- **managed_object**
Identifies the object to which the ACTION command is sent.
- **scope**
Contains a fixed value 0 in this version of the protocol.
- **action_type**
Identifies the specific method that should be called (and thus the data type that is appended to this data structure).
NOM_ACT_POLL_MDIB_DATA is used for a Single Poll Data Request.
NOM_ACT_POLL_MDIB_DATA_EXT is used for an Extended Poll Data Request
- **length**
This field defines the remaining number of bytes in the message (which is the size of the method specific data appended to this data structure).

Method-specific data is appended to the data type.

Action Result

The Action Result command is used as a response message to the Action message. It is appended to the Operation Result message or an Operation Linked Result message (if the size of the returned data exceeds a maximum message size). The *command_type* is set to CMD_CONFIRMED_ACTION.

The Action Result uses the following data structure:

```
typedef struct {
    ManagedObjectId managed_object;
    OIDType          action_type; /* identification of method */
    u_16             length;      /* size of appended data */
} ActionResult;
```

- **managed_object**
Identifies the object that responds to the ACTION command (usually mirrored from ACTION command).
- **action_type**
Identifies the specific method that was called (and thus the data type that is appended to this data structure).

- length
This field defines the remaining number of bytes in the message (which is the size of the method specific result data appended to this data structure).
- Method-specific data is appended to the data type.

Get

The Get command (CMD_GET) specifies attributes that should be returned. It is appended to an Operation Invoke message.

The Get command uses the following data structure:

```
typedef struct {
    ManagedObjectId    managed_object;
    u_32               scope;
    AttributeIdList    attributeIdList;
} GetArgument;
```

- managed_object
Identifies the object to which the Get command is sent.
- scope
Contains a fixed value 0 in this version of the protocol.
- attributeIdList
Contains the list of attribute identifiers.

```
typedef struct {
    u_16               count;
    u_16               length;
    OIDType            value[1];
} AttributeIdList;
```

Get Result

The Get Result is returned in response to the Get command. It is appended to an Operation Result or Operation Linked Result message.

The Get Result uses the following data structure:

```
typedef struct {
    ManagedObjectId    managed_object;
    AttributeList      attributeList;
} GetResult;
```

- managed_object
Identifies the object that responds to the Get command.
- attributeList
Contains the requested attributes.

Set

The Set command (CMD_SET) or Confirmed Set command (CMD_CONFIRMED_SET) specifies attributes that should be added, replaced, or removed. It is appended to an Operation Invoke message.

The Set command uses the following data structures:

```
typedef struct {
    ManagedObjectId    managed_object;
    u_32               scope;
```

```

        ModificationList      modificationList;
    } SetArgument;

```

- **managed_object**
Identifies the object to which the Get command is sent.
- **scope**
Contains a fixed value 0 in this version of the protocol.
- **modificationList**
Contains the attribute ids and values to be modified.

```

typedef struct {
    u_16          count;
    u_16          length;
    AttributeModEntry value[1];
} ModificationList;

typedef struct {
    ModifyOperator      modifyOperator;
    AVAType             attribute;
} AttributeModEntry;

typedef u_16          ModifyOperator;
#define REPLACE                0
#define ADD_VALUES             1
#define REMOVE_VALUES         2
#define SET_TO_DEFAULT        3

```

Set Result

The Set Result is returned in response to the Confirmed Set command. It is appended to an Operation Result or Operation Linked Result message.

The Set Result uses the following data structure:

```

typedef struct {
    ManagedObjectId      managed_object;
    AttributeList        attributeList;
} SetResult;

```

- **managed_object**
Identifies the object that responds to the Set command.
- **attributeList**
Contains all modified attributes.

Command Structure Summary

The following diagram shows how the different generic Protocol Command command structures are built from the different data type definitions that were introduced in this section.

SPpdu						
ROapdus						
ROIvapdu			RORSapdu ROLSapdu			ROERapdu
Event Report Argument	Action Argument	Get Argument Set Argument	Event Report Result	Action Result	Get Result Set Result	Error Data
Event Data	Action Data		Event Result Data	Action Result Data		

From this generic message structure the specific Protocol Command messages introduced in “Protocol Dialog” on page 25 are derived by:

- Defining identifier codes for the supported specific Event Report and Action types. These identifier codes are the values of the *event_type* and *action_type* fields.
- Defining the specific Event Data and Action Data data types for these Event Report and Action types.

Protocol Commands

This section describes the actual commands as constructed from the building blocks. Consult the “Command Structure Summary” on page 51 as a reference.

Notation

The Protocol Commands are constructed from the data types previously defined. A generic protocol machine must parse the individual elements of a command message separately, so in this chapter a special notation is used to define how the command messages are constructed (rather than defining composite C data type definitions).

Example:

```
MDSCreateEventReport ::=
  <SPpdu>
  <ROapdus (ro_type := ROIV_APDU)>
  <ROIvapdu (command_type := CMD_CONFIRMED_EVENT_REPORT) >
  <EventReportArgument (event_type := NOM_NOTI_MDS_CREAT)>
  <MDSCreateInfo>
```

This notation means that an MDS Create Event Report Command message is constructed from the individual data types listed in the < > brackets, which are C data types. Some elements of these data types have specific values. E.g. the *ro_type* field in the *ROapdus* data type has the value *ROIV_APDU*.

Additional data structures for appended event specific or method specific data are defined in the usual C type definition notation.

Most of the elements of the command messages contain length fields. You must take care to correctly set and parse these fields so that the message can be correctly parsed.

Device Discovery Messages

The Device Discovery messages lets the client locate new IntelliVue monitor devices in the network without prior knowledge of their IP address. The IntelliVue monitor only *sends* a Device Discovery on the LAN interface. This message is not available on the MIB/RS232 interface.

CONNECT INDICATION EVENT

The Connect Indication Event message is a sub-net-wide broadcast message in the normal Event Report format. It is sent to the port 24005.

The IntelliVue monitor resends the Connect Indication message as long as no logical connection to a central station has been established. The connection of a Data Export Computer Client does not stop the transmission of Connect Indication messages.

The IntelliVue monitor uses the retransmit strategy described in RFC 951. The initial resend period is 4 seconds, and this is doubled with each resend. The maximum resend period is approximately 64 seconds. The actual resend period contains a random component to avoid network congestion, e.g., after a power failure.

The UDP checksum in the Connect Indication message may be set to 0, indicating that no checksum has been calculated.

The Connect Indication message has the following structure:

```
ConnectIndication ::=
  <Nomenclature>
  <ROapdus (ro_type := ROIVapdu)>
  <ROIVapdu (command_type := CMD_EVENT_REPORT)>
  <EventReportArgument
    (managed_object := {NOM_MOC_MDS_COMPOS_SINGLE_BED, 0, 0},
     event_type := NOM_NOTI_MDS_CONNECT_INDIC)>
  <ConnectIndInfo>

typedef u_32 Nomenclature;
```

The nomenclature starts with two bytes 0x0, followed by one byte major and one byte minor version.

```
typedef AttributeList ConnectIndInfo;
```

See the section “Connect Indication Attributes” on page 107 for a list of attributes contained in the appended attribute list.

The Computer Client should parse the *ConnectIndInfo* to find out about the port for the Data Export protocol. The Computer Client must send requests to the port that is specified for the Data Export protocol.

The Computer Client application can run on any free local port, but must not change the port during the association (refer to “Definition of the Association Control Protocol” on page 65 for more information).

Connection Startup

After the logical connection has been established between the IntelliVue monitor and the Computer Client, the IntelliVue monitor sends the MDS Create Event message to announce version and status information.

MDS CREATE EVENT

The MDS Create Event describes the software and hardware configuration of the IntelliVue monitor. The Computer Client should parse this message to learn about the system configuration.

The MDS Create Event message has the following structure:

```
MDSCreateEventReport ::=
  <SPpdu>
  <ROapdus (ro_type := ROIV_APDU)>
  <ROIVapdu (command_type := CMD_CONFIRMED_EVENT_REPORT)>
  <EventReportArgument
    (managed_object := {NOM_MOC_VMS_MDS, 0, 0},
     event_type := NOM_NOTI_MDS_CREAT)>
  <MDSCreateInfo>
```

The MDS Create Information uses the following C type definition:

```
typedef struct {
    ManagedObjectId    managed_object;
    AttributeList      attribute_list;
} MdsCreateInfo;
```

- **managed_object**

Identifies the MDS object. Contents is the same as in the *managed_object* field in the Event Report structure.

- **attribute_list**

The attached *attribute_list* contains the IntelliVue monitor MDS attributes from the System Identification and from the System Application Attribute Group. See “Wave Objects” on page 82 for a list of all attributes.

Depending on the protocol and the protocol options which were negotiated when the association was established, the IntelliVue monitor may map its internal data representation to a representation which is supported by the negotiated protocol. Hence, the Connect Indication message may describe the system differently from the MDS Create Event message. In the case of differences, the MDS Create Event is the relevant message.

The MDS Create Event message contains both the "Date and Time" and the "Relative Time" attributes. The Computer Client can use this data to make a mapping from the relative time to the absolute time of the IntelliVue monitor. The Computer Client should regularly check if the mapping is still valid by sending a Single Poll Data Request for the MDS attributes (“SINGLE POLL DATA REQUEST” on page 55).

If the size of the Event Report (Event Report Result + Event Result Data) exceeds the size of a maximum message (MTU - Maximum Transmit Unit), multiple messages are sent. Each of these messages is sent as a single Event Report.

The Computer Client must confirm the MDS CREATE EVENT with a MDS CREATE EVENT RESULT message, otherwise the association will be aborted by the IntelliVue monitor. The MDS CREATE EVENT message is resent with a period of about 3 seconds. The association is aborted if the Event message has been sent 3 times without receiving a confirmation.

When the MDS Create Event message is resent, it has the same invoke ID as the original message.

MDS CREATE EVENT RESULT

As the MDS Create Event Report is a confirmed operation, the Computer Client must send a MDS Create Event Result message to confirm it.

The reply message has the following structure:

```
MDSCreateEventResult ::=
  <SPpdu>
  <ROapdus (ro_type := RORS_APDU)>
  <RORSapdu
    (invoke_id := mirrored from event report,
     command_type := CMD_CONFIRMED_EVENT_REPORT)>
  <EventReportResult
    (managed_object := mirrored from event report,
     event_type := NOM_NOTI_MDS_CREAT)
    length := 0 >
```

As the MDS Create Event Result message does not contain any appended additional information, the length of the appended information is set to 0.

The result message must have the same *invoke_id* as the event message.

Specific Data Access Commands

The following protocol commands are used to access the different types of data in the IntelliVue monitor.

SINGLE POLL DATA REQUEST

This message can be sent as soon as the logical connection is established and the MDS Create Event/Reply message sequence is finished. The message calls a method that returns IntelliVue monitor device data in a single response message.

The message has the following structure:

```
MDSPollAction ::=
  <SPpdu>
  <ROapdus (ro_type := ROIV_APDU)>
  <ROIVapdu (command_type := CMD_CONFIRMED_ACTION)>
  <ActionArgument
    (managed_object := {NOM_MOC_VMS_MDS, 0, 0},
     action_type := NOM_ACT_POLL_MDIB_DATA)>
  <PollMdibDataReq>
```

The *managed_object* must be the same as the *managed_object* in the MDS Create Event message. This is the top level object which actually implements the Data Export protocol.

The appended *PollMdibDataRequest* has the following data type:

```
typedef struct{
  u_16          poll_number;
  TYPE          polled_obj_type;
  OIDType       polled_attr_grp;
} PollMdibDataReq;
```

- **poll_number**
This field will be sent back in the response message. It is recommended to use this field as a counter.
- **polled_obj_type**
Defines which objects (Numerics or Alarms or MDS or Patient Demographics) is polled.

The following is a list of supported objects and their corresponding TYPE values:

NUMERICS:	partition:	0x0001
	code:	NOM_MOC_VMO_METRIC_NU
WAVES:	partition:	0x0001
	code:	NOM_MOC_VMO_METRIC_SA_RT
ALERTS:	partition:	0x0001
	code:	NOM_MOC_VMO_AL_MON
Pat.Demog:	partition:	0x0001
	code:	NOM_MOC_PT_DEMOG
MDS:	partition:	0x0001
	code:	NOM_MOC_VMS_MDS

The codes are taken from the Object Oriented Elements partition of the nomenclature (see “Object Classes” on page 111).

- `polled_attr_grp`

Defines which set of attributes is polled. For more information on the supported attribute groups and their contents, please refer to the section “Attribute Data Types and Constants Used” on page 75.

The IntelliVue monitor specifies limits on the maximum frequency for incoming SINGLE POLL DATA REQUEST messages. If the Computer Client sends messages with a frequency above the limit, some of the messages will be ignored (no response is sent). Separate limits are calculated for each object.

The IntelliVue monitor will process a maximum of one POLL DATA REQUEST messages for each object type per second. An additional POLL DATA REQUEST for Numeric Observed Values is allowed.

SINGLE POLL DATA RESULT

This message is sent by the IntelliVue monitor in response to the Single Poll Data Request.

The message has the following structure:

```
MDSPollActionResult ::=
  <SPpdu>
  <ROapdus (ro_type := RORS_APDU)>
  <RORSapdu (invoke_id := "mirrored from request message"
    command_type := CMD_CONFIRMED_ACTION)>
  <ActionResult
    (managed_object := {NOM_MOC_VMS_MDS, 0, 0},
    action_type := NOM_ACT_POLL_MDIB_DATA)>
  <PollMdibDataReply>
```

The appended *PollMdibDataReply* is constructed from the following data types:

- The *PollMdibDataReply* structure is the top level data structure returned in the Single Poll Data Result message. It contains the following fields:

```
typedef struct {
  u_16                poll_number;
  RelativeTime        rel_time_stamp;
  AbsoluteTime        abs_time_stamp;
  TYPE                polled_obj_type;
  OIDType             polled_attr_grp;
  PollInfoList        poll_info_list;
} PollMdibDataReply;
```

- `poll_number`
The poll number field contains the value of the same field in the Poll Request message.
- `rel_time_stamp`
The Relative Time Stamp is a high resolution time stamp that represents the system time when the event message is sent by the IntelliVue monitor.
For Numerics, the Relative Time Stamp denotes the time, when the Numeric measurement was generated. It may contain 0 if no measurement has been made yet.
- `abs_time_stamp`
The IntelliVue monitor does not support Absolute Time Stamps in the Poll Data Result. All fields contain 0xff. If the Computer Client needs Absolute Time Stamps, it should use the corresponding MDS attributes ("Relative Time" and "Date and Time" to map the `rel_time_stamp` to an *abs_time_stamp*.
- `polled_obj_type`
Defines for which objects (Numerics or Alarms or MDS or Patient Demographics) data is returned in the Poll Result message.
- `polled_attr_grp`
Defines which set of attributes is returned in the Poll Result message.
- `poll_info_list`
This structure contains the attribute values of the objects included in the poll.
The Poll Info List is an array structure where each *SingleContextPoll* element contains the poll result data of one naming context.

```
typedef struct {
    u_16          count;
    u_16          length;
    SingleContextPoll value[1];
} PollInfoList;
```

- `count`
Number of Single Context Poll structures that are appended.
- `length`
Size in bytes of the appended Single Context Poll structures.
- `value`
This field is a placeholder field only. It represents the specified number of appended Single Context Poll structures.

The Single Context Poll structure contains polled data of all object instances within one unique naming context (IntelliVue monitor supports multiple naming contexts). It contains the following fields:

```
typedef struct {
    MdsContext          context_id;
    struct {
        u_16            count;
        u_16            length;
        ObservationPoll  value[1];
    } poll_info;
} SingleContextPoll;
```

- `context_id`
The *context_id* field is used when the sourcing device represents multiple physical devices, so that the Handle attribute would not allow a unique identification of the object instance.
- `poll_info.count`
This field contains number of appended Observation Poll structures.
- `poll_info.length`
This field contains the length in bytes of the appended list of Observation Poll structures.
- `poll_info.value`
This field is a placeholder field only. It represents the specified number of appended Observation Poll structures.

The *ObservationPoll* represents the polled data of one object instance. It contains the following fields:

```
typedef struct {
    Handle          obj_handle;
    AttributeList    attributes;
} ObservationPoll;
```

- `obj_handle`
The handle identifies the object instance. It is used to identify the object in different Poll Reply Messages.
- `attributes`
The attributes field is a list structured field containing the values of the polled object attributes. For a list of supported object attributes, see the chapter on “Attribute Data Types and Constants Used” on page 75.

If the size of data returned for a Poll Result (Action Result + Action Result Data) exceeds the size of a maximum message (MTU - Maximum Transmit Unit), multiple messages are returned. These messages use the Remote Operation Linked Result mechanism (“Remote Operation Linked Result” on page 44). This means that in all result messages except the last result message the *ROLSapdu* is used instead of the *RORSapdu*.

When the Linked Result mechanism is used, the IntelliVue monitor may send the terminating Remote Operation Result message with an empty *PollInfoList* (count and length fields of the *PollInfoList* set to 0). It also may send Linked Result messages with one empty *SingleContextPoll* (count and length field of the *SingleContextPoll* set to 0).

EXTENDED POLL DATA REQUEST

The Extended Poll Data Request allows the following extensions of the Single Poll Data Request:

- Access 12 second, 1 minute and 5 minute averaged Numerics.
- Access wave data
- Request periodic Poll Replies without sending a Poll Request every time.
- Request that only a limited number of objects is encoded within a Poll Result

The Extended Poll Data Request message is only allowed, if the Poll Profile Extensions optional package has been negotiated during the association phase. For more information on the negotiation of optional packages see the sections “Association Request Message” on page 67 and “Association Response Message” on page 73.

The message has the following structure:

```
MDSPollAction ::=
  <SPpdu>
  <ROapdu (ro_type := ROIV_APDU)>
  <ROIVapdu (command_type := CMD_CONFIRMED_ACTION)>
  <ActionArgument
    (managed_object := {NOM_MOC_VMS_MDS, 0, 0},
     action_type := NOM_ACT_POLL_MDIB_DATA_EXT)>
  <PollMdibDataReqExt>
```

The appended *PollMdibDataRequestExt* has the following data type:

```
typedef struct{
  u_16          poll_number;
  TYPE          polled_obj_type;
  OIDType       polled_attr_grp;
  AttributeList poll_ext_attr;
} PollMdibDataReqExt;
```

- *poll_number*
This field will be sent back in the response message. It is recommended to use this field as a counter. See also the section “EXTENDED POLL DATA RESULT” on page 61 for more information about the handling of *poll_number*.
- *polled_obj_type*
Defines for which objects data is returned in the Poll Result message. The Extended Poll Data Request message only allows the polling of Numerics, Waves and the Alert Monitor.
- *polled_attr_grp*
Defines which set of attributes is returned in the Poll Result message.
- *poll_ext_attr*
The appended *AttributeList* allows to define additional options.

Accessing 12 second, 1 minute and 5 minute averaged Numerics

Within the Poll Profile Extensions optional package, the Computer Client and the IntelliVue monitor have negotiated which data source (real-time or averaged) is used to obtain the Numeric data (refer to the chapter “Definition of the Association Control Protocol” on page 65 for more information on how to negotiate optional packages). Currently, the IntelliVue monitor allows the specification of one data source for Numeric data.

The IntelliVue monitor responds to an Extended Poll Data Request message with an Extended Poll Data Result message, which contains the Numeric data from the source specified in the Poll Profile Extensions optional package.

The normal Poll Data Request message always returns data from real-time measurements. If another data source has been negotiated in the Poll Profile Extensions optional package, the Poll Data Request message will fail, if no data from real-time measurements is available.

The *poll_ext_attr AttributeList* in the Extended Poll Data Request message allows to specify additional options. Currently, the following attributes are supported:

Attribute: Time Periodic Data Poll

The Time Periodic Data Poll attribute allows to request periodic Poll Replies for a given time.

```
Attribute ID:      NOM_ATTR_TIME_PD_POLL
Attribute Type:    PollDataReqPeriod
Attribute Groups:  -
Availability:      Optional
```

The PollDataReqPeriod is defined as follows:

```
typedef struct {
    RelativeTime    active_period;
} PollDataReqPeriod;
```

The *active_period* specifies the time for which the IntelliVue monitor will send periodic Poll Replies.

The *AttributeList* Structure may contain additional attributes, e.g. in future releases.

If the Computer Client adds the Time Periodic Data Poll attribute to the Extended Poll Data Request message, the IntelliVue monitor sends periodic Extended Poll Data Result messages for the time specified in the attribute.

Data Source	Result Period
real-time waves	256ms
real-time measurements	1s
12 second averaged data	6s
1 minute averaged data	30s
5 minute averaged data	150s
alert data	1s

When the IntelliVue monitor receives an Extended Poll Data Request message, the first result message is sent immediately as a confirmation. It has the sequence number zero (see below). This allows the Computer Client to detect that its request was successful. The following messages are sent with the period specified in the table above.

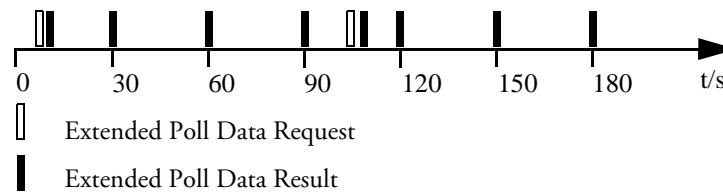
Figure 2 Period of Extended Poll Replies

The Computer Client should send a new Extended Poll Request before the time specified in the Time Periodic Data Poll attribute has expired. Each new Extended Poll Request is confirmed with an immediate Extended Poll Result message. However, the basic period of the replies is continued as illustrated in Figure 2 for 1 minute averaged data.

If the Computer Client uses the Extended Poll Request to access Realtime Numerics, it may happen that the IntelliVue monitor needs more than 1 second to encode all the data for the numerics (e.g. systems with a huge number of measurement modules). In this case the Poll Results will be sent at the highest possible frequency.

Limiting the Number of Objects in the Poll Result

In some cases, a Computer Client may want to limit the number of objects which are contained in a Poll Result. If the IntelliVue monitor is connected to a large number of measurement modules, a Poll Request for numerics will result in a large amount of data being sent from the IntelliVue monitor to the Computer Client.

Attribute: Number of Prioritized Objects

The attribute Number of Prioritized Objects specifies the maximum number of objects which will be encoded in the Poll Result.

```
Attribute ID:      NOM_ATTR_POLL_OBJ_PRIO_NUM
Attribute Type:    u_16
Attribute Groups:  -
Availability:      Optional
```

Based on an internal priority table, the IntelliVue monitor determines which objects will be added to the Poll Result. The priority table is constructed in the background, if the system configuration changes, it may take up to two minutes until the table has been updated. During this transition phase, the Poll Results sent by the monitor may contain less than the requested number of objects.

EXTENDED POLL DATA RESULT

When the IntelliVue monitor receives an Extended Poll Data Request message, it responds with a single or periodic Extended Poll Data Result messages.

The message has the following structure:

```

MDSPollActionResultExt ::=
  <SPpdu>
  <ROapdus (ro_type := RORS_APDU)>
  <RORSapdu (invoke_id := "mirrored from request message"
             command_type := CMD_CONFIRMED_ACTION)>
  <ActionResult
    (managed_object := {NOM_MOC_VMS_MDS, 0, 0},
     action_type := NOM_ACT_POLL_MDIB_DATA_EXT)>
  <PollMdibDataReplyExt>

```

The *PollMdibDataReplyExt* is defined as follows:

```

typedef struct PollMdibDataReplyExt {
  u_16          poll_number;
  u_16          sequence_no;
  RelativeTime  rel_time_stamp;
  AbsoluteTime  abs_time_stamp;
  TYPE          polled_obj_type;
  OIDType       polled_attr_grp;
  PollInfoList  poll_info_list;
} PollMdibDataReplyExt;

```

The *PollMdibDataReplyExt* structure is the top level data structure returned in the Extended Poll Data Result message. The appended data has the same structure as for the Single Poll Data Result.

The *PollMdibDataReplyExt* structure contains the following fields:

- **poll_number**
The poll number field contains the value of the same field in the Extended Poll Request message.
- **sequence_number**
The *sequence_number* is set to 0 when a new Extended Poll Data Request message is received. The IntelliVue monitor increases it with each periodic result message. This field allows the Computer Client to verify the sequence of the received result messages.
- **rel_time_stamp**
The Relative Time Stamp is a high resolution time stamp that represents the system time when the event message is sent by the IntelliVue monitor.
For Numerics, the Relative Time Stamp denotes the time when the Numeric measurement was generated. It may contain 0 if no measurement has been made yet.
For Waves, the Relative Time Stamp denotes the beginning of the 256ms result period for real-time waves.
- **abs_time_stamp**
The IntelliVue monitor does not support Absolute Time Stamps in the Poll Data Result. All fields contain 0xff. If the Computer Client needs Absolute Time Stamps, it should use the corresponding MDS attributes ("Relative Time" and "Date and Time" to map the *rel_time_stamp* to an *abs_time_stamp*).
- **polled_obj_type**
Defines for which objects (Numerics or Alarms or MDS or Patient Demographic) data is returned in the Poll Result message.
- **polled_attr_grp**
For more information on the supported attribute groups and their contents, please refer to the section "Attribute Data Types and Constants Used" on page 75.

- poll_info_list

This structure contains the attribute values of the objects included in the poll.

Keep Alive Message

The IntelliVue monitor closes an association if it does not receive any protocol commands within a specified time (see “Definition of the Association Control Protocol” on page 65 to learn how the limit for a timeout is negotiated). If the Computer Client sends messages with a very low frequency (e.g. when using the extended poll mechanism) it must send a keep alive message to prevent the IntelliVue monitor from closing the association.

It is suggested that the Computer Client sends a Poll Data Request message for this purpose. This has the advantage that the message is confirmed and the Computer Client can detect a possible loss of the message. The Computer Client should chose a Poll Request which results in as little processing overhead as possible.

A suitable keep alive message would be a Poll Request for the Alert Monitor object, requesting the VMO Static Context Attribute group. The associated Poll Result sent by the IntelliVue monitor is a short message.

Specify Objects in the Poll Result

The Get and Set operations can be used to specify wave objects or numeric objects to be reported within the Poll Results.

There is a default priority list which depends on an internal priority table and the current system configuration. For wave objects and numeric objects, the default list can be replaced by a user defined priority list.

Due to the high amount of data it is always recommended to specify the required wave objects before requesting wave data.

NOTE Software versions < E.0 may have limited support of this command.

GET PRIORITY LIST REQUEST

The message has the following structure:

```
MDSGetPriorityList ::=
  <SPpdu>
  <ROapdus (ro_type := ROIV_APDU)>
  <ROIVapdu (command_type := CMD_GET)>
  <GetArgument
    (managed_object := {NOM_MOC_VMS_MDS, 0, 0})>
```

The Get argument's *AttributeIdList* specifies the attribute identifiers:

- NOM_ATTR_POLL_RTSA_PRIO_LIST
Wave object priority list.
- NOM_ATTR_POLL_NU_PRIO_LIST
Numeric object priority list.

GET PRIORITY LIST RESULT

This message is sent in response to the Get Priority List Request.

The message has the following structure:

```
MDSGetPriorityListResult ::=
  <SPpdu>
  <ROapdus (ro_type := RORS_APDU)>
  <RORSapdu
    (invoke_id := "mirrored from request message",
     command_type := CMD_GET)>
  <GetResult
    (managed_object := {NOM_MOC_VMS_MDS, 0, 0})>
```

The Get result's *AttributeList* contains the requested attribute identifiers and values. The *TextIdList* structure is used to define the wave object priority list:

```
typedef struct {
    u_16 count;
    u_16 length;
    TextId value[1];
} TextIdList;
```

The array of *TextIds* specifies the objects by their label, as returned in the dynamic context.

SET PRIORITY LIST REQUEST

The message has the following structure:

```
MDSSetPriorityList ::=
  <SPpdu>
  <ROapdus (ro_type := ROIV_APDU)>
  <ROIVapdu (command_type := CMD_CONFIRMED_SET)>
  <SetArgument
    (managed_object := {NOM_MOC_VMS_MDS, 0, 0})>
```

The Set argument's *ModificationList* specifies the modify operations, attribute identifiers, and new values (if needed).

For the REPLACE operation, a wave object priority list attribute with modified *TextIdList* structure is attached.

For the SET_TO_DEFAULT operation, there is an empty attribute (*length* is 0) attached.

The ADD_VALUES and REMOVE_VALUES operations are not supported.

SET PRIORITY LIST RESULT

This message is sent in response to the Set Priority List Request.

The message has the following structure:

```
MDSSetPriorityListResult ::=
  <SPpdu>
  <ROapdus (ro_type := RORS_APDU)>
  <RORSapdu
    (invoke_id := "mirrored from request message",
     command_type := CMD_CONFIRMED_SET)>
  <SetResult
    (managed_object := {NOM_MOC_VMS_MDS, 0, 0})>
```

The Set result returns the modified *AttributeList*, as defined above.

Definition of the Association Control Protocol

Protocol Command Structure

The Protocol messages to establish the logical connection (association) between the IntelliVue monitor and a Computer Client follow the definitions of the ACSE Standard (ISO/IEC 8649 and ISO/IEC 8650), with some proprietary extensions.

All Association Control Commands share a common structure as shown here:

Session Header
Session Data
Presentation Header
User Data
Presentation Trailer

Figure 3 Protocol Commands for Association Control

For some messages, the Session Data and the User Data block may be empty.

A Computer Client can use the pre-defined building blocks for the Session Data, Presentation Header, and Presentation Trailer listed in the appendix to conveniently build valid messages (“Association Control Protocol Examples” on page 298 for a list of building blocks). Only the User Data block of the Association Request must be filled with Computer Client-specific data.

Protocol Commands

Protocol Command messages as defined in this section are the data structures that are transported within the transport layer messages.

The following commands are used to manage a logical connection between a Computer Client and a IntelliVue monitor:

- Association Request Message
- Association Response Message
- Refuse Message
- Release Request Message

- Release Response Message
- Abort Message

The Association Request message is sent from the Computer Client to the IntelliVue monitor when it wants to establish a new association. The *AssocReqUserData* contains information about the requested protocol and protocol options.

```
AssociationRequestMessage ::=
    <AssocReqSessionHeader>
    <AssocReqSessionData>
    <AssocReqPresentationHeader>
    <AssocReqUserData>
    <AssocReqPresentationTrailer>
```

The Association Response message is sent by the IntelliVue monitor if an Association Request message was parsed successfully and the association is accepted.

```
AssociationResponseMessage ::=
    <AssocRespSessionHeader>
    <AssocRespSessionData>
    <AssocRespPresentationHeader>
    <AssocRespUserData>
    <AssocRespPresentationTrailer>
```

If the Association Request message is corrupt, or if the association cannot be accepted (e.g. there is already another association), the IntelliVue monitor sends a Refuse message.

```
RefuseMessage ::=
    <RefuseSessionHeader>
    <RefuseSessionData>
    <RefusePresentationData>
    <RefuseUserData>
    <RefusePresentationTrailer>
```

When the Computer Client wants to terminate an association, it can send a Release Request message.

```
ReleaseRequestMessage ::=
    <ReleaseReqSessionHeader>
    <ReleaseReqSessionData>
    <ReleaseReqPresentationHeader>
    <ReleaseReqUserData>
    <ReleaseReqPresentationTrailer>
```

When the IntelliVue monitor receives a Release Request message, it sends a Release Response message as confirmation. The Release Response message indicates that the association has been terminated.

```
ReleaseRespMessage ::=
    <ReleaseRespSessionHeader>
    <ReleaseRespSessionData>
    <ReleaseRespPresentationHeader>
    <ReleaseRespUserData>
    <ReleaseRespPresentationTrailer>
```

The Abort message terminates an association without further confirmation. For example, the IntelliVue monitor sends an Abort message if an association is timed out (no communication from the Computer Client).

```
AbortMessage ::=
    <AbortSessionHeader>
    <AbortSessionData>
    <AbortPresentationHeader>
    <AbortUserData>
    <AbortPresentationTrailer>
```

Session Headers

The Session Headers can be used to identify the protocol commands. Each Session Header type maps to one protocol command.

The Session Header occupies the first bytes of the message. It is defined as follows:

```
typedef struct {
    u_8      type;
    #   define CN_SPDU_SI          0x0D
    #   define AC_SPDU_SI          0x0E
    #   define RF_SPDU_SI          0x0C
    #   define FN_SPDU_SI          0x09
    #   define DN_SPDU_SI          0x0A
    #   define AB_SPDU_SI          0x19
    LI      length;
} SessionHeader;
```

The *type* has the following meaning:

CN_SPDU_SI: A Session Connect header. The message contains an Association Request.

AC_SPDU_SI: A Session Accept header. The message contains an Association Response, indicating that the association has been established.

RF_SPDU_SI: A Session Refuse header. An association could not be established.

FN_SPDU_SI: A Session Finish header. The message contains a Release Request, indicating that the association should be terminated.

DN_SPDU_SI: A Session Disconnect header. The message contains a Release Response, indicating that the association has been terminated.

AB_SPDU_SI: A Session Abort header. The message contains an Abort message, indicating the immediate termination of the association.

If the first byte is 0xE1, the message is a Data Export Protocol command message (see “Definition of the Data Export Protocol” on page 35).

The *LI* field contains the length of the appended data (including all presentation data). The length encoding uses the following rules:

- If the length is smaller or equal 254 bytes, *LI* is one byte containing the actual length.
- If the length is greater than 254 bytes, *LI* is three bytes, the first being 0xff, the following two bytes containing the actual length.

Examples:

L = 15 is encoded as 0x0f

L = 256 is encoded as {0xff,0x01,0x00}

Message Encoding

The following section describes how a Computer Client can use the building blocks in the section “Association Control Protocol Examples” on page 298 to format correct Association Control messages.

Association Request Message

For the Association Request message, only the Session Header and the User Data must be filled out individually, as they contain variable data.

When using the building blocks, the presentation context ID for the Data Export Protocol is set to 2. This ID is sent in the SPpdu of all Data Export Protocol Commands.

The Session Header of the Association Request Message is defined as follows:

```
AssocReqSessionHeader ::=
    <SessionHead (type := CN_SPDU_SI)>
```

The length field in the Session Header must be set to the total length of the all appended data (including the presentation trailer).

Also the length field of the Presentation Header must be set to the total length of the appended message after this field. The field starts at the 2nd byte of the Presentation Header. It has the same format as the length field in the Session Header.

The User Data contains a specification of the requested protocol and protocol options. It is defined as follows:

```
AssocReqUserData ::=
    <ASNLength>
    <MDSEUserInfoStd>
```

The *ASNLength* contains the length of the *MDSEUserInfoStd*. It uses the following encoding rules:

- if the length is less or equal to 127, *ASNLength* is one byte, containing the actual length.
- if the length is greater than 127, *ASNLength* is several bytes long. The most significant bit (bit 0) of the first byte is set to 1, the bits 1 to 7 indicate the number of bytes which are appended to encode the actual length.

Examples:

L = 15 is encoded as 0x0f

L = 256 is encoded as {0x82,0x01,0x00}

The *MDSEUserInfoStd* is defined as follows:

```
typedef struct MDSEUserInfoStd {
    ProtocolVersion    protocol_version;
    NomenclatureVersion nomenclature_version;
    FunctionalUnits    functional_units;
    SystemType         system_type;
    StartupMode        startup_mode;
    AttributeList      option_list;
    AttributeList      supported_aprofiles;
} MDSEUserInfoStd;
```

The Computer Client must fill out the *MDSEUserInfoStd* data structure. It specifies the protocol versions and options the Computer Client supports. The IntelliVue monitor parses the *MDSEUserInfoStd* and constructs an Association Response message, which also contains a *MDSEUserInfoStd* data structure. The Association Response specifies which protocol versions and options will be used for the session.

The *ProtocolVersion* is a bit field containing the supported versions of the Data Export protocol. The Computer Client must set the bits for each version it supports. The IntelliVue monitor checks the supported versions and returns the bit for the highest commonly supported protocol version. If no matching version is found, the Association Request is refused.

```
typedef u_32    ProtocolVersion;
#define         MDDL_VERSION1  0x80000000
```

The *NomenclatureVersion* is a bit field containing the revision of the nomenclature which is used to name objects and their attributes. The Computer Client must set the bits for each version it supports. The IntelliVue monitor checks the supported versions and returns the bit for the highest commonly supported nomenclature version. If no matching version is found, the Association Request is refused.

```
typedef u_32 NomenclatureVersion;
#define         NOMEN_VERSION  0x40000000;
```

The *FunctionalUnits* is used to activate additional protocol functions. The Computer Client must set the bit for each functional unit it supports. The IntelliVue monitor checks the supported functional units and returns the bits for all commonly supported units (bitwise AND). No additional protocol functions have been defined yet.

```
typedef u_32      FunctionalUnits;
```

The *SystemType* is a bit field indicating whether the device is a Computer Client or a server. The Computer Client must set the SYST_CLIENT bit and the IntelliVue monitor will return the SYST_SERVER bit. If the SYST_CLIENT bit is not set in the Association Request, the association is refused.

```
typedef u_32 SystemType;
#define      SYST_CLIENT      0x80000000
#define      SYST_SERVER      0x00800000
```

The *StartupMode* is used to indicate the startup mode of the Computer Client and the IntelliVue monitor respectively. The IntelliVue monitor sets the bit for the startup mode which was used for the last reboot.

```
typedef u_32      StartupMode;
#define      HOT_START      0x80000000
#define      WARM_START      0x40000000
#define      COLD_START      0x20000000
```

If the IntelliVue monitor performs a COLD_START, all device settings are reset to the factory defaults. The configurations of the measurements might have changed and the patient data is lost.

The startup mode WARM_START and HOT_START indicate that configuration was not reset during the last restart.

The *option_list* can be used to negotiate additional protocol options in the form of an *AttributeList*. Currently, no further options are supported.

The *option_list* has a variable length. The offset of the *supported_aprofiles* field depends on the length of the *option_list*.

The *supported_aprofiles AttributeList* is used to define the available application profiles. An application profile specifies a set of protocol commands that is supported by the system. The Computer Client must add an entry for each supported profile to this list. The IntelliVue monitor parses the *supported_aprofiles* and returns the first profile in the list that is supported. If none of the profiles is supported, the Association Request is refused. The IntelliVue monitor supports the following profile:

Attribute: Poll Profile Support

The Poll Profile Support attribute contains the specification of the polling profile supported by the system.

```
Attribute ID:      NOM_POLL_PROFILE_SUPPORT
Attribute Type:    PollProfileSupport
Attribute Groups:  -
```

The *PollProfileSupport* is defined as follows:

```
typedef struct PollProfileSupport {
    PollProfileRevision    poll_profile_revision;
    RelativeTime            min_poll_period;
    u_32                    max_mtu_rx;
    u_32                    max_mtu_tx;
    u_32                    max_bw_tx;
    PollProfileOptions      options;
    AttributeList            optional_packages;
} PollProfileSupport;
```

The *PollProfileRevision* is a bit field containing the supported versions of the Polling Profile. The Computer Client must set the bits for each version it supports. The IntelliVue monitor checks the supported versions and returns the bit for the highest commonly supported profile version. If no matching version is found, the profile is not supported.

```
typedef u_32 PollProfileRevision;
#define POLL_PROFILE_REV_0          0x80000000
```

The *min_poll_period* specifies the minimum period with which the Computer Client wants to poll. If the IntelliVue monitor supports the requested poll period, it will return the value, otherwise it will return the minimum poll period it supports. The Computer Client should not send poll requests with a higher period than the negotiated value. For more information on poll periods, refer to the section “SINGLE POLL DATA REQUEST” on page 55.

The *min_poll_period* is also used to specify association time-outs. If the IntelliVue monitor does not receive any messages from the Computer Client within a given time, it sends an Abort message and terminates the association. The time-out periods depend on the negotiated *min_poll_period*, they are listed in the table below.

<i>min_poll_period</i>	Association Time out
< 3.3s	10s
3.3s ... 43s	3* <i>min_poll_period</i>
> 43s	130s

The *max_mtu_rx* and *max_mtu_tx* fields contain the maximum size (MTU - Maximum Transport Unit) for protocol commands (the size of the protocol command is the size of the data appended after the Remote Operation Header).

The MTU negotiation uses the following procedure:

- The Computer Client determines the maximum size of a protocol command it can send and receive.
- The Computer Client sets *max_mtu_tx* to the maximum size it can transmit (i.e. the IntelliVue monitor should provide receive capabilities for messages of this size) and the *max_mtu_rx* to the maximum size it can receive (i.e. the IntelliVue monitor should not send larger commands).
- The IntelliVue monitor determines the maximum size of a protocol command it can send and receive.
- The IntelliVue monitor sets *max_mtu_tx* to the maximum size the Computer Client is allowed to transmit (this is the minimum of the *max_mtu_tx* the Computer Client requested and the message size the IntelliVue monitor *can receive*). The IntelliVue monitor sets *max_mtu_rx* to the maximum size the client must be able to receive (this is the minimum of the *max_mtu_rx* the Computer Client requested and the message size the IntelliVue monitor *can send*).

Example:

- The Computer Client can send 800 bytes and receive 500 bytes of user data in one message.
- The Computer Client sets *max_mtu_tx* to 800 and *max_mtu_rx* to 500.
- The IntelliVue monitor can send 700 bytes and receive 600 bytes in one message.
- The IntelliVue monitor sets *max_mtu_tx* to 600 bytes (the IntelliVue monitor cannot receive larger messages) and *max_mtu_rx* to 500 bytes (the Computer Client can not receive more than 500 bytes in a message).

The IntelliVue monitor requires that the Computer Client can receive protocol commands of at least 300 bytes. Otherwise the profile is not supported. Smaller command sizes would lead to a considerable communication overhead. The largest negotiable MTU is 1364 bytes for the LAN interface and 1000 Bytes for the MIB/RS232 interface. The resulting size of the data packets may be larger than the MTU, because the MTU covers only the size of the Command Header and the Command Specific Data.

It is recommended that the Computer Client uses a large MTU. This reduces processing overhead and in most cases avoids splitting of messages.

For wave data export, the Computer Client needs to be able to receive observed values with 256 ms of wave data in one message. The MTU should be at least 500 bytes (700 bytes with multiplexed context).

The *max_bw_tx* contains the estimated maximum transmit bandwidth which will be used. The IntelliVue monitor fills in the maximum transmit bandwidth it uses, the value 0xffffffff indicates that no estimation is possible (this is the default). The current software does not support bandwidth estimation.

The *PollProfileOptions* bit field is used to set additional profile options. The IntelliVue monitor sets the P_OPT_DYN_CREATE_OBJECTS and P_OPT_DYN_DELETE_OBJECTS bits to indicate that the number of internal objects (e.g. the number of Numerics) may change dynamically. The *PollProfileOptions* is defined as follows:

```
typedef u_32 PollProfileOptions;
#define P_OPT_DYN_CREATE_OBJECTS      0x40000000
#define P_OPT_DYN_DELETE_OBJECTS     0x20000000
```

The *optional_packages AttributeList* allows the definition of additional options supported in the profile. The Computer Client must add an entry for each optional package it requests. The IntelliVue monitor checks the packages and adds an entry for each package it supports in the Association Response.

An attribute constitutes an optional package. The Poll Profile Extension is an optional package available for use.

Attribute: Poll Profile Extensions

The Poll Profile Extensions attribute specifies some extensions for the standard polling profile. For more information on how to use these extensions refer to the section “EXTENDED POLL DATA REQUEST” on page 59.

```
Attribute ID:      NOM_ATTR_POLL_PROFILE_EXT
Attribute Type:    PollProfileExt
Attribute Groups:  -
```

The *PollProfileExt* is defined as follows:

```
typedef struct {
    PollProfileExtOptions options;
    AttributeList          ext_attr;
} PollProfileExt;

typedef u_32 PollProfileExtOptions;
#define POLL_EXT_PERIOD_NU_1SEC      0x80000000
#define POLL_EXT_PERIOD_NU_AVG_12SEC 0x40000000
#define POLL_EXT_PERIOD_NU_AVG_60SEC 0x20000000
#define POLL_EXT_PERIOD_NU_AVG_300SEC 0x10000000
#define POLL_EXT_PERIOD_RTSA         0x08000000
#define POLL_EXT_ENUM                 0x04000000
#define POLL_EXT_NU_PRIO_LIST         0x02000000
#define POLL_EXT_DYN_MODALITIES       0x01000000
```

The *PollProfileExtOptions* bit field defines available options for the Poll Profile Extensions package.

If the POLL_EXT_PERIOD_NU_1SEC bit is set, the Computer Client requests real-time measurements as source for Numeric data.

If the POLL_EXT_PERIOD_NU_AVG_12SEC bit is set, the Computer Client requests 12 second averaged data as source for Numeric data.

If the POLL_EXT_PERIOD_NU_AVG_60SEC bit is set, the Computer Client requests 1 minute averaged data as source for Numeric data.

If the POLL_EXT_PERIOD_NU_AVG_300SEC bit is set, the Computer Client requests 5 minute averaged data as source for Numeric data.

The Computer Client must set at least one of the bits for the numeric period, otherwise the optional package is ignored. Currently, the IntelliVue monitor supports only one source for an association. If more than one of the bits is set, the source with the smallest measurement period is selected. The IntelliVue monitor sets the corresponding bit in the Association Response message.

There may be only one active numeric source at a given time. If there is an active association on the LAN interface which has requested realtime numerics, it is not possible to establish another association on the MIB/RS232 interface which requests 1 minute averaged data. In this case, the association request would result in a refuse message.

If the POLL_EXT_PERIOD_RTSA bit is set, the computer client requests wave data. The patient monitor sets the corresponding bit in its response message to indicate wave data export capability.

The Computer Client must parse the Association Response message to find out whether the requested options have been accepted by the IntelliVue monitor.

If the POLL_EXT_ENUM bit is set, the computer client is allowed to request Enumeration objects.

If the POLL_NU_PRIO_LIST bit is set, the computer client is allowed to set the numeric priority list.

If the POLL_EXT_DYN_MODALITIES bit is set, the computer client gets all timestamps for metrics with dynamic modalities. They are not exported otherwise for compatibility reasons.

The *ext_attr AttributeList* is reserved for future extensions.

Release Request Message

The Release Request message does not contain variable data. It is sufficient for the Computer Client to use the building blocks listed in the section “Association Control Protocol Examples” on page 298.

Abort Message

The Abort message does not contain variable data. It is sufficient for the Computer Client to use the building blocks listed in the section “Association Control Protocol Examples” on page 298.

Message Parsing

In most cases, it is sufficient for the Computer Client to check the first byte of the association control message. The first byte defines the Session Layer header, which can be mapped to an Association Control command.

Association Response Message

The IntelliVue monitor sends the Association Response message if an association has been established successfully. The Computer Client must parse the User Data within this message to find out which protocol options have been negotiated.

The Computer Client should not assume that the same Association Request message will always lead to the same Association Response message. The internal state of the IntelliVue monitor might lead to different responses.

The Association Response message is identified by its Session Header:

```
AssocRespSessionHeader ::=
  <SessionHead (type := AC_SPDU_SI)>
```

When parsing the Association Response message, the Computer Client must find the beginning of the User Data. This can be done by identifying the following byte sequence within the message;

```
0xBE 0x80 0x28 0x80 0x81
```

or

```
0xBE 0x80 0x28 0x80 0x02 0x01 0x02 0x81
```

The User Data is defined as follows;

```
AssocRespUserData ::=
  <ASNLength>
  <MDSEUserInfoStd>
```

The last byte of the User Data must be followed by 16 bytes 0x00.

The *MDSEUserInfo* follows the same definitions as described above for the Association Request Message.

Refuse

The IntelliVue monitor sends a Refuse message if an Association Request message was not accepted, because it was formatted incorrectly or because the requested protocol and protocol options are not supported by the IntelliVue monitor.

A Refuse message is also sent, if the maximum number of concurrent associations has been reached. Currently, the IntelliVue monitor only supports one active association.

The Refuse messages is identified by its Session Header:

```
RefuseSessionHeader ::=
  <SessionHead (type := RF_SPDU_SI)>
```

Release Response

It is sufficient to check the Session Header to detect a Release Response message. The Session Header is defined as follows:

```
ReleaseRespSessionHeader ::=
  <SessionHead (type := DN_SPDU_SI)>
```

Attribute Data Types and Constants Used

The data types in this chapter are based on the data types introduced in the chapter “Definition of the Data Export Protocol” on page 35. Refer to this chapter for more information about the base data types.

All data types used in this guide assume that elements of structures are aligned on 2 byte boundaries. Many compilers use different alignment modes by default. Make sure that the compiler uses the right alignments when parsing and formatting protocol messages.

The Poll Reply messages may contain attributes which are not documented here. A Computer Client should ignore all unknown attributes.

With IntelliVue release G the nomenclature of some numeric and wave labels have been changed. The labels that previously resided in the namespace NOM_EMFC are now merged into the NOM_SCADA namespace and the new defined NOM_SETTING namespace.

At the end of this chapter is a mapping table to guide you through the transition.

Numeric Objects

Numeric Object Attributes

This section defines the attributes of the Numeric object, together with the attribute identifier codes and attribute data types.

Attribute: Handle

The Handle attribute identifies the Numeric object in the form of a numeric value. The Handle is unique within a device context (see Common Data Type - Global Handle). The actual value of the Handle attribute does not have a meaning. It is used for reference and relation purposes (e.g. Alert Monitor entries reference the Numeric object instance by means of the Handle).

Attribute ID:	NOM_ATTR_ID_HANDLE
Attribute Type:	Handle (see Definitions Shared by Protocols)
Attribute Groups:	VMO Static Context Group
Availability:	Mandatory

Attribute: Type

The Type attribute contains an identification of the object type.

```
Attribute ID:      NOM_ATTR_ID_TYPE
Attribute Type:    TYPE (see Definitions Shared by Protocols)
Attribute Groups:  VMO Static Context Group
Availability:      Mandatory
```

Attribute: Numeric Observed Value

The Numeric Observed Value attribute represents the (measured) value, along with state and identification data.

```
Attribute ID:      NOM_ATTR_NU_VAL_OBS
Attribute Type:    NuObsValue (see below)
Attribute Groups:  Metric Observed Value Group
Availability:      Conditional (either NuObsValue or
                   NuObsValueCmp must be present)
```

The *NuObsValue* data type is defined as follows:

```
typedef struct {
    OIDType      physio_id;
    MeasurementState state;
    OIDType      unit_code;
    FLOATType    value;
} NuObsValue;
```

The *physio_id* (physiological identifier) field contains a nomenclature code from the SCADA partition that identifies the represented value (typically a physiological measurement).

The *unit_code* field contains a nomenclature code from the dimension nomenclature partition. It identifies the units of measure.

The *value* field is a floating point number with the actual value. Before interpreting the numeric value, the *state* must be checked. Only if *state* indicates a valid measurement, should the *value* field be interpreted.

The *state* field is a bit field structure (multiple bits can be set simultaneously) defined as follows:

```
typedef u_16      MeasurementState;
#define INVALID           0x8000
#define QUESTIONABLE     0x4000
#define UNAVAILABLE      0x2000
#define CALIBRATION_ONGOING 0x1000
#define TEST_DATA        0x0800
#define DEMO_DATA         0x0400
#define VALIDATED_DATA    0x0080
#define EARLY_INDICATION  0x0040
#define MSMT_ONGOING      0x0020
#define MSMT_STATE_IN_ALARM 0x0002
#define MSMT_STATE_AL_INHIBITED 0x0001
```

The bits have the following meaning:

INVALID: The source detects a sufficient degradation to render the data meaningless.

QUESTIONABLE: A problem exists, but it is still appropriate to present the data. This occurs when (1) either the degradation in the data is marginal or (2) the source cannot make a definite judgement on the reliability of the data.

UNAVAILABLE: The signal does not permit derivation of the numeric in question. This could be a transient state (e.g. first breath detected after an apnea -> no rate available), or a continuous state (no etCO₂ detection possible on a flat CO₂ wave).

CALIBRATION_ONGOING: Parameter is currently being calibrated.

TEST_DATA: The signal is an automatically generated test signal only and is not a valid patient signal. If this bit is set, the value is not suitable for patient diagnosis.

DEMO_DATA: The IntelliVue monitor runs in demonstration mode, the signal is automatically generated and is not a valid patient signal. If this bit is set, the value is not suitable for patient diagnosis.

VALIDATED_DATA: The value has been manually validated.

EARLY_INDICATION: The value represents an early estimate of the actual signal (the Non-Invasive Blood Pressure measurement e.g. sets this bit as soon as it has derived a systolic value, even if mean and diastolic values are still missing).

MSMT_ONGOING: A new aperiodic measurement is currently ongoing.

MSMT_STATE_IN_ALARM: Indicates that the numeric has an active alarm condition

MSMT_STATE_AL_INHIBITED: Alarms are switched off for the numeric (crossed bell)

The measurement is valid if the first octet of the state is all 0.

Attribute: Compound Numeric Observed Value

The Compound Numeric Observed Value attribute represents multiple (measured) values modelled in one Numeric object, along with state and identification data.

The Compound Numeric Observed Value is e.g. used to represent Blood Pressure measurements. For these measurements, systolic, diastolic and mean values are represented by a single Numeric object.

Attribute ID:	NOM_ATTR_NU_CMPD_OBS_VAL
Attribute Type:	NuObsValCmp (see below)
Attribute Groups:	Metric Observed Value Group
Availability:	Conditional (either NuObsValue or NuObsValueCmp must be present)

The *NuObsValueCmp* data type is defined as follows:

```
typedef struct {
    u_16      count;
    u_16      length;
    NuObsValue value[1];
} NuObsValueCmp;
```

The count field defines the number of *NuObsValue* elements in the structure. Note that the count field is variable, the number of elements may change over time. For a Blood Pressure measurement e.g there can be 3 values (systolic, diastolic, mean) or a single value only (mean only).

The length field defines the size of the array of *NuObsValue* structures in bytes.

The value field is a place holder for parsing.

Attribute: Absolute Time Stamp

The Absolute Time Stamp attribute is used to define a time tag for the current Numeric value. In the IntelliVue monitor, the attribute is used for aperiodic measurements only.

```
Attribute ID:      NOM_ATTR_TIME_STAMP_ABS
Attribute Type:    AbsoluteTime(see Definitions Shared by Protocols)
Attribute Groups:  Metric Observed Value Group
Availability:      Optional
```

Attribute: Relative Time Stamp

The Relative Time Stamp attribute is used to define a high resolution time tag for the current Numeric value.

```
Attribute ID:      NOM_ATTR_TIME_STAMP_REL
Attribute Type:    RelativeTime(see Definitions Shared by Protocols)
Attribute Groups:  Metric Observed Value Group
Availability:      Optional
```

Attribute: Label

The Label attribute is a 32 bit wide ID which represents the Numeric label string. The Label is unique for all numerics in the system.

```
Attribute ID:      NOM_ATTR_ID_LABEL
Attribute Type:    TextId
                  (see Protocol Common Definitions)
Attribute Group:   VMO Dynamic Context Group
Availability:      Optional
```

Attribute: Label String

The Label String attribute is a unicode string which contains the label string for a Numeric.

```
Attribute ID:      NOM_ATTR_ID_LABEL_STRING
Attribute Type:    String
                  (see Protocol Common Definitions)
Attribute Group:   VMO Dynamic Context Group
Availability:      Optional
```

The Label String does not contain the asterisk prefix displayed by the monitor to indicate whether a numeric has been manually entered. The client has to check the MetricCategory field of the MetricSpec attribute and add this prefix to obtain the identical label string that is displayed by the monitor.

Attribute: Display Resolution

The Display Resolution attribute is present if the resolution of the numeric shown on the display must be different from the resolution communicated in the Numeric Observed Value attribute. E.g. a Temperature is displayed with a resolution of 1/10, but the Observed Value is sent with a precision of 1/100 to get the necessary accuracy for differential temperatures. The Display Resolution attribute describes the format in which the value of a numeric is displayed on the screen.

```
Attribute ID:      NOM_ATTR_DISP_RES
Attribute Type:    DispResolution
Attribute Group:   VMO Dynamic Context Group
Availability:      Optional
```

The *DispResolution* is defined as follows:

```
typedef struct
{
    u_8    pre_point;
    u_8    post_point;
} DispResolution;
```

The value of *pre_point* denotes the number of digits before the decimal point. The value of *post_point* denotes the number of digits after the decimal point.

Attribute: Color

The Color attribute describes the color in which a numeric is displayed on the screen.

Attribute ID:	NOM_ATTR_COLOR
Attribute Type:	SimpleColour
Attribute Group:	VMO Dynamic Context Group
Availability:	Optional

The *SimpleColour* is defined as follows:

```
typedef u_16 SimpleColour;
#define COL_BLACK      0
#define COL_RED        1
#define COL_GREEN      2
#define COL_YELLOW     3
#define COL_BLUE       4
#define COL_MAGENTA    5
#define COL_CYAN       6
#define COL_WHITE      7
#define COL_PINK       20
#define COL_ORANGE     35
#define COL_LIGHT_GREEN 50
#define COL_LIGHT_RED  65
```

Attribute: Metric Specification

The Metric Specification attribute describes static properties of a numeric.

Attribute ID:	NOM_ATTR_METRIC_SPECN
Attribute Type:	MetricSpec
Attribute Group:	VMO Static Context Group
Availability:	Mandatory

The *MetricSpec* is defined as follows:

```
typedef struct
{
    RelativeTime    update_period;
    MetricCategory  category;
    MetricAccess    access;
    MetricStructure structure;
    MetricRelevance relevance;
} MetricSpec;
```

The *update_period* is the minimum time between changes of the observed value.

The *MetricCategory* is defined as follows:

```
typedef u_16      MetricCategory;
#define MCAT_UNSPEC      0
#define AUTO_MEASUREMENT 1
#define MANUAL_MEASUREMENT 2
#define AUTO_SETTING     3
#define MANUAL_SETTING   4
#define AUTO_CALCULATION 5
#define MANUAL_CALCULATION 6
#define MULTI_DYNAMIC_CAPABILITIES 50
#define AUTO_ADJUST_PAT_TEMP 128
#define MANUAL_ADJUST_PAT_TEMP 129
#define AUTO_ALARM_LIMIT_SETTING 130
```

It allows to distinguish between measurements, calculations and settings. The values have the following meaning:

MCAT_UNSPEC: not specified

AUTO_MEASUREMENT: automatic measurement

MANUAL_MEASUREMENT: manual measurement

AUTO_SETTING: automatic setting

MANUAL_SETTING: manual setting

AUTO_CALCULATION: automatic calculation, e.g. differential temperature

MANUAL_CALCULATION: manual calculation

MULTI_DYNAMIC_CAPABILITIES: this measurement may change its category during operation or may be used in various modes.

AUTO_ADJUST_PAT_TEMP: measurement is automatically adjusted for patient temperature

MANUAL_ADJUST_PAT_TEMP: measurement manually adjusted for patient temperature

AUTO_ALARM_LIMIT_SETTING: this is not a measurement, but an alarm limit setting

The *MetricAccess* bit field provides info on how the metric value can be accessed and when a measurement is available.

```
typedef u_16      MetricAccess;
#define AVAIL_INTERMITTEND 0x8000
#define UPD_PERIODIC      0x4000
#define UPD_EPISODIC      0x2000
#define MSMT_NONCONTINUOUS 0x1000
```

The values have the following meaning:

AVAIL_INTERMITTEND: The intermitted availability bit is set, if the observed values not always available (e.g. only if a measurement is explicitly started).

UPD_PERIODIC: observed value is updated periodically

UPD_EPISODIC: observed value is updated episodically (exactly one update mode (UPD_) must be set

MSMT_NONCONTINUOUS: indicates that the measurement is non continuous (this is different from the update mode)

The *MetricStructure* describes if the object represents a single measurement or multiple related measurements (an invasive blood pressure could be compound when it represents

a pulsatile pressure like ABP and derives systolic, diastolic, mean values)

```
typedef struct MetricStructure {
    u_8          ms_struct;
    u_8          ms_comp_no;
} MetricStructure;
```

ms_struct describes the structure of the object, 0 means simple, 1 means compound object.

ms_comp_no contains the maximum number of components in the compound, it contains 0 for simple objects.

The MetricRelevance is a 16 bit wide field for internal use only.

```
typedef u_16 MetricRelevance;
```

Attribute MetricModality

The MetricModality attribute describes metric properties of a numeric that may depend on the usage of the measurement device to obtain a measurement.

Attribute ID:	NOM_ATTR_METRIC_MODALITY
Attribute Type:	MetricModality
Attribute Group:	Metric Observed Value Group
Availability:	Mandatory

The *MetricModality* is defined as follows:

```
typedef u_16 MetricModality;

#define METRIC_MODALITY_MANUAL      0x4000
#define METRIC_MODALITY_APERIODIC  0x2000
#define METRIC_MODALITY_VERIFIED   0x1000
```

Attribute Groups

The attributes of the Numeric object are arranged in the following attribute groups:

Attribute Group:	VMO Static Context Group
Group ID:	NOM_ATTR_GRP_VMO_STATIC
Description:	Static context of the object
Attributes:	Type, Handle, Metric Specification
Attribute Group:	VMO Dynamic Context Group
Group ID:	NOM_ATTR_GRP_VMO_DYN
Description:	Dynamic context of the object
Attributes:	Label, Label String, Color, Display Resolution
Attribute Group:	Metric Observed Value Group
Group ID:	NOM_ATTR_GRP_METRIC_VAL_OBS
Description:	Observed values of the object
Attributes:	Nu Observed Value, Compound Nu Observed Value, Absolute Time Stamp, Relative Time Stamp, MetricModality

Dynamic Context Changes

Internally, the IntelliVue monitor uses two different communication channels for attributes from the VMO Dynamic Context Group and the Metric Observed Value Group. This can lead to possible inconsistencies between these two attribute groups. Imagine that a Computer Client is polling all attribute groups. If the user changes the Label of a numeric (VMO Dynamic Context Group), the *physio_id* in the Nu Observed Value (Metric Observed Value Group) may be updated a short period later.

For real-time Numerics, this inconsistency is typically resolved after less than one second with the periodic update of the Observed Values. For averaged Numerics, the update of the Observed Values depends on the averaging period. It may be 12 seconds, 1 minute or 5 minutes.

Wave Objects

Wave Object Attributes

This section defines the attributes of the Wave object, together with the attribute identifier codes and attribute data types.

Attribute: Handle

The Handle attribute contains an identification of the wave object in the form of a numeric value. The actual value of the Handle attribute does not have a meaning. It is used for reference and relation purposes.

Attribute ID:	NOM_ATTR_ID_HANDLE
Attribute Type:	Handle (see Definitions Shared by Protocols)
Attribute Groups:	VMO Static Context Group
Availability:	Mandatory

Attribute: Type

The Type attribute contains an identification of the object type.

Attribute ID:	NOM_ATTR_ID_TYPE
Attribute Type:	TYPE (see Definitions Shared by Protocols)
Attribute Groups:	VMO Static Context Group
Availability:	Mandatory

Attribute: Metric Specification

The Metric Specification describes static properties of a metric object.

Attribute ID:	NOM_ATTR_METRIC_SPECN
Attribute Type:	MetricSpec
Attribute Groups:	VMO Static Context Group
Availability:	Mandatory

The *MetricSpec* is defined as follows:

```
typedef struct {
    RelativeTime      update_period;
    MetricCategory    category;
    MetricAccess      access;
    MetricStructure    structure;
    MetricRelevance    relevance;
} MetricSpec;
```

The *update_period* specifies the time between observed values.

MetricCategory, *MetricAccess*, *MetricStructure*, and *MetricRelevance* are already defined for the Numeric object.

Attribute: Sample Array Specification

The Sample Array Specification describes static properties of a wave object.

Attribute ID:	NOM_ATTR_SA_SPECN
---------------	-------------------

```

Attribute Type:      SaSpec
Attribute Groups:    VMO Static Context Group
Availability:        Mandatory

```

The *SaSpec* is defined as follows:

```

typedef struct {
    u_16          array_size;
    SampleType     sample_type;
    SaFlags        flags;
} SaSpec;

```

The *array_size* specifies the maximum number of samples in one observed value.

The *SampleType* is defined as follows:

```

typedef struct {
    u_8          sample_size;
    u_8          significant_bits;
} SampleType;

```

The *sample_size* specifies the number of bits used to encode one wave sample.

The number of *significant_bits* is less or equal *sample_size*. To get the actual sample value, non-significant bits must be masked if indicated in the flags value.

The *SaFlags* is defined as follows:

```

typedef u_16      SaFlags;
#define SMOOTH_CURVE          0x8000
#define DELAYED_CURVE         0x4000
#define STATIC_SCALE          0x2000
#define SA_EXT_VAL_RANGE      0x1000

```

The values have the following meaning:

SMOOTH_CURVE, DELAYED_CURVE : used for wave presentation

STATIC_SCALE: Scale and range specification does not change.

SA_EXT_VAL_RANGE: The non-significant bits in the sample value must be masked.

Attribute: Sample Array Fixed Value Specification

The Sample Array Fixed Value Specification defines a list of fixed sample values or bit masks that indicate specific conditions.

```

Attribute ID:      NOM_ATTR_SA_FIXED_VAL_SPECN
Attribute Type:    SaFixedValSpec16
Attribute Groups:  VMO Static Context Group
Availability:      Optional

```

The *SaFixedValSpec16* is a sequence of *SaFixedValSpecEntry16* elements:

```

typedef struct {
    u_16          count;
    u_16          length;
    SaFixedValSpecEntry16 value[1];
} SaFixedValSpec16;
typedef struct {
    SaFixedValId  sa_fixed_val_id;
    u_16          sa_fixed_val;
} SaFixedValSpecEntry16;

```

The *SaFixedValId* is defined as follows:

```

typedef u_16      SaFixedValId;
#define SA_FIX_UNSPEC          0
#define SA_FIX_INVALID_MASK   1
#define SA_FIX_PACER_MASK     2

```

```
#define SA_FIX_DEFIB_MARKER_MASK    3
#define SA_FIX_SATURATION           4
#define SA_FIX_QRS_MASK             5
```

The values have the following meaning:

SA_FIX_UNSPEC: Not specified.

SA_FIX_INVALID_MASK: Invalid sample mask.

SA_FIX_PACER_MASK: Pace pulse detected.

SA_FIX_DEFIB_MARKER_MASK: Defib marker in this sample.

SA_FIX_SATURATION: Indicates saturation condition in this sample.

(Note: despite the name, this is a mask as well.)

SA_FIX_QRS_MASK: Indicates QRS trigger around this sample.

The *sa_fixed_val* may be a value or a bit mask, as indicated in the *sa_fixed_val_id*.

Attribute: Sample Period

The Sample Period specifies the sample rate.

```
Attribute ID:      NOM_ATTR_TIME_PD_SAMP
Attribute Type:    RelativeTime (see Definitions Shared by Protocols)
Attribute Groups:  VMO Static Context Group
Availability:      Mandatory
```

Attribute: Label

The Label attribute contains a 32 bit wide ID which represents the wave label string. The Label is unique for all waves in the system.

```
Attribute ID:      NOM_ATTR_ID_LABEL
Attribute Type:    TextId (see Definitions Shared by Protocols)
Attribute Groups:  VMO Dynamic Context Group
Availability:      Optional
```

Attribute: Label String

The Label String is a unicode string which contains the label string for a wave.

```
Attribute ID:      NOM_ATTR_ID_LABEL_STRING
Attribute Type:    String (see Definitions Shared by Protocols)
Attribute Groups:  VMO Dynamic Context Group
Availability:      Optional
```

Attribute: Metric State

The Metric State attribute indicates metric on or off state.

```
Attribute ID:      NOM_ATTR_METRIC_STAT
Attribute Type:    MetricState
Attribute Groups:  VMO Dynamic Context Group
Availability:      Optional
```

The *MetricState* is a bit field defined as follows:

```
typedef u_16      MetricState;
#define METRIC_OFF      0x8000
```

Attribute: Unit Code

The Unit Code attribute contains a nomenclature code from the dimension partition. It identifies the units of measure.

Attribute ID:	NOM_ATTR_UNIT_CODE
Attribute Type:	OIDType (see Definitions Shared by Protocols)
Attribute Groups:	VMO Dynamic Context Group
Availability:	Optional

Attribute: Color

The Color attribute describes the color in which a wave is displayed on the screen.

Attribute ID:	NOM_ATTR_COLOR
Attribute Type:	SimpleColour
Attribute Groups:	VMO Dynamic Context Group
Availability:	Optional

The *SimpleColour* is already defined for the Numeric object.

Attribute: Measure Mode

The Measure Mode attribute defines specific measurement modes.

Attribute ID:	NOM_ATTR_MODE_MSMT
Attribute Type:	MeasureMode
Attribute Groups:	VMO Dynamic Context Group
Availability:	Optional

For wave objects, the following *MeasureMode* bits are defined:

```
typedef u_16      MeasureMode;
#define CO2_SIDESTREAM      0x0400
#define ECG_PACED           0x0200
#define ECG_NONPACED       0x0100
#define ECG_DIAG            0x0080
#define ECG_MONITOR        0x0040
#define ECG_FILTER          0x0020
#define ECG_MODE_EASI       0x0008
#define ECG_LEAD_PRIMARY    0x0004
```

The values have the following meaning:

CO2_SIDESTREAM: CO₂ sidestream.

ECG_PACED, ECG_NONPACED: Paced mode setting.

ECG_DIAG, ECG_MONITOR, ECG_FILTER: ECG filter setting.

ECG_MODE_EASI: EASI derived lead.

ECG_LEAD_PRIMARY: ECG primary lead.

Attribute: Metric Info Label

The Metric Info Label allows to specify an additional dynamic text (32 bit ID).

Attribute ID:	NOM_ATTR_METRIC_INFO_LABEL
Attribute Type:	TextId (see Definitions Shared by Protocols)
Attribute Groups:	VMO Dynamic Context Group
Availability:	Optional

Attribute: Metric Info Label String

The Metric Info Label String allows to specify an additional dynamic text (unicode string).

Attribute ID:	NOM_ATTR_METRIC_INFO_LABEL_STR
---------------	--------------------------------

Attribute Type:	String (see Definitions Shared by Protocols)
Attribute Groups:	VMO Dynamic Context Group
Availability:	Optional

Attribute: Scale and Range Specification

The Scale and Range Specification describes a relation between scaled values and absolute values and also defines the range of the measured values and samples.

Attribute ID:	NOM_ATTR_SCALE_SPECN_I16
Attribute Type:	ScaleRangeSpec16
Attribute Groups:	VMO Dynamic Context Group
Availability:	Mandatory

The *ScaleRangeSpec16* is defined as follows:

```
typedef struct {
    FLOATType    lower_absolute_value;
    FLOATType    upper_absolute_value;
    u_16         lower_scaled_value;
    u_16         upper_scaled_value;
} ScaleRangeSpec16;
```

The scaled values refer to the wave samples in the observed values.

If the wave does not represent any absolute value, the absolute value fields must be *NaN* (Not a Number).

Attribute: Sample Array Physiological Range

The Sample Array Physiological Range is used for display scaling.

Attribute ID:	NOM_ATTR_SA_RANGE_PHYS_I16
Attribute Type:	ScaledRange16
Attribute Groups:	VMO Dynamic Context Group
Availability:	Optional

The *ScaledRange16* is defined as follows:

```
typedef struct {
    u_16         lower_scaled_value;
    u_16         upper_scaled_value;
} ScaledRange16;
```

Attribute: Visual Grid

The Visual Grid attribute allows to define grid lines.

Attribute ID:	NOM_ATTR_GRID_VIS_I16
Attribute Type:	SaVisualGrid16
Attribute Groups:	VMO Dynamic Context Group
Availability:	Optional

The *SaVisualGrid16* is defined as follows:

```
typedef struct {
    u_16         count;
    u_16         length;
    SaGridEntry16 value[1];
} SaVisualGrid16;
typedef struct {
    FLOATType    absolute_value;
    u_16         scaled_value;
    u_16         level;
} SaGridEntry16;
```

Different *levels* define relative importance of grid lines. 0 is the first (most important) level.

Attribute: Sample Array Calibration Specification

The Sample Array Calibration Specification allows to define the presence of a calibration bar or calibration stair.

Attribute ID:	NOM_ATTR_SA_CALIB_I16
Attribute Type:	SaCalData16
Attribute Groups:	VMO Dynamic Context Group
Availability:	Optional

The *SaCalData16* is defined as follows:

```
typedef struct {
    FLOATType      lower_absolute_value;
    FLOATType      upper_absolute_value;
    u_16           lower_scaled_value;
    u_16           upper_scaled_value;
    u_16           increment;
    u_16           cal_type;
#define BAR        0
#define STAIR      1
} SaCalData16;
```

Attribute: Sample Array Observed Value

The Sample Array Observed Value attribute represents the wave samples, along with state and identification data.

Attribute ID:	NOM_ATTR_SA_VAL_OBS
Attribute Type:	SaObsValue
Attribute Groups:	Metric Observed Value Group
Availability:	Conditional (either SaObsValue or SaObsValueCmp is present)

The *SaObsValue* data type is defined as follows:

```
typedef struct {
    OIDType      physio_id;
    MeasurementState state;
    struct {
        u_16      length;
        u_8      value[1];
    } array;
} SaObsValue;
```

The *physio_id* (physiological identifier) field contains a nomenclature code from the SCADA partition that identifies the represented wave (typically a physiological measurement).

The *state* indicates measurement validity. Refer to the Numeric object for a definition of the bit field. The measurement is valid if the first octet of the *state* is all 0.

Attribute: Compound Sample Array Observed Value

The Compound Sample Array Observed Value attribute represents multiple waves modelled in one Wave object, along with state and identification data.

Compound Sample Array Observed Values are used to provide 250 samples/s ECG waves with common context.

Attribute ID:	NOM_ATTR_SA_CMPD_VAL_OBS
Attribute Type:	SaObsValueCmp
Attribute Groups:	Metric Observed Value Group

Availability: Conditional (either *SaObsValue* or *SaObsValueCmp* is present)

The *SaObsValueCmp* data type is defined as follows:

```
typedef struct {
    u_16          count;
    u_16          length;
    SaObsValue    value[1];
} SaObsValueCmp;
```

The *count* field defines the number of *SaObsValue* elements in the structure.

The *length* field defines the size of the array of *SaObsValue* structures in bytes.

The *SaObsValue* data type is defined above. The elements in a compound observed value can be identified by their *physio_id*.

Attributes Groups

The attributes of the Wave object are arranged in the following attribute groups:

Attribute Group:	VMO Static Context Group
Group ID:	NOM_ATTR_GRP_VMO_STATIC
Description:	Static context of the object
Attributes:	Handle, Type, Metric Specification, Sample Array Specification, Sample Array Fixed Value Specification, Sample Period
Attribute Group:	VMO Dynamic Context Group
Group ID:	NOM_ATTR_GRP_VMO_DYN
Description:	Dynamic context of the object
Attributes:	Label, Label String, Metric State, Unit Code, Color, Measure Mode, Metric Info Label, Metric Info Label String, Scale and Range Specification, Sample Array Physiological Range, Visual Grid, Sample Array Calibration Specification
Attribute Group:	Metric Observed Value Group
Group ID:	NOM_ATTR_GRP_METR_VAL_OBS
Description:	Observed values of the object
Attributes:	Sample Array Observed Value, Compound Sample Array Observed Value

Enumeration Objects

Enumeration Object Attributes

This section defines the attributes of the enumeration objects, together with the attribute identifier codes and attribute data types.

Enumeration Objects are not available in software revisions below E.0.

Attribute: Handle

The Handle attribute identifies the enumeration object in the form of a numeral value. The Handle is unique within a device context (see Common Data Type - Global Handle). The actual value of the Handle attribute does not have a meaning. It is used for reference and relation purposes.

Attribute ID:	NOM_ATTR_ID_HANDLE
Attribute Type:	Handle (see Definitions Shared by Protocols)
Attribute Groups:	-
Availability:	Mandatory

Attribute: Type

The Type attribute contains an identification of the object type.

Attribute ID:	NOM_ATTR_ID_TYPE
Attribute Type:	TYPE (see Definitions Shared by Protocols)
Attribute Groups:	VMO Static Context Group
Availability:	Mandatory

Attribute: Metric Specification

The Metric Specification describes static properties of a metric object.

Attribute ID:	NOM_ATTR_METRIC_SPECN
Attribute Type:	MetricSpec
Attribute Groups:	VMO Static Context Group
Availability:	Mandatory

The *MetricSpec* is defined as follows:

```
typedef struct {
    RelativeTime      update_period;
    MetricCategory    category;
    MetricAccess      access;
    MetricStructure    structure;
    MetricRelevance   relevance;
} MetricSpec;
```

The *update_period* specifies the time between observed values.

MetricCategory, *MetricAccess*, *MetricStructure*, and *MetricRelevance* are already defined for the Numeric object.

Attribute: Label

The Label attribute is a 32 bit wide ID which represents the enumeration label string.

Attribute ID:	NOM_ATTR_ID_LABEL
Attribute Type:	TextId (see Protocol Common Definitions)
Attribute Group:	VMO Dynamic Context Group
Availability:	Optional

Attribute: Label String

The Label String attribute is a unicode string which contains the label string for an enumeration.

Attribute ID:	NOM_ATTR_ID_LABEL_STRING
Attribute Type:	String (see Protocol Common Definitions)
Attribute Group:	VMO Dynamic Context Group
Availability:	Optional

Attribute: Color

The Color attribute describes the color in which an enumeration is displayed on the screen.

Attribute ID:	NOM_ATTR_COLOR
Attribute Type:	SimpleColour
Attribute Group:	VMO Dynamic Context Group
Availability:	Optional

Attribute: Enum-Observed-Value

The Enum-Observed-Value attribute describes the current state of the enumeration object.

```

Attribute ID:      NOM_ATTR_VAL_ENUM_OBS
Attribute Type:    EnumObsVal
Attribute Group:   VMO Observed Value Group
Availability:      Optional

typedef struct {
    OIDType      physio_id;
    MeasurementState state;
    EnumVal      value;
} EnumObsVal;

typedef struct {
    OIDType      obj_id;
    FLOATType    num_val;
    OIDType      unit_code;
} EnumObjIdVal;

typedef struct {
    u_16          choice;
    u_16          length;
    union {
        OIDType      enum_obj_id;
        EnumObjIdVal enum_obj_id_val;
    } u;
} EnumVal;

```

The field choice of structure EnumVal defines the valid structure of union u. Its values are defined as follows:

```

#define ENUM_OBJ_ID_CHOSEN      1
#define ENUM_OBJ_ID_VAL_CHOSEN 4

```

Attribute Groups

The attributes of the enumeration object are arranged in the following attribute groups:

```

Attribute Group:  VMO Static Context Group
Group ID:         NOM_ATTR_GRP_VMO_STATIC
Description:      Static context of the object
Attributes:       Type, Handle, Metric Specification

Attribute Group:  VMO Dynamic Context Group
Group ID:         NOM_ATTR_GRP_VMO_DYN
Description:      Dynamic context of the object
Attributes:       Label, Label String, Color

Attribute Group:  Metric Observed Value Group
Group ID:         NOM_ATTR_GRP_METRIC_VAL_OBS
Description:      Observed values of the object
Attributes:       Enumeration Observed Value,
                  Absolute Time Stamp, Relative Time Stamp

```

System Objects

System Objects Attributes

This section defines the attributes of the Medical Device System (MDS) object, together with the attribute identifier codes and attribute data types.

Attribute: Handle

The Handle attribute identifies the MDS object in the form of a numeral value. The Handle is unique within a device context (see Common Data Type - Global Handle). The actual value of the Handle attribute does not have a meaning. It is used for reference and relation purposes.

```
Attribute ID:      NOM_ATTR_ID_HANDLE
Attribute Type:    Handle (see Definitions Shared by Protocols)
Attribute Groups:  -
Availability:      Mandatory
```

Attribute: System Type

The System Type attribute contains an identification of the device type identified with the MDS object (e.g. monitor)

```
Attribute ID:      NOM_ATTR_SYS_TYPE
Attribute Type:    TYPE (see Definitions Shared by Protocols)
Attribute Groups:  System Identification Attribute Group
Availability:      Mandatory
```

For the MDS object, the OBJ nomenclature partition is used. The code value is a static identification.

Attribute: System Model

The System Model attribute contains a manufacturer ID and a manufacturer-specific model number for the device.

```
Attribute ID:      NOM_ATTR_ID_MODEL
Attribute Type:    SystemModel
Attribute Groups:  System Identification Attribute Group
Availability:      Mandatory
```

The *SystemModel* is defined as follows:

```
typedef struct {
    VariableLabel manufacturer;
    VariableLabel model_number;
} SystemModel;
```

The *manufacturer* field is of variable length, hence the offset of *model_number* depends on the length of *manufacturer*. Currently, the IntelliVue monitor uses 4 characters for the *manufacturer* and 6 characters for the *model_number* (including the terminating '\0').

Attribute: System ID

The Sytem ID attribute contains a unique identifier for the device.

```
Attribute ID:      NOM_ATTR_SYS_ID
Attribute Type:    VariableLabel
                  (see Definitions Shared by Protocols)
Attribute Groups:  System Identification Attribute Group
Availability:      Mandatory
```

The IntelliVue monitor uses the 6 byte MAC address as identifier. Future versions might use an 8 byte EUI-64 identifier.

Attribute: Nomenclature Version

The Nomenclature Version attribute contains the version of the nomenclature used by the device.

```
Attribute ID:      NOM_ATTR_NOM_VERS
Attribute Type:    u_32
Attribute Groups:  System Identification Attribute Group
Availability:      Mandatory
```

The Nomenclature Version is composed of 16 bit major and 16 bit minor version number. The IntelliVue monitor currently uses the Nomenclature Version 1.0.

Attribute: System Localization

The System Localization attribute contains information about the language version used by the device.

```
Attribute ID:      NOM_ATTR_LOCALIZN
Attribute Type:    SystemLocal
Attribute Groups:  System Identification Attribute Group
Availability:      Optional
```

The *SystemLocal* is defined as follows:

```
typedef struct {
    u_32          text_catalog_revision;
    Language      language;
    StringFormat  format;
} SystemLocal;
```

The *text_catalog_revision* contains revision information about the texts used by the monitor. The two most significant bytes contain the version of the text catalog (one byte major, one byte minor revision). The text catalog defines the possible values for Attributes of the type *TextId*. A client which depends on a *TextId* having a specific value can use this information for revision control.

The lower two bytes of the *text_catalog_revision* are used for a language revision (one byte major, one byte minor revision). The language revision denotes the mapping from a *TextId* to an actual string in the monitor language.

The *Language* describes the language used by the monitor. It is defined as follows:

```
typedef u_16 Language;
#define LANGUAGE_UNSPEC      0
#define ENGLISH              1
#define GERMAN                2
#define FRENCH                3
#define ITALIAN               4
#define SPANISH               5
#define DUTCH                 6
#define SWEDISH               7
#define FINNISH               8
#define NORWEG                9
#define DANISH                10
#define JAPANESE              11
#define REP_OF_CHINA           12
#define PEOPLE_REP_CHINA      13
#define PORTUGUESE            14
#define RUSSIAN                15
#define BYELORUSSIAN          16
#define UKRAINIAN              17
#define CROATIAN               18
#define SERBIAN                19
#define MACEDONIAN            20
#define BULGARIAN              21
#define GREEK                  22
#define POLISH                 23
#define CZECH                  24
#define SLOVAK                 25
#define SLOVENIAN              26
#define HUNGARIAN              27
#define ROMANIAN               28
#define TURKISH                29
#define LATVIAN                30
#define LITHUANIAN             31
#define ESTONIAN               32
#define KOREAN                 33
```

The *StringFormat* describes how strings are encoded. The IntelliVue monitor uses unicode encoding.

```
typedef u_16 StringFormat;
#define STRFMT_UNICODE_NT    11
```

Attribute: System Specification

The System Specification attribute contains a set of functional components supported by the system.

```
Attribute ID:      NOM_ATTR_SYS_SPECN
Attribute Type:    SystemSpec
Attribute Groups:  System Application Attribute Group
Availability:      Optional
```

The *SystemSpec* is defined as follows:

```
typedef struct {
    u_16      count;
    u_16      length;
    SystemSpecEntry value[1];
} SystemSpec;

typedef struct {
    PrivateOid component_capab_id;
    u_16      length;
    u_16      value[1];
} SystemSpecEntry;
```

The supported components are:

```
Component ID:      NOM_MDIB_OBJ_SUPPORT
Component Type:    MdibObjectSupport
Availability:      Mandatory
```

The *MdibObjectSupport* is defined as follows:

```
typedef struct {
    u_16    count;
    u_16    length;
    MdibObjectSupportEntry    value[1];
} MdibObjectSupport;

typedef struct {
    TYPE    object_type;
    u_32    max_inst;
} MdibObjectSupportEntry;
```

The *MdibObjectSupport* contains a list of all object classes supported by the system and the maximum number of instances per class. If *max_inst* contains 0xffffffff, it is not defined.

Attribute: Mds General System Info

The Mds General System Info attribute contains global information about the monitor and its configuration.

```
Attribute ID:      NOM_ATTR_MDS_GEN_INFO
Attribute Type:    MdsGenSystemInfo
Attribute Group:   System Application Attribute Group
Availability:      Optional
```

The *MdsGenSystemInfo* is defined as follows:

```
typedef struct
{
    u_16    count;
    u_16    length;
    MdsGenSystemInfoEntry    value[1];
} MdsGenSystemInfo;
```

The *MdsGenSystemInfoEntry* allows to encode generic system information. It has the following structure:

```
typedef struct
{
    u_16    choice;
#define MDS_GEN_SYSTEM_INFO_SYSTEM_PULSE_CHOSEN    1
    u_16    length;
    u_8     value[1]; /* placeholder for appended data */
} MdsGenSystemInfoEntry;
```

One *MdsGenSystemInfoEntry* is used to encode the System Pulse information. The monitor can generate a pulse rate from several sources.

```
Choice:      MDS_GEN_SYSTEM_INFO_SYSTEM_PULSE_CHOSEN 1
Type:        SystemPulseInfo
Availability: Optional
```

The *SystemPulseInfo* is defined as follows:

```
typedef struct
{
    ManagedObjectId    system_pulse;
    ManagedObjectId    alarm_source;
} SystemPulseInfo;
```

It enfolds the *ManagedObjectIds* of the object instances selected as system-pulse respectively alarm-source.

Attribute: Production Specification

The Production Specification attribute contains a list of component revisions and serial numbers within the system.

```
Attribute ID:      NOM_ATTR_ID_PROD_SPECN
Attribute Type:    ProductionSpec
Attribute Groups:  System Production Attribute Group
Availability:      Optional
```

The *ProductionSpec* is defined as follows:

```
typedef struct {
    u_16      count;
    u_16      length;
    ProdSpecEntry value[1];
} ProductionSpec;

typedef struct {
    u_16      spec_type;
#define UNSPECIFIED      0
#define SERIAL_NUMBER    1
#define PART_NUMBER      2
#define HW_REVISION      3
#define SW_REVISION      4
#define FW_REVISION      5
#define PROTOCOL_REVISION 6
    PrivateOid      component_id;
    VariableLabel    prod_spec;
} ProdSpecEntry;
```

The current IntelliVue monitor uses 10 characters for a serial number, 14 characters for part numbers and 8 characters for revision strings. The strings are not null-terminated.

The supported components are:

```
Component ID:      ID_COMP_PRODUCT
Description:        Overall product specification

Component ID:      ID_COMP_CONFIG
Description:        Specific system configuration

Component ID:      ID_COMP_BOOT
Description:        Boot code specification

Component ID:      ID_COMP_MAIN_BD
Description:        Mainboard hardware specification

Component ID:      ID_COMP_APPL_SW
Description:        Application software specification
```

See the section “Component IDs” on page 8-194 for the values of the *component_id*. The *ProductionSpec* may contain additional private entries.

To retrieve the IntelliVue monitor software revision, read the ProductSpecEntry with the Component ID “ID_COMP_APPL_SW”. Its prod_spec attribute contains a string of the form “J.00.00” describing the running software revision.

Attribute: MDS Status

The MDS Status attribute describes the device state.

```
Attribute ID:      NOM_ATTR_VMS_MDS_STAT
Attribute Type:    MDSStatus
Attribute Groups:  System Application Attribute Group
Availability:      Mandatory
```

The *MDSStatus* is defined as follows:

```
typedef u_16      MDSStatus;
#define DISCONNECTED 0
#define UNASSOCIATED 1
#define OPERATING    6
```

The *MDSStatus* values have the following meaning:

DISCONNECTED: The IntelliVue monitor is not connected to the network.

UNASSOCIATED: The IntelliVue monitor is connected to the network, but no association is currently active.

OPERATING: The IntelliVue monitor has an association with a Computer Client.

Currently, a Computer Client will only see the MDS Status OPERATING, if the MDS has another Status, there is no association with a Computer Client.

Attribute: Bed Label

The Bed Label attribute contains a printable string identifying the system location.

```
Attribute ID:      NOM_ATTR_ID_BED_LABEL
Attribute Type:    String
                  (see Definitions Shared by Protocols)
Attribute Groups:  System Application Attribute Group
Availability:      Optional
```

The Bed Label can be entered in the Admit/Discharge dialog. It uses 16 bit unicode character encoding. Currently, the Bed Label is 17 characters (including terminating '\0'). If the actual label is shorter, the string is filled with '\0' characters.

Attribute: Operating Mode

The Operating Mode attribute identifies the current operating mode of the device.

```
Attribute ID:      NOM_ATTR_MODE_OP
Attribute Type:    PrivateOID
Attribute Groups:  System Application Attribute Group
Availability:      Optional
```

The Operating Mode is defined as a bit field. The following mode bits are defined:

```
#define OPMODE_UNSPEC 0x8000
#define MONITORING    0x4000
#define DEMO          0x2000
#define SERVICE       0x1000
#define OPMODE_STANDBY 0x0002
#define CONFIG        0x0001
```

The values have the following meaning:

OPMODE_UNSPEC: The Operating Mode is not specified.

MONITORING: Device is configured to monitor patient data (the default mode).

DEMO: Demonstration Mode with simulated patient data.

SERVICE: Device is in Service Mode.

STANDBY: Standby and Power Safe Mode.

CONFIG: Device is in Configuration Mode.

Exactly one of the bit out of the bits 0 - 4 must be set, bits 14 and 15 (the stand-by and config mode bits) can be set optionally.

Attribute: Application Area

The Application Area attribute describes the intended application area for the device.

```
Attribute ID:      NOM_ATTR_AREA_APPL
Attribute Type:    ApplicationArea
Attribute Groups:  System Application Attribute Group
Availability:      Optional
```

The *ApplicationArea* is defined as follows:

```
typedef u_16      ApplicationArea;
#define AREA_UNSPEC      0
#define AREA_OPERATING_ROOM  1
#define AREA_INTENSIVE_CARE  2
#define AREA_NEONATAL_INTENSIVE_CARE 3
#define AREA_CARDIOLOGY_CARE 4
```

The values have the following meaning:

AREA_UNSPEC: The application area has not been specified.

AREA_OPERATING_ROOM: The application area has been specified as an operating room.

AREA_INTENSIVE_CARE: The application area has been specified as an intensive care unit.

AREA_NEONATAL_INTENSIVE_CARE: The application area has been specified as a neonatal intensive care unit.

AREA_CARDIOLOGY_CARE: The application area has been specified as a cardiology care unit.

Attribute: Date and Time

The Date and Time attribute contains the current device time.

```
Attribute ID:      NOM_ATTR_TIME_ABS
Attribute Type:    AbsoluteTime
                  (see Definitions Shared by Protocols)
Attribute Groups:  System Application Attribute Group
Availability:      Optional
```

Attribute: Relative Time

The Relative Time attribute contains the current device relative time.

```
Attribute ID:      NOM_ATTR_TIME_REL
Attribute Type:    RelativeTime
                  (see Definitions Shared by Protocols)
Attribute Groups:  System Application Attribute Group
Availability:      Optional
```

The Relative Time is set to zero after each power cycle.

Attribute: Altitude

The Altitude attribute contains the system altitude above or below sea level.

```
Attribute ID:      NOM_ATTR_ALTITUDE
Attribute Type:    i_16
Attribute Groups:  System Application Attribute Group
Availability:      Optional
```

Attribute: Line Frequency

The Line Frequency attribute describes the frequency of the main power supply in Hz.

```
Attribute ID:      NOM_ATTR_LINE_FREQ
Attribute Type:    LineFrequency
Attribute Groups:  System Application Attribute Group
Availability:      Optional
```

The *LineFrequency* is defined as follows:

```
typedef u_16      LineFrequency;
#define LINE_F_UNSPEC 0
#define LINE_F_50HZ  1
#define LINE_F_60HZ  2
```

Attribute: Association Invoke ID

The Association Invoke ID attribute is a counter for the number of associations. It is incremented with each new association.

```
Attribute ID:      NOM_ATTR_ID_ASSOC_NO
Attribute Type:    u_16
Attribute Groups:  System Identification Attribute Group
Availability:      Optional
```

Attribute Groups

The attributes of the Medical Device System object are arranged in the following attribute groups:

```
Attribute Group:  System Identification Attribute Group
Group ID:         NOM_ATTR_GRP_SYS_ID
Description:      Identification of the system
Attributes:       System Type, System Model, System Id,
                  Nomenclature Version, System Localization, Association Invoke Id
```

```
Attribute Group:  System Application Attribute Group
Group ID:         NOM_ATTR_GRP_SYS_APPL
Description:      System Capabilities and Settings
Attributes:       System Specification, MDS Status, Bed Label,
                  Operating Mode, Application Area, Data and
                  Time, Relative Time, Altitude, Line
                  Frequency, Mds General System Info
```

```
Attribute Group:  System Production Attribute Group
Group ID:         NOM_ATTR_GRP_SYS_PROD
Description:      HW and SW configuration
Attributes:       Production Specification
```

Alert Monitor Object

Attributes of the Alert Monitor Object

This section defines the attributes of the Alert Monitor object, together with the attribute identifier codes and attribute data types.

The Alert Monitor object represents the overall device alert condition. It contains a global alert status and a list of active technical and patient alerts.

Attribute: Handle

The Handle attribute identifies the Alert Monitor object in the form of a numeral value. The Handle is unique within a device context (see Common Data Type - Global Handle). The actual value of the Handle attribute does not have a meaning. It is used for reference and relation purposes.

```
Attribute ID:      NOM_ATTR_ID_HANDLE
Attribute Type:    Handle (see Definitions Shared by Protocols)
Attribute Groups:  VMO Static Context Group
Availability:      Mandatory
```

Attribute: Type

The Type attribute contains an identification of the object type represented by the Alert Monitor.

```
Attribute ID:      NOM_ATTR_ID_TYPE
Attribute Type:    TYPE (see Definitions Shared by Protocols)
Attribute Groups:  VMO Static Context Group
Availability:      Mandatory
```

Attribute: Device Alert Condition

The Device Alert Condition attribute contains global device alert status information.

```
Attribute ID:      NOM_ATTR_DEV_AL_COND
Attribute Type:    DeviceAlertCondition
Attribute Groups:  Alert Monitor Group
Availability:      Mandatory
```

The *DeviceAlertCondition* is defined as follows:

```
typedef struct {
    AlertState    device_alert_state;
    u_16          al_stat_chg_cnt;
    AlertType     max_p_alarm;
    AlertType     max_t_alarm;
    AlertType     max_aud_alarm;
} DeviceAlertCondition;
```

The *AlertState* is a bit field defined as follows:

```
typedef u_16 AlertState;
#define AL_INHIBITED          0x8000
#define AL_SUSPENDED         0x4000
#define AL_LATCHED           0x2000
#define AL_SILENCED_RESET    0x1000
#define AL_DEV_IN_TEST_MODE   0x0400
#define AL_DEV_IN_STANDBY     0x0200
#define AL_DEV_IN_DEMO_MODE   0x0100
#define AL_NEW_ALERT          0x0008
```

The *AlertState* is used for the overall device alert state and for the specific state of each alert. The bits in *AlertState* have the following meaning:

AL_INHIBITED: Alert is switched off.

AL_SUSPENDED: Alert inactivated temporarily, alert condition is acknowledged.

AL_LATCHED: Alert condition is not active but latched, note that technical alarms are never latching.

AL_SILENCED_RESET: Alert condition stopped but alarming re-enabled (only for *DeviceAlertCondition*).

AL_DEV_IN_TEST_MODE: Device is in a temporary test mode.

AL_DEV_IN_STANDBY: Device is in standby mode.

AL_DEV_IN_DEMO_MODE: Indicates that the device is in demo mode.

AL_NEW_ALERT: Indicate a new alarm (not in *DeviceAlertCondition*). A Computer Client might not see this bit if it does not poll fast enough or other delays occur.

The *al_stat_chg_cnt* is an internal change counter. A Computer Client should not interpret this field, because it can not be guaranteed that no internal message is missed.

The *AlertType* is a bit field defined as follows:

```
typedef u_16 AlertType;
#define NO_ALERT 0
#define LOW_PRI_T_AL 1
#define MED_PRI_T_AL 2
#define HI_PRI_T_AL 4
#define LOW_PRI_P_AL 256
#define MED_PRI_P_AL 512
#define HI_PRI_P_AL 1024
```

Intellivue monitors with software revision E.0 or higher allow changing of the inop severity for various inop alarms. These changes are reflected in the *AlertType* bitfield.

The bits have the following meaning:

NO_ALERT: No alert active.

LOW_PRI_T_AL: Low priority technical alarm (soft inop). These inops are generated after a signal analysis (e.g. “Noisy ECG”).

MED_PRI_T_AL: Medium priority technical alarm (hard inop). These inops are generated during inoperable parameter measurement because of hardware faults or no transducer connected (e.g. “Leads Off”, “ABP No Transducer”).

HI_PRI_T_AL: High priority technical alarm (severe inop).

LOW_PRI_P_AL: Awareness Condition (short yellow alarm): These alarms are marked with a “**” in the alarm string and a specific short yellow alarm sound is issued. Today short yellow alarms are generated only from arrhythmia computer.

MED_PRI_P_AL: Medium priority patient alarm (yellow alarm): These alarms are marked with a “***” in the alarm string. They indicate a less critical patient condition usually due to violation of user defined criteria (e.g. limit violation alarm).

HI_PRI_P_AL: High priority patient alarm (red alarm): These alarms are marked with a “****” in the alarm string. These alarms indicate a life threatening patient condition.

Attribute: Device T-Alarm List

The Device T-Alarm List attribute contains the active technical alarms (inops) in the system.

```
Attribute ID:      NOM_ATTR_AL_MON_T_AL_LIST
Attribute Type:    DevAlarmList
Attribute Groups:  Alert Monitor Group
Availability:      Mandatory
```

The *DevAlarmList* is defined as follows:

```
typedef struct {
    u_16          count;
    u_16          length;
    DevAlarmEntry value[1];
} DevAlarmList;

typedef struct {
    OIDType      al_source;
    OIDType      al_code;
    AlertType     al_type;
    AlertState    al_state;
    ManagedObjectId object;
    PrivateOid    alert_info_id;
#define GEN_ALMON_INFO 513
#define STR_ALMON_INFO 516
    u_16          length;
} DevAlarmEntry;
```

The *al_source* is taken from the Object Oriented or the SCADA partition (depending on *al_code*). It identifies the origin of the alert (e.g. temperature).

The *al_code* is taken from the Events partition and describes the reason for the alert (e.g. high alarm). The least significant bit is used to define the nomenclature partition for *al_source*. Last bit 0 means SCADA partition, last bit 1 means Object Oriented partition.

The definitions for *AlertType* and *AlertState* can be found in the paragraph about the Device Alert Condition.

The *object* field contains a reference to the object which generated the alert. The object may not be known to the Computer Client, if the Data Export protocol does not allow accessing the specific object.

If the *alert_info_id* is set to GEN_ALMON_INFO, an *AlMonGenInfo* structure is appended:

```
typedef struct {
    u_16          al_inst_no;
    TextId        al_text;
    AlertPriority  priority;
    AlertFlags     flags;
} AlMonGenInfo;
```

If the *alert_info_id* is set to STR_ALMON_INFO, an *StrAlMonInfo* structure is appended:

```
typedef struct {
    u_16          al_inst_no;
    TextId        al_text;
    AlertPriority  priority;
    AlertFlags     flags;
    String        string;
} StrAlMonInfo;
```

Currently, the IntelliVue monitor only supports the *StrAlMonInfo* data type.

The *al_inst_no* is a private ID.

The *al_text* is a private ID.

The *AlertPriority* is defined as follows:

```
typedef u_16      AlertPriority;
```

The *AlertPriority* only allows prioritization within a group of alarms. A Computer Client application should use the *AlertType* to distinguish low and high priority alarms.

The *AlertFlags* type is defined as follows:

```
typedef u_16      AlertFlags;
#define BEDSIDE_AUDIBLE      0x4000
#define CENTRAL_AUDIBLE      0x2000
#define VISUAL_LATCHING      0x1000
#define AUDIBLE_LATCHING      0x0800
#define SHORT_YELLOW_EXTENSION 0x0400
#define DERIVED              0x0200
```

The bits in the *AlertFlag* have the following meaning:

BEDSIDE_AUDIBLE: Alert sound at the bedside

CENTRAL_AUDIBLE: Alert sound at the central station

VISUAL_LATCHING: Alert is visible after the alarm condition has ceased. The alarm indication will exist until a specific action is taken by a user (e.g. Silence/Reset).

AUDIBLE_LATCHING: Alert is sound issued after the alarm condition has ceased. The alarm indication will exist until a specific action is taken by a user (e.g. Silence/Reset).

SHORT_YELLOW_EXTENSION: Alarm is not active but artificially extended for short yellow behavior.

DERIVED: Derived alarm.

The *String* contains the a description of the alarm in the language supported by the IntelliVue monitor. *Strings* for patient alarms are prefixed with two "***" or three "****" alarm stars (see “Definitions Shared by Protocols” on page 6-35 for UNICODE character encoding). Currently, the String is 19 characters long, including the terminating '\0'.

Attribute: Device P-Alarm List

The Device P-Alarm List attribute contains the active patient alarm in the system.

```
Attribute ID:      NOM_ATTR_AL_MON_P_AL_LIST
Attribute Type:    DevAlarmList
Attribute Groups:  Alert Monitor Group
Availability:      Mandatory
```

The *DevAlarmList* data type is the same as for the Device T-Alarm List.

The data in a Device T-Alarm List or Device P-Alarm List might be too large to fit in a single message. In this case the Remote Operation Linked Result message will be used (see “Remote Operation Linked Result” on page 6-44). In this case each message will contain a correctly formatted Alarm list and the Computer Client must merge the lists to get the complete Device T-Alarm List or Device P-Alarm List.

Attribute Groups

The attributes of the Alert Monitor object are arranged in the following attribute groups:

Attribute Group:	VMO Static Context Group
Group ID:	NOM_ATTR_GRP_VMO_STATIC
Description:	Static context of the object
Attributes:	TYPE, Handle
Attribute Group:	Alert Monitor Group
Group ID:	NOM_ATTR_GRP_AL_MON
Description:	Alarm related attributes
Attributes:	Device Alert Condition, Device P-Alarm List, Device T-Alarm List

Patient Demographics Object

Attributes of the Patient Demographic Object

This section defines the attributes of the Patient Demographics object, together with the attribute identifier codes and attribute data types.

The Patient Demographics object contains the patient information present in the system.

Attribute: Handle

The Handle attribute identifies the Patient Demographics object in the form of a numeral value. The Handle is unique within a device context (see Common Data Type - Global Handle). The actual value of the Handle attribute does not have a meaning. It is used for reference and relation purposes.

Attribute ID:	NOM_ATTR_ID_HANDLE
Attribute Type:	Handle (see Definitions Shared by Protocols)
Attribute Groups:	Patient Demographics Attribute Group
Availability:	Mandatory

Attribute: Pat Demo State

The Pat Demo State attribute describes the current state of the Patient Demographics object.

Attribute ID:	NOM_ATTR_PT_DEMOG_ST
Attribute Type:	PatDemoState
Attribute Groups:	Patient Demographics Attribute Group
Availability:	Mandatory

The *PatDemoState* is defined as follows:

```
typedef u_16      PatDmgState;
#define EMPTY      0
#define PRE_ADMITTED 1
#define ADMITTED   2
#define DISCHARGED 8
```

The values have the following meaning:

EMPTY: No patient information present.

PRE_ADMITTED: Currently not used.

ADMITTED: Patient information is present and valid.

DISCHARGED: Data is still available, but patient is no longer assigned to device.

Attribute: Patient Type

The Patient Type attribute describes the type of patient admitted to the system.

```
Attribute ID:      NOM_ATTR_PT_TYPE
Attribute Type:    PatientType
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

The Patient Type is defined as follows:

```
typedef u_16      PatientType;
#define PAT_TYPE_UNSPECIFIED 0
#define ADULT 1
#define PEDIATRIC 2
#define NEONATAL 3
```

The Patient Type can be set by the user in the Admit/Discharge dialog (Patient Cat.).

Attribute: Patient Paced Mode

The Patient Paced Mode attribute indicates whether the patient is paced or not.

```
Attribute ID:      NOM_ATTR_PT_PACED_MOD
Attribute Type:    PatPacedMode
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

The *PatPacedMode* is defined as follows:

```
typedef u_16      PatPacedMode;
#define PAT_NOT_PACED 0
#define PAT_PACED_GEN 1
```

Values greater one are reserved to indicate special paced modes. The Computer Client should test for "==" 0" or "!= 0".

Attribute: Given Name

The Given Name attribute contains the first name of the patient.

```
Attribute ID:      NOM_ATTR_PT_NAME_GIVEN
Attribute Type:    String
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Currently, the Given Name can be up to 19 characters long, including the terminating '\0'.

Attribute: Middle Name

The Middle Name attribute contains the middle name of the patient.

```
Attribute ID:      NOM_ATTR_PT_NAME_MIDDLE
Attribute Type:    String
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Currently, the Middle Name can be up to 19 characters long, including terminating '\0'.

Attribute: Family Name

The Family Name attribute contains the last name of the patient.

```
Attribute ID:      NOM_ATTR_PT_NAME_FAMILY
Attribute Type:    String
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Currently, the Family Name can be up to 19 characters long, including terminating '\0'.

Attribute: Patient ID

The Patient ID attribute contains the ID of the patient.

```
Attribute ID:      NOM_ATTR_PT_LIFETIME_ID (identical to previous NOM_ATTR_PT_ID)
Attribute Type:    String
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Currently, the Patient ID (Medical Record Number - MRN) can be up to 17 characters long, including the terminating '\0'.

Since Rev. G Intellivue the Patient ID is called Lifetime ID.

Attribute: Encounter ID

The Encounter ID attribute contains the ID of the current visit of the patient.

```
Attribute ID:      NOM_ATTR_PT_ENCOUNTER_ID
Attribute Type:    String
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Currently, the Encounter ID can be up to 17 characters long, including the terminating '\0'.

Attribute: Patient Sex

The Patient Sex attribute contains the sex of the patient.

```
Attribute ID:      NOM_ATTR_PT_SEX
Attribute Type:    PatientSex
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

The *PatientSex* is described as follows:

```
typedef u_16      PatientSex;
#define SEX_UNKNOWN      0
#define MALE             1
#define FEMALE           2
#define SEX_UNSPECIFIED  9
```

The values have the following meaning:

SEX_UNKNOWN: Patient sex is not known

MALE: Patient is male

FEMALE: Patient is female

SEX_UNSPECIFIED: Patient sex is not specified

Attribute: Date of Birth

The Date of Birth attribute contains the Date of Birth of the patient.

```
Attribute ID:      NOM_ATTR_PT_DOB
Attribute Type:    AbsoluteTime
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Attribute: Patient Height

The Patient Height attribute contains the height of the patient.

```
Attribute ID:      NOM_ATTR_PT_HEIGHT
Attribute Type:    PatMeasure
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

The *PatMeasure* is defined as follows:

```
typedef struct {
    FLOATType    value;
    OIDType      m_unit;
} PatMeasure;
```

The *value* contains the actual value of the attribute and the *m_units* indicates the unit of measurement for the *value*.

Attribute: Patient Weight

The Patient Height attribute contains the weight of the patient.

```
Attribute ID:      NOM_ATTR_PT_WEIGHT
Attribute Type:    PatMeasure
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Attribute: Patient Age

The Patient Age attribute contains the age of the patient.

```
Attribute ID:      NOM_ATTR_PT_AGE
Attribute Type:    PatMeasure
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Attribute: Patient BSA

The Patient BSA attribute contains the body surface area of the patient.

```
Attribute ID:      NOM_ATTR_PT_BSA
Attribute Type:    PatMeasure
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Attribute: Patient BSA Formula

The Patient BSA Formula attribute describes the formula which is used for the calculation of the patient body surface area.

```
Attribute ID:      NOM_ATTR_PT_BSA_FORMULA
Attribute Type:    PatBsaFormula
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

The *PatBsaFormula* is described as follows:

```
typedef u_16      PtBsaFormula;
#define BSA_FORMULA_UNSPEC    0
#define BSA_FORMULA_BOYD      1
#define BSA_FORMULA_DUBOIS    2
```

The values have the following meaning:

BSA_FORMULA_UNSPEC: Formula not specified

BSA_FORMULA_BOYD: BSA calculation according to Boyd

BDA_FORMULA_DUBOIS: BSA calculation according to Dubois

Attribute: Notes1

The Notes1 attribute provides additional information about the patient.

```
Attribute ID:      NOM_ATTR_PT_NOTES1
Attribute Type:    String
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Currently, the Notes1 field can be up to 31 characters long, including the terminating '\0'.

Attribute: Notes2

The Notes2 attribute provides additional information about the patient.

```
Attribute ID:      NOM_ATTR_PT_NOTES2
Attribute Type:    String
Attribute Groups:  Patient Demographics Attribute Group
Availability:      Optional
```

Currently, the Notes2 field can be up to 31 characters long, including the terminating '\0'.

Attribute Groups

The attributes of the Patient Demographics object are arranged in the following attribute groups:

```
Attribute Group:  Patient Demographics Attribute Group
Group ID:         NOM_ATTR_GRP_PT_DEMOG
Description:      Attributes containing patient information
Attributes:       all attributes
```

Patient Conflict Handling

The patient information is stored in the monitor, the measurement server and the central station (if present). This can lead to patient conflicts when the patient information in these locations differ. If the IntelliVue monitor detects a patient conflict, it will display a “Patient Selection” window which allows the user to resolve the conflict.

In the case of a patient conflict, the behavior of the Data Export software is as follows:

- If the Patient Type or Patient Paced Mode attribute is different, the data from the measurement server is considered as relevant.
- If the patient is different (devices have been disconnected and a new patient has been admitted), the Patient Type and Patient Paced Mode information from the measurement server is exported. The other attributes are cleared and the Family Name attribute is set to “???”.

Connect Indication Attributes

This section describes the attributes contained in the Connect Indication Message.

Attribute: System Type

The System Type attribute describes the type of the system (e.g. Monitor).

```
Attribute ID:      NOM_ATTR_SYS_TYPE
Attribute Type:    TYPE (see Definitions Shared by Protocols)
Attribute Groups:  -
Availability:      -
```

Attribute: Protocol Support

The Protocol Support contains an entry for each protocol supported on the network interface.

```
Attribute ID:      NOM_ATTR_PCOL_SUPPORT
Attribute Type:    ProtoSupport
Attribute Groups:  -
Availability:      -
```

The *ProtoSupport* is defined as follows:

```
typedef struct {
    u_16          count;
    u_16          length;
    ProtoSupportEntry value[1];
} ProtoSupport;

typedef struct {
    ApplProtoId    appl_proto;
    TransProtoId   trans_proto;
    u_16           port_number;
    ProtoOptions    options;
} ProtoSupportEntry;

typedef u_16 ApplProtoId;
#define AP_ID_ACSE 1
#define AP_ID_DATA_OUT 5

typedef u_16 TransProtoId;
#define TP_ID_UDP 1

typedef u_16 ProtoOptions;
#define P_OPT_WIRELESS 0x8000
```

The Computer Client should parse the available protocols and search for the AP_ID_DATA_OUT. This entry specifies the port for the Data Export Protocol. The corresponding Association Control Protocol runs on the same port.

The Computer Client must only send requests to the port specified for the Data Export Protocol.

Attribute: System Localization

The System Localization attribute describes the handling of natural language items.

```
Attribute ID:      NOM_ATTR_LOCALIZN
Attribute Type:    SystemLocal
Attribute Groups:  -
Availability:      -
```

The *SystemLocal* is defined as follows:

```
typedef struct {
    u_32          syslocal_revision;
    Language       language;
    StringFormat   format;
} SystemLocal;
```

The *syslocal_revision* contains the revision of the text catalog used for internal texts.

The *Language* describes the language used in any String type. It is defined as follows:

```
typedef u_16          Language;
#define LANGUAGE_UNSPEC 0
#define ENGLISH        1
#define GERMAN         2
#define FRENCH         3
#define ITALIAN        4
#define SPANISH        5
#define DUTCH          6
#define SWEDISH        7
#define FINNISH        8
#define NORWEG         9
#define DANISH        10
#define JAPANESE       11
#define REP_OF_CHINA   12
#define PEOPLE_REP_CHINA 13
#define PORTUGUESE     14
#define RUSSIAN        15
#define BYELORUSSIAN   16
#define UKRAINIAN      17
#define CROATIAN       18
#define SERBIAN        19
#define MACEDONIAN     20
#define BULGARIAN      21
#define GREEK          22
#define POLISH         23
#define CZECH          24
#define SLOVAK         25
#define SLOVENIAN      26
#define HUNGARIAN      27
#define ROMANIAN       28
#define TURKISH        29
#define LATVIAN        30
#define LITHUANIAN     31
#define ESTONIAN       32
#define KOREAN         33
```

The *StringFormat* defines the format used for the *String* data type. The IntelliVue monitor uses 16bit Unicode characters.

```
typedef u_16          StringFormat;
#define STRFMT_UNICODE_NT 11
```

Attribute: IP Address Information

The IP Address Information attribute identifies the network interface of the IntelliVue monitor.

```
Attribute ID:      NOM_ATTR_NET_ADDR_INFO
Attribute Type:    IpAddressInfo
Attribute Groups:  -
Availability:      -
```

The *IpAddressInfo* is defined as follows:

```
typedef struct IpAddressInfo {
    MACAddress    mac_address;
    IPAddress     ip_address;
    IPAddress     subnet_mask;
} IpAddressInfo;

typedef struct MACAddress {
    u_8           value[6];
} MACAddress;

typedef struct IPAddress {
    u_8           value[4];
} IPAddress;
```

Partition IDs

The following sections contain a list of identifiers which are used within the IntelliVue monitor. Each identifier is unique within a given partition.

```
#define NOM_PART_OBJ                1
    /* Object Oriented Elements */
#define NOM_PART_SCADA              2
    /* Physiological Measurements */
#define NOM_PART_EVT                3
    /* Events for Alerts */
#define NOM_PART_DIM                4
    /* Units of Measurement */
#define NOM_PART_PGRP               6
    /* Identification of Parameter Groups */
#define NOM_PART_INFRASTRUCT        8
    /* Infrastructure Elements */
#define NOM_PART_EMFC               1025
    /* EMFC */
#define NOM_PART_SETTINGS           1026
    /* Settings */
```

Object Classes

The following IDs identify object types. They are taken from the Object Oriented Elements partition. These objects may be the source of alerts (see “Alert Monitor Object” on page 99).

NOM_MOC_VMO	1
VMO	
NOM_MOC_VMO_METRIC_NU	6
Numeric	
NOM_MOC_VMO_METRIC_SA_RT	9
Realtime Sample Array	
NOM_MOC_VMS_MDS	33
MDS	
NOM_MOC_VMS_MDS_COMPOS_SINGLE_BED	35
Composit Single Bed MDS	
NOM_MOC_VMS_MDS_SIMP	37
Simple MDS	
NOM_MOC_BATT	41
Battery	
NOM_MOC_PT_DEMOG	42
Patient Demographics	
NOM_MOC_VMO_AL_MON	54
Alert Monitor	
NOM_ACT_POLL_MDIB_DATA	3094
Poll Action	
NOM_NOTI_MDS_CREAT	3334
MDS Create	
NOM_NOTI_CONN_INDIC	3351
Connect Indication	
NOM_DEV_METER_CONC_SKIN_GAS	4264
Skin Gas	
NOM_DEV_METER_FLOW_BLD	4284
Blood Flow	
NOM_DEV_ANALY_CONC_GAS_MULTI_PARAM_MDS	4113
Gas Analyzer	
NOM_DEV_ANALY_CONC_GAS_MULTI_PARAM_VMD	4114
Gas	
NOM_DEV_METER_CONC_SKIN_GAS_MDS	4265
Skin Gas	
NOM_DEV_MON_PHYSIO_MULTI_PARAM_MDS	4429
Multi-Param	
NOM_DEV_PUMP_INFUS_MDS	4449
Pump Infus	
NOM_DEV_SYS_PT_VENT_MDS	4465
Ventilator	
NOM_DEV_SYS_PT_VENT_VMD	4466
Ventilator	
NOM_DEV_SYS_MULTI_MODAL_MDS	4493
Multi-Modal MDS	
NOM_DEV_SYS_MULTI_MODAL_VMD	4494
Multi-Modal	
NOM_DEV_SYS_VS_CONFIG_MDS	5209
config MDS	
NOM_DEV_SYS_VS_UNCONFIG_MDS	5213
unconfig MDS	
NOM_DEV_ANALY_SAT_O2_VMD	4106
sat O2	
NOM_DEV_ANALY_CONC_GAS_MULTI_PARAM_VMD	4114
Gas	
NOM_DEV_ANALY_FLOW_AWAY_VMD	4130
Flow Away	
NOM_DEV_ANALY_CARD_OUTPUT_VMD	4134
C.O.	

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Physiological Identifier

A Physiological Identifier denotes the origin of a physiological measurement. The identifiers are located in the SCADA partition. The Physiological Identifier is transmitted as part of the numeric or wave observed value. The Physiological Identifier may not be unique. However, it is guaranteed that the Label ID is unique. The Label ID is mapped to a Label String based on the text catalogue (see “Attribute: System Localization” on page 92). Note that the mapping listed below may not be complete and is subject to changes and additions, due to revision changes and additions from additional interfaced devices. The table below should be viewed as an example.

The list below shows the numerics and waves which are supported by the monitor. The numerics and waves are sorted according to their internal priority, i.e. numerics or waves with a higher priority are listed first. This information depends heavily on the software revision of the monitor and the connected devices. Especially data coming from a VueLink module depends on the version of the VueLink driver and the specification of the connected external device. The list contains the possible unit codes for the numerics and waves. The unit codes for numerics/waves acquired through data import interfaces (e.g. VueLink) are not documented, because this data depends on the implementation of the specific data import driver.

For a given software revision, the IntelliVue monitor may not export all of the numerics specified below. The IntelliVue monitor may export numerics, which are not specified here. If a numeric is exported also depends on the configuration of the monitor. In general, a numeric will only be available if the required measurement module is connected and if the specific measurement is activated. Some measurements require the presents of more than one measurement module or special configuration steps may be necessary to activate the measurement.

Numerics

HR	Heart Rate	
	Label:	
	NLS_NOM_ECG_CARD_BEAT_RATE	0x00024182
	Observed Value:	
	NOM_ECG_CARD_BEAT_RATE	0x4182
btbHR	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Cardiac Beat-to-Beat Rate	
	Label:	
	NLS_NOM_ECG_CARD_BEAT_RATE_BTBT	0x0002418A
PVC	Observed Value:	
	NOM_ECG_CARD_BEAT_RATE_BTBT	0x418A
	Premature Ventricular Contractions	
	Label:	
	NLS_NOM_ECG_V_P_C_CNT	0x00024261
ST	Observed Value:	
	NOM_ECG_V_P_C_CNT	0x4261
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	ST generic label	
	Label:	
	NLS_NOM_ECG_AMPL_ST	0x00020300
	Compound Observed Value:	
	NOM_ECG_AMPL_ST_I	0x0301
	NOM_ECG_AMPL_ST_II	0x0302
	NOM_ECG_AMPL_ST_III	0x033D
	NOM_ECG_AMPL_ST_AVR	0x033E
	NOM_ECG_AMPL_ST_AVL	0x033F
	NOM_ECG_AMPL_ST_AVF	0x0340
	NOM_ECG_AMPL_ST_V	0x0343
	NOM_ECG_AMPL_ST_MCL	0x034B

	NOM_ECG_AMPL_ST_V1	0x0303
	NOM_ECG_AMPL_ST_V2	0x0304
	NOM_ECG_AMPL_ST_V3	0x0305
	NOM_ECG_AMPL_ST_V4	0x0306
	NOM_ECG_AMPL_ST_V5	0x0307
	NOM_ECG_AMPL_ST_V6	0x0308
	NOM_ECG_AMPL_ST_AS	0x0365
	NOM_ECG_AMPL_ST_ES	0x0364
	NOM_ECG_AMPL_ST_AI	0x0366
	Units:	
	NOM_DIM_MILLI_M	0x0512
STindx	ST Index	
	Label:	
	NLS_NOM_ECG_AMPL_ST_INDEX	0x0002F03D
	Observed Value:	
	NOM_ECG_AMPL_ST_INDEX	0xF03D
	Units:	
	NOM_DIM_MILLI_M	0x0512
QTc	Label:	
	NLS_NOM_ECG_TIME_PD_QTc	0x00023F24
	Observed Value:	
	NOM_ECG_TIME_PD_QTc	0x3F24
	Units:	
	NOM_DIM_MILLI_SEC	0x0892
DeltaQTc	Label:	
	NLS_NOM_ECG_TIME_PD_QTc_DELTA	0x0002F156
	Observed Value:	
	NOM_ECG_TIME_PD_QTc_DELTA	0xF156
	Units:	
	NOM_DIM_MILLI_SEC	0x0892
QT	Label:	
	NLS_NOM_ECG_TIME_PD_QT_GL	0x00023F20
	Observed Value:	
	NOM_ECG_TIME_PD_QT_GL	0x3F20
	Units:	
	NOM_DIM_MILLI_SEC	0x0892
QT-HR	QT HEARTRATE	
	Label:	
	NLS_NOM_ECG_TIME_PD_QT_HEART_RATE	0x0002F154
	Observed Value:	
	NOM_ECG_TIME_PD_QT_HEART_RATE	0xF154
	Units:	
	NOM_DIM_BEATS_PER_MIN	
QT Bsl	Label:	
	NLS_NOM_ECG_TIME_PD_QT_BASELINE	0x0002F155
	Observed Value:	
	NOM_ECG_TIME_PD_QT_BASELINE	0xF155
	Units:	
	NOM_DIM_MILLI_SEC	0x0892
QTHRBl	QT BASELINE HEARTRATE	
	Label:	
	NLS_NOM_ECG_TIME_PD_QT_BASELINE_HEART_RATE	0x0002F157
	Observed Value:	
	NOM_ECG_TIME_PD_QT_BASELINE_HEART_RATE	0xF157
	Units:	
	NOM_DIM_MILLI_SEC	0x0892
Pulse	Pulse Rate	
	Label:	
	NLS_NOM_PULS_RATE	0x0002480A
	Observed Value:	
	NOM_PULS_RATE	0x480A

SpO2	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Arterial Oxygen Saturation	
	Label:	
	NLS_NOM_PULS_OXIM_SAT_O2	0x00024BB8
	Observed Value:	
	NOM_PULS_OXIM_SAT_O2	0x4BB8
	Units:	
Pulse	NOM_DIM_PERCENT	0x0220
	Pulse Rate from Plethysmogram	
	Label:	
	NLS_NOM_PULS_OXIM_PULS_RATE	0x00024822
	Observed Value:	
	NOM_PLETH_PULS_RATE	0x4822
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
SpO2pr	Oxygen Saturation	
	Label:	
	NLS_NOM_PULS_OXIM_SAT_O2_PRE_DUCTAL	0x0002F1C0
	Observed Value:	
	NOM_PULS_OXIM_SAT_O2_PRE_DUCTAL	0xF1C0
	Units:	
	NOM_DIM_PERCENT	0x0220
	Pulse Rate from Plethysmogram (pre ductal)	
Pulse	Label:	
	NLS_SPO2_NAMES_PULS_OXIM_PULS_RATE_PRE_DUCTAL	0x8015543D
	Observed Value:	
	NOM_PLETH_PULS_RATE	0x4822
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Oxygen Saturation	
	Label:	
SpO2po	NLS_NOM_PULS_OXIM_SAT_O2_POST_DUCTAL	0x0002F1D4
	Observed Value:	
	NOM_PULS_OXIM_SAT_O2_POST_DUCTAL	0xF1D4
	Units:	
	NOM_DIM_PERCENT	0x0220
	Pulse Rate from Plethysmogram (post ductal)	
	Label:	
	NLS_SPO2_NAMES_PULS_OXIM_PULS_RATE_POST_DUCTAL	0x80155440
Pulse	Observed Value:	
	NOM_PLETH_PULS_RATE	0x4822
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
%SpO2T	SpO2 parameter label as sourced by the Telemetry system	
	Label:	
	NLS_NOM_PULS_OXIM_SAT_O2_TELE	0x0002F09C
	Observed Value:	
	NOM_PULS_OXIM_SAT_O2_TELE	0xF09C
	Units:	
	NOM_DIM_PERCENT	0x0220
	Pulse parameter label as sourced by the Telemetry system	
PulseT	Label:	
	NLS_NOM_PULS_OXIM_PULS_RATE_TELE	0x0002F09D
	Observed Value:	
	NOM_PULS_OXIM_PULS_RATE_TELE	0xF09D
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Arterial Oxygen Saturation (right)	
	Label:	
SpO2 r	NLS_NOM_PULS_OXIM_SAT_O2_ART_RIGHT	0x00024BCC
	Observed Value:	
	NOM_PULS_OXIM_SAT_O2_ART_RIGHT	0x4BCC
	Units:	
	NOM_DIM_PERCENT	0x0220

Pulse	Pulse Rate from Plethysmogram (right)	
	Label:	
	NLS_SPO2_NAMES_PULS_OXIM_PULS_RATE_RIGHT	0x80155402
	Observed Value:	
	NOM_PLETH_PULS_RATE	0x4822
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
SpO2 l	Arterial Oxygen Saturation (left)	
	Label:	
	NLS_NOM_PULS_OXIM_SAT_O2_ART_LEFT	0x00024BC8
	Observed Value:	
	NOM_PULS_OXIM_SAT_O2_ART_LEFT	0x4BC8
	Units:	
	NOM_DIM_PERCENT	0x0220
Pulse	Pulse Rate from Plethysmogram (left)	
	Label:	
	NLS_SPO2_NAMES_PULS_OXIM_PULS_RATE_LEFT	0x80155401
	Observed Value:	
	NOM_PLETH_PULS_RATE	0x4822
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
DeltaSpO2	Difference between two SpO2 Values (like Left - Right)	
	Label:	
	NLS_NOM_PULS_OXIM_SAT_O2_DIFF	0x00024BC4
	Observed Value:	
	NOM_PULS_OXIM_SAT_O2_DIFF	0x4BC4
	Units:	
	NOM_DIM_PERCENT	0x0220
Perf	Perfusion Indicator	
	Label:	
	NLS_NOM_PULS_OXIM_PERF_REL	0x00024BB0
	Observed Value:	
	NOM_PULS_OXIM_PERF_REL	0x4BB0
	Units:	
	NOM_DIM_DIMLESS	0x0200
PerfPr	Relative Perfusion Left	
	Label:	
	NLS_NOM_PULS_OXIM_PERF_REL_PRE_DUCTAL	0x0002F22C
	Observed Value:	
	NOM_PULS_OXIM_PERF_REL_PRE_DUCTAL	0xF22C
	Units:	
	NOM_DIM_DIMLESS	0x0200
PerfPo	Relative Perfusion Left	
	Label:	
	NLS_NOM_PULS_OXIM_PERF_REL_POST_DUCTAL	0x0002F1DC
	Observed Value:	
	NOM_PULS_OXIM_PERF_REL_POST_DUCTAL	0xF1DC
	Units:	
	NOM_DIM_DIMLESS	0x0200
Perf T	Perf from Telemetry	
	Label:	
	NLS_NOM_PULS_OXIM_PERF_REL_TELE	0x0002F12C
	Observed Value:	
	NOM_PULS_OXIM_PERF_REL_TELE	0xF12C
	Units:	
	NOM_DIM_DIMLESS	0x0200
Perf r	Relative Perfusion Right label	
	Label:	
	NLS_NOM_PULS_OXIM_PERF_REL_RIGHT	0x0002F08B
	Observed Value:	
	NOM_PULS_OXIM_PERF_REL_RIGHT	0xF08B
	Units:	
	NOM_DIM_DIMLESS	0x0200
Perf l	Relative Perfusion Left	
	Label:	

	NLS_NOM_PULS_OXIM_PERF_REL_LEFT	0x0002F08A
	Observed Value:	
	NOM_PULS_OXIM_PERF_REL_LEFT	0xF08A
	Units:	
	NOM_DIM_DIMLESS	0x0200
NBP	non-invasive blood pressure	
	Label:	
	NLS_NOM_PRESS_BLD_NONINV	0x00024A04
	Observed Value (from VueLink):	
	NOM_PRESS_BLD_NONINV	0x4A04
	Compound Observed Value:	
	NOM_PRESS_BLD_NONINV_SYS	0x4A05
	NOM_PRESS_BLD_NONINV_DIA	0x4A06
	NOM_PRESS_BLD_NONINV_MEAN	0x4A07
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse from NBP	
	Label:	
	NLS_NOM_PRESS_BLD_NONINV_PULS_RATE	0x0002F0E5
	Observed Value:	
	NOM_PRESS_BLD_NONINV_PULS_RATE	0xF0E5
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
ABP	Arterial Blood Pressure (ABP)	
	Label:	
	NLS_NOM_PRESS_BLD_ART_ABP	0x00024A14
	Observed Value (from VueLink):	
	NOM_PRESS_BLD_ART_ABP	0x4A14
	Compound Observed Value:	
	NOM_PRESS_BLD_ART_ABP_SYS	0x4A15
	NOM_PRESS_BLD_ART_ABP_DIA	0x4A16
	NOM_PRESS_BLD_ART_ABP_MEAN	0x4A17
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from ABP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_ABP	0x80035402
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
ART	Arterial Blood Pressure (ART)	
	Label:	
	NLS_NOM_PRESS_BLD_ART	0x00024A10
	Observed Value (from VueLink):	
	NOM_PRESS_BLD_ART	0x4A10
	Compound Observed Value:	
	NOM_PRESS_BLD_ART_SYS	0x4A11
	NOM_PRESS_BLD_ART_DIA	0x4A12
	NOM_PRESS_BLD_ART_MEAN	0x4A13
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from ART	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_ART	0x80035403
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
Ao	Arterial Blood Pressure in the Aorta (Ao)	
	Label:	
	NLS_NOM_PRESS_BLD_AORT	0x00024A0C

	Observed Value (from VueLink):	
	NOM_PRESS_BLD_AORT	0x4A0C
	Compound Observed Value:	
	NOM_PRESS_BLD_AORT_SYS	0x4A0D
	NOM_PRESS_BLD_AORT_DIA	0x4A0E
	NOM_PRESS_BLD_AORT_MEAN	0x4A0F
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from Ao	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_AO	0x80035404
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
PAP	Pulmonary Arterial Pressure (PAP)	
	Label:	
	NLS_NOM_PRESS_BLD_ART_PULM	0x00024A1C
	Observed Value (from VueLink):	
	NOM_PRESS_BLD_ART_PULM	0x4A1C
	Compound Observed Value:	
	NOM_PRESS_BLD_ART_PULM_SYS	0x4A1D
	NOM_PRESS_BLD_ART_PULM_DIA	0x4A1E
	NOM_PRESS_BLD_ART_PULM_MEAN	0x4A1F
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from PAP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_PAP	0x80035405
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
CVP	Central Venous Pressure (CVP)	
	Label:	
	NLS_NOM_PRESS_BLD_VEN_CENT	0x00024A44
	Observed Value (from VueLink):	
	NOM_PRESS_BLD_VEN_CENT	0x4A44
	Compound Observed Value:	
	NOM_PRESS_BLD_VEN_CENT_SYS	0x4A45
	NOM_PRESS_BLD_VEN_CENT_DIA	0x4A46
	NOM_PRESS_BLD_VEN_CENT_MEAN	0x4A47
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from CVP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_CVP	0x80035406
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
RAP	Right Atrial Pressure (RAP)	
	Label:	
	NLS_NOM_PRESS_BLD_ATR_RIGHT	0x00024A34
	Observed Value (from VueLink):	
	NOM_PRESS_BLD_ATR_RIGHT	0x4A34
	Compound Observed Value:	
	NOM_PRESS_BLD_ATR_RIGHT_SYS	0x4A35
	NOM_PRESS_BLD_ATR_RIGHT_DIA	0x4A36
	NOM_PRESS_BLD_ATR_RIGHT_MEAN	0x4A37
	Units:	
	NOM_DIM_MMHG	0x0F20

Pulse	NOM_DIM_KILO_PASCAL	0x0F03
	Pulse derived from RAP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_RAP	0x80035407
	Observed Value:	
LAP	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Left Atrial Pressure (LAP)	
	Label:	
	NLS_NOM_PRESS_BLD_ATR_LEFT	0x00024A30
	Observed Value (from VueLink):	
	NOM_PRESS_BLD_ATR_LEFT	0x4A30
	Compound Observed Value:	
	NOM_PRESS_BLD_ATR_LEFT_SYS	0x4A31
	NOM_PRESS_BLD_ATR_LEFT_DIA	0x4A32
	NOM_PRESS_BLD_ATR_LEFT_MEAN	0x4A33
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from LAP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_LAP	0x80035408
	Observed Value:	
	NOM_PULS_RATE	0x480A
ICP	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Intra-cranial Pressure (ICP)	
	Label:	
	NLS_NOM_PRESS_INTRA_CRAN	0x00025808
	Observed Value (from VueLink):	
	NOM_PRESS_INTRA_CRAN	0x5808
	Compound Observed Value:	
	NOM_PRESS_INTRA_CRAN_SYS	0x5809
	NOM_PRESS_INTRA_CRAN_DIA	0x580A
	NOM_PRESS_INTRA_CRAN_MEAN	0x580B
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from ICP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_ICP	0x80035409
	Observed Value:	
	NOM_PULS_RATE	0x480A
UAP	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Umbilical Arterial Pressure (UAP)	
	Label:	
	NLS_NOM_PRESS_BLD_ART_UMB	0x00024A28
	Observed Value (from VueLink):	
	NOM_PRESS_BLD_ART_UMB	0x4A28
	Compound Observed Value:	
	NOM_PRESS_BLD_ART_UMB_SYS	0x4A29
	NOM_PRESS_BLD_ART_UMB_DIA	0x4A2A
	NOM_PRESS_BLD_ART_UMB_MEAN	0x4A2B
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from UAP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_UAP	0x8003540A
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0

UVP	Umbilical Venous Pressure (UVP)	
	Label:	
	NLS_NOM_PRESS_BLD_VEN_UMB	0x00024A48
	Observed Value (from VueLink):	
	NOM_PRESS_BLD_VEN_UMB	0x4A48
	Compound Observed Value:	
	NOM_PRESS_BLD_VEN_UMB_SYS	0x4A49
	NOM_PRESS_BLD_VEN_UMB_DIA	0x4A4A
	NOM_PRESS_BLD_VEN_UMB_MEAN	0x4A4B
	Units:	
Pulse	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
	Pulse derived from UVP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_UVP	0x8003540B
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
FAP	Femoral Arterial Pressure (FAP)	
	Label:	
	NLS_NOM_PRESS_BLD_ART_FEMORAL	0x0002F0BC
	Compound Observed Value:	
	NOM_PRESS_BLD_ART_FEMORAL_SYS	0xF0BD
	NOM_PRESS_BLD_ART_FEMORAL_DIA	0xF0BE
	NOM_PRESS_BLD_ART_FEMORAL_MEAN	0xF0BF
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from FAP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_FAP	0x80035434
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
BAP	Brachial Arterial Blood Pressure (BAP)	
	Label:	
	NLS_NOM_PRESS_BLD_ART_BRACHIAL	0x0002F0C0
	Compound Observed Value:	
	NOM_PRESS_BLD_ART_BRACHIAL_SYS	0xF0C1
	NOM_PRESS_BLD_ART_BRACHIAL_DIA	0xF0C2
	NOM_PRESS_BLD_ART_BRACHIAL_MEAN	0xF0C3
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from BAP	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_BAP	0x80035437
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
IC1	Intracranial Pressure 1 (IC1)	
	Label:	
	NLS_NOM_PRESS_INTRA_CRAN_1	0x0002F0B4
	Compound Observed Value:	
	NOM_PRESS_INTRA_CRAN_1_DIA	0xF0B6
	NOM_PRESS_INTRA_CRAN_1_SYS	0xF0B5
	NOM_PRESS_INTRA_CRAN_1_MEAN	0xF0B7
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from IC1	
	Label:	

IC2	NLS_PRESS_NAMES_PULSE_FROM_IC1	0x8003542E
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Intracranial Pressure 2 (IC2)	
	Label:	
	NLS_NOM_PRESS_INTRA_CRAN_2	0x0002F0B8
	Compound Observed Value:	
	NOM_PRESS_INTRA_CRAN_2_SYS	0xF0B9
Pulse	NOM_PRESS_INTRA_CRAN_2_DIA	0xF0BA
	NOM_PRESS_INTRA_CRAN_2_MEAN	0xF0BB
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
	Pulse derived from IC2	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_IC2	0x80035431
	Observed Value:	
	NOM_PULS_RATE	0x480A
P	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
	unspecific pressure	
	Label:	
	NLS_NOM_PRESS_BLD	0x00024A00
	Observed Value (from VueLink):	
	NOM_PRESS_BLD	0x4A00
	Compound Observed Value:	
	NOM_PRESS_BLD_SYS	0x4A01
	NOM_PRESS_BLD_DIA	0x4A02
Pulse	NOM_PRESS_BLD_MEAN	0x4A03
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
	Pulse derived from unspecific Pressure	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_P	0x80035401
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
P1	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Generic Pressure 1 (P1)	
	Label:	
	NLS_NOM_PRESS_GEN_1	0x0002F0A4
	Observed Value (from VueLink):	
	NOM_PRESS_GEN_1	0xF0A4
	Compound Observed Value:	
	NOM_PRESS_GEN_1_SYS	0xF0A5
	NOM_PRESS_GEN_1_DIA	0xF0A6
	NOM_PRESS_GEN_1_MEAN	0xF0A7
Pulse	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
	Pulse derived from P1	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_P1	0x80035422
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
	NOM_DIM_BEAT_PER_MIN	0x0AA0
P2	Generic Pressure 2 (P2)	
	Label:	
	NLS_NOM_PRESS_GEN_2	0x0002F0A8
	Observed Value (from VueLink):	
	NOM_PRESS_GEN_2	0xF0A8

	Compound Observed Value:	
	NOM_PRESS_GEN_2_SYS	0xF0A9
	NOM_PRESS_GEN_2_DIA	0xF0AA
	NOM_PRESS_GEN_2_MEAN	0xF0AB
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from P2	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_P2	0x80035425
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
P3	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Generic Pressure 3 (P3)	
	Label:	
	NLS_NOM_PRESS_GEN_3	0x0002F0AC
	Observed Value (from VueLink):	
	NOM_PRESS_GEN_3	0xF0AC
	Compound Observed Value:	
	NOM_PRESS_GEN_3_SYS	0xF0AD
	NOM_PRESS_GEN_3	0xF0AC
	NOM_PRESS_GEN_3_MEAN	0xF0AF
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from P3	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_P3	0x80035428
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
P4	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Generic Pressure 4 (P4)	
	Label:	
	NLS_NOM_PRESS_GEN_4	0x0002F0B0
	Observed Value (from VueLink):	
	NOM_PRESS_GEN_4	0xF0B0
	Compound Observed Value:	
	NOM_PRESS_GEN_4_SYS	0xF0B1
	NOM_PRESS_GEN_4_DIA	0xF0B2
	NOM_PRESS_GEN_4_MEAN	0xF0B3
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
Pulse	Pulse derived from P4	
	Label:	
	NLS_PRESS_NAMES_PULSE_FROM_P4	0x8003542B
	Observed Value:	
	NOM_PULS_RATE	0x480A
	Units:	
IUP	NOM_DIM_BEAT_PER_MIN	0x0AA0
	Intra-Uterine Pressure	
	Label:	
	NLS_NOM_PRESS_INTRA_UTERIAL	0x0002F0D8
	Observed Value:	
	NOM_PRESS_BLD	0x4A00
PAWP	Pulmonary Artery Wedge Pressure	
	Label:	
	NLS_NOM_PRESS_BLD_ART_PULM_WEDGE	0x00024A24
	Observed Value:	
	NOM_PRESS_BLD_ART_PULM_WEDGE	0x4A24
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03

CPP	Cerebral Perfusion Pressure	
	Label:	
	NLS_NOM_PRESS_CEREB_PERF	0x00025804
	Observed Value:	
	NOM_PRESS_CEREB_PERF	0x5804
PPV	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
	Pulse Pressure Variation	
	Label:	
CCO	NLS_NOM_PULS_PRESS_VAR	0x0002F0E3
	Observed Value:	
	NOM_PULS_PRESS_VAR	0xF0E3
	Continuous Cardiac Output	
	Label:	
CCI	NLS_NOM_OUTPUT_CARD_CTS	0x00024BDC
	Observed Value:	
	NOM_OUTPUT_CARD_CTS	0x4BDC
	Units:	
	NOM_DIM_X_L_PER_MIN	0x0C00
SV	Continuous Cardiac Output Index	
	Label:	
	NLS_NOM_OUTPUT_CARD_INDEX_CTS	0x0002F047
	Observed Value:	
	NOM_OUTPUT_CARD_INDEX_CTS	0xF047
SV	Units:	
	NOM_DIM_X_L_PER_MIN_PER_M_SQ	0x0B20
	Stroke Volume	
	Label:	
	NLS_NOM_VOL_BLD_STROKE	0x00024B84
SI	Observed Value:	
	NOM_VOL_BLD_STROKE	0x4B84
	Units:	
	NOM_DIM_MILLI_L	0x0652
	Stroke Index	
SVV	Label:	
	NLS_NOM_VOL_BLD_STROKE_INDEX	0x0002F048
	Observed Value:	
	NOM_VOL_BLD_STROKE_INDEX	0xF048
	Units:	
dPmax	NOM_DIM_MILLI_L_PER_M_SQ	0x0592
	Stroke Volume Variation	
	Label:	
	NLS_NOM_VOL_BLD_STROKE_VAR	0x0002F049
	Observed Value:	
C.O.	NOM_VOL_BLD_STROKE_VAR	0xF049
	Units:	
	NOM_DIM_PERCENT	0x0220
	Index of Left Ventricular Contractility	
	Label:	
C.I.	NLS_NOM_GRAD_PRESS_BLD_AORT_POS_MAX	0x00024C25
	Observed Value:	
	NOM_GRAD_PRESS_BLD_AORT_POS_MAX	0x4C25
	Cardiac Output	
	Label:	
	NLS_NOM_OUTPUT_CARD	0x00024B04
	Observed Value:	
	NOM_OUTPUT_CARD	0x4B04
	Units:	
	NOM_DIM_X_L_PER_MIN	0x0C00
	Cardiac Index	
	Label:	
	NLS_NOM_OUTPUT_CARD_INDEX	0x0002490C
	Observed Value:	
	NOM_OUTPUT_CARD_INDEX	0x490C

ITBV	Units:	
	NOM_DIM_X_L_PER_MIN_PER_M_SQ	0x0B20
	Intrathoracic Blood Volume	
	Label:	
ITBVI	NLS_NOM_VOL_BLD_INTRA_THOR	0x0002F040
	Observed Value:	
	NOM_VOL_BLD_INTRA_THOR	0xF040
	Units:	
EVLW	NOM_DIM_MILLI_L	0x0652
	Intrathoracic Blood Volume Index	
	Label:	
	NLS_NOM_VOL_BLD_INTRA_THOR_INDEX	0x0002F041
EVLWI	Observed Value:	
	NOM_VOL_BLD_INTRA_THOR_INDEX	0xF041
	Units:	
	NOM_DIM_MILLI_L_PER_M_SQ	0x0592
GEDV	Extravascular Lung Water	
	Label:	
	NLS_NOM_VOL_LUNG_WATER_EXTRA_VASC	0x0002F042
	Observed Value:	
GEDVI	NOM_VOL_LUNG_WATER_EXTRA_VASC	0xF042
	Units:	
	NOM_DIM_MILLI_L	0x0652
	Extravascular Lung Water Index	
CFI	Label:	
	NLS_NOM_VOL_LUNG_WATER_EXTRA_VASC_INDEX	0x0002F043
	Observed Value:	
	NOM_VOL_LUNG_WATER_EXTRA_VASC_INDEX	0xF043
PVPI	Units:	
	NOM_DIM_MILLI_L_PER_KG	0x0C72
	Global End Diastolic Volume	
	Label:	
GEF	NLS_NOM_VOL_GLOBAL_END_DIA	0x0002F044
	Observed Value:	
	NOM_VOL_GLOBAL_END_DIA	0xF044
	Units:	
SNR	NOM_DIM_MILLI_L	0x0652
	Global End Diastolic Volume Index	
	Label:	
	NLS_NOM_VOL_GLOBAL_END_DIA_INDEX	0x0002F045
PVPI	Observed Value:	
	NOM_VOL_GLOBAL_END_DIA_INDEX	0xF045
	Units:	
	NOM_DIM_MILLI_L_PER_M_SQ	0x0592
CFI	Cardiac Function Index	
	Label:	
	NLS_NOM_CARD_FUNC_INDEX	0x0002F046
	Observed Value:	
PVPI	NOM_CARD_FUNC_INDEX	0xF046
	Units:	
	NOM_DIM_DIMLESS	0x0200
	Pulmonary Vascular Permeability Index	
GEF	Label:	
	NLS_NOM_PERM_VASC_PULM_INDEX	0x0002F106
	Observed Value:	
	NOM_PERM_VASC_PULM_INDEX	0xF106
SNR	Global Ejection Fraction	
	Label:	
	NLS_NOM_FRACT_EJECT	0x0002F105
	Observed Value:	
SNR	NOM_FRACT_EJECT	0xF105
	Signal to Noise ratio	
	Label:	
	NLS_NOM_SNR	0x0002F101
	Observed Value:	

RLShnt	NOM_SNR	0xF101
	Right-to-Left Heart Shunt	
	Label:	
	NLS_NOM_SHUNT_RIGHT_LEFT	0x0002F14A
	Observed Value:	
SaO2	NOM_SHUNT_RIGHT_LEFT	0xF14A
	Units:	
	NOM_DIM_MILLI_SECOND	
	Oxygen Saturation	
	Label:	
SvO2	NLS_NOM_SAT_O2_ART	0x00024B34
	Observed Value:	
	NOM_SAT_O2_ART	0x4B34
	Mixed Venous Oxygen Saturation	
	Label:	
ScvO2	NLS_NOM_SAT_O2_VEN	0x00024B3C
	Observed Value:	
	NOM_SAT_O2_VEN	0x4B3C
	Units:	
	NOM_DIM_PERCENT	0x0220
SO2	Central Venous Oxygen Saturation	
	Label:	
	NLS_NOM_SAT_O2_VEN_CENT	0x0002F100
	Observed Value:	
	NOM_SAT_O2_VEN_CENT	0xF100
SO2 1	O2 Saturation	
	Label:	
	NLS_NOM_SAT_O2	0x00024B2C
	Observed Value:	
	NOM_SAT_O2	0x4B2C
SO2 r	Units:	
	NOM_DIM_PERCENT	
	Oxygen Saturation Left Side	
	Label:	
	NLS_NOM_SAT_O2_LEFT	0x0002F89D
SO2 1	Observed Value:	
	NOM_SAT_O2_ART	0x4B34
	Oxygen Saturation Right Side	
	Label:	
	NLS_NOM_SAT_O2_RIGHT	0x0002F89E
SO2 2	Observed Value:	
	NOM_SAT_O2_ART	0x4B34
	O2 Saturation 1 (generic)	
	Label:	
	NLS_NOM_SAT_O2_GEN_1	0x0002F962
SO2 3	Observed Value:	
	NOM_SAT_O2_GEN_1	0xF962
	Units:	
	NOM_DIM_PERCENT	0x0220
	O2 Saturation 2 (generic)	
SO2 4	Label:	
	NLS_NOM_SAT_O2_GEN_2	0x0002F963
	Observed Value:	
	NOM_SAT_O2_GEN_2	0xF963
	Units:	
SO2 5	NOM_DIM_PERCENT	0x0220
	O2 Saturation 3 (generic)	
	Label:	
	NLS_NOM_SAT_O2_GEN_3	0x0002F964
	Observed Value:	
SO2 6	NOM_SAT_O2_GEN_3	0xF964
	Units:	
	NOM_DIM_PERCENT	0x0220
	O2 Saturation 4 (generic)	
	Label:	

	NLS_NOM_SAT_O2_GEN_4	0x0002F965
	Observed Value:	
	NOM_SAT_O2_GEN_4	0xF965
	Units:	
	NOM_DIM_PERCENT	0x0220
LI	Light Intenisty. SvO2	
	Label:	
	NLS_NOM_INTENS_LIGHT	0x0002F072
	Observed Value:	
	NOM_INTENS_LIGHT	0xF072
DO2	Oxygen Availability DO2	
	Label:	
	NLS_NOM_SAT_O2_DELIVER	0x0002F06D
	Observed Value:	
	NOM_SAT_O2_DELIVER	0xF06D
DO2I	Oxygen Availability Index	
	Label:	
	NLS_NOM_SAT_O2_DELIVER_INDEX	0x0002F06E
	Observed Value:	
	NOM_SAT_O2_DELIVER_INDEX	0xF06E
O2ER	Oxygen Extraction Ratio	
	Label:	
	NLS_NOM_RATIO_SAT_O2_CONSUMP_DELIVER	0x0002F06F
	Observed Value:	
	NOM_RATIO_SAT_O2_CONSUMP_DELIVER	0xF06F
Qs/Qt	Percent Alveolarvenous Shunt Qs/Qt	
	Label:	
	NLS_NOM_RATIO_ART_VEN_SHUNT	0x0002F070
	Observed Value:	
	NOM_RATIO_ART_VEN_SHUNT	0xF070
AaDO2	Alveolar- Arterial Oxygen Difference	
	Label:	
	NLS_NOM_SAT_DIFF_O2_ART_ALV	0x00024B40
	Observed Value:	
	NOM_SAT_DIFF_O2_ART_ALV	0x4B40
Sp-vO2	Difference between Spo2 and SvO2	
	Label:	
	NLS_NOM_SAT_DIFF_O2_ART_VEN	0x0002F06C
	Observed Value:	
	NOM_SAT_DIFF_O2_ART_VEN	0xF06C
tcGas	Generic Term for the Transcutaneous Gases	
	Label:	
	NLS_NOM_GAS_TCUT	0x0002F051
	Observed Value:	
	NOM_GAS_TCUT	0xF051
tcpO2	Transcutaneous Oxygen Partial Pressure	
	Label:	
	NLS_NOM_O2_TCUT	0x000250D0
	Observed Value:	
	NOM_O2_TCUT	0x50D0
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
tcpCO2	Transcutaneous Carbon Dioxide Partial Pressure	
	Label:	
	NLS_NOM_CO2_TCUT	0x000250CC
	Observed Value:	
	NOM_CO2_TCUT	0x50CC
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_KILO_PASCAL	0x0F03
SitTim	NOM_DIM_MIN	
	Label:	
	NLS_NOM_TIME_TCUT_SENSOR	0x0002F03E
	Observed Value:	

	NOM_TIME_TCUT_SENSOR	0xF03E
SensrT	Sensor Temperature	
	Label:	
	NLS_NOM_TEMP_TCUT_SENSOR	0x0002F03F
	Observed Value:	
	NOM_TEMP_TCUT_SENSOR	0xF03F
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
HeatPw	NOM_DIM_MILLI_WATT	
	Label:	
	NLS_NOM_HEATING_PWR_TCUT_SENSOR	0x0002F076
	Observed Value:	
	NOM_HEATING_PWR_TCUT_SENSOR	0xF076
CO2	CO2 concentration	
	Label:	
	NLS_NOM_AWAY_CO2	0x000250AC
	Observed Value (from VueLink):	
	NOM_AWAY_CO2	0x50AC
	Compound Observed Value:	
	NOM_AWAY_CO2_ET	0x50B0
	NOM_AWAY_CO2_INSP_MIN	0x50BA
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_PERCENT	0x0220
	NOM_DIM_KILO_PASCAL	0x0F03
RRspir	Respiration Rate from Spirometry	
	Label:	
	NLS_NOM_AWAY_RESP_RATE_SPIRO	0x0002F0E2
	Observed Value:	
	NOM_AWAY_RESP_RATE_SPIRO	0xF0E2
	Units:	
	NOM_DIM_RESP_PER_MIN	0x0AE0
awRR	Airway Respiration Rate	
	Label:	
	NLS_NOM_AWAY_RESP_RATE	0x00025012
	Observed Value:	
	NOM_AWAY_RESP_RATE	0x5012
	Units:	
	NOM_DIM_RESP_PER_MIN	0x0AE0
O2	Generic oxygen measurement label	
	Label:	
	NLS_NOM_CONC_AWAY_O2	0x00025164
	Observed Value (from VueLink):	
	NOM_CONC_AWAY_O2	0x5164
	Compound Observed Value:	
	NOM_CONC_AWAY_O2_ET	0x5378
	NOM_CONC_AWAY_O2_INSP	0x5284
	Units:	
	NOM_DIM_MMHG	0x0F20
	NOM_DIM_PERCENT	0x0220
	NOM_DIM_KILO_PASCAL	0x0F03
FIO2	Fractional Inspired Oxygen FIO2	
	Label:	
	NLS_NOM_VENT_CONC_AWAY_O2_INSP	0x00027498
	Observed Value:	
	NOM_VENT_CONC_AWAY_O2_INSP	0x7498
RR	Respiration Rate	
	Label:	
	NLS_NOM_RESP_RATE	0x0002500A
	Observed Value:	
	NOM_RESP_RATE	0x500A
	Units:	
	NOM_DIM_RESP_PER_MIN	0x0AE0
T.I.	Transthoracic Impedance	

VCO2	Label:	
	NLS_NOM_IMPED_TTHOR	0x000250E4
	Observed Value:	
VCO2ti	NOM_IMPED_TTHOR	0x50E4
	CO2 Production	
	Label:	
VCO2ti	NLS_NOM_FLOW_CO2_PROD_RESP	0x000250E0
	Observed Value:	
	NOM_FLOW_CO2_PROD_RESP	0x50E0
Pplat	CO2 Tidal Production	
	Label:	
	NLS_NOM_FLOW_CO2_PROD_RESP_TIDAL	0x0002F882
Pplat	Observed Value:	
	NOM_FLOW_CO2_PROD_RESP_TIDAL	0xF882
	Plateau Pressure	
AWP	Label:	
	NLS_NOM_PRESS_RESP_PLAT	0x000250E8
	Observed Value:	
AWP	NOM_PRESS_RESP_PLAT	0x50E8
	Airway Pressure Wave	
	Label:	
AWPmin	NLS_NOM_PRESS_AWAY	0x000250F0
	Observed Value:	
	NOM_PRESS_AWAY	0x50F0
CPAP	Airway Pressure Minimum	
	Label:	
	NLS_NOM_PRESS_AWAY_MIN	0x000250F2
CPAP	Observed Value:	
	NOM_PRESS_AWAY_MIN	0x50F2
	Continuous Positive Airway Pressure	
iPEEP	Label:	
	NLS_NOM_PRESS_AWAY_CTS_POS	0x000250F4
	Observed Value:	
iPEEP	NOM_PRESS_AWAY_CTS_POS	0x50F4
	Intrinsic PEEP Breathing Pressure	
	Label:	
AWPin	NLS_NOM_PRESS_AWAY_END_EXP_POS_INTRINSIC	0x00025100
	Observed Value:	
	NOM_PRESS_AWAY_END_EXP_POS_INTRINSIC	0x5100
AWPin	Airway Pressure Wave - measured in the inspiratory path	
	Label:	
	NLS_NOM_PRESS_AWAY_INSP	0x00025108
PIP	Observed Value:	
	NOM_PRESS_AWAY_INSP	0x5108
	Positive Inspiratory ressure	
MnAwP	Label:	
	NLS_NOM_PRESS_AWAY_INSP_MAX	0x00025109
	Observed Value:	
MnAwP	NOM_PRESS_AWAY_INSP_MAX	0x5109
	Mean Airway Pressure. Printer Context	
	Label:	
I:E 1:	NLS_NOM_PRESS_AWAY_INSP_MEAN	0x0002510B
	Observed Value:	
	NOM_PRESS_AWAY_INSP_MEAN	0x510B
Vd/Vt	Inpired:Expired Ratio	
	Label:	
	NLS_NOM_RATIO_IE	0x00025118
Vd/Vt	Observed Value:	
	NOM_RATIO_IE	0x5118
	Ratio of Deadspace to Tidal Volume Vd/Vt	
Raw	Label:	
	NLS_NOM_RATIO_AWAY_DEADSP_TIDAL	0x0002511C
	Observed Value:	
Raw	NOM_RATIO_AWAY_DEADSP_TIDAL	0x511C
	Static Lung Resistance	

TV	Label:	
	NLS_NOM_RES_AWAY	0x00025120
	Observed Value:	
TVexp	NOM_RES_AWAY	0x5120
	Tidal Volume	
	Label:	
TVexp	NLS_NOM_VOL_AWAY_TIDAL	0x0002513C
	Observed Value (from VueLink):	
	NOM_VOL_AWAY_TIDAL	0x513C
TVin	Compound Observed Value:	
	expired Tidal Volume	
	Label:	
TVin	NLS_NOM_VOL_AWAY_EXP_TIDAL	0x0002F0E1
	Observed Value:	
	NOM_VOL_AWAY_EXP_TIDAL	0xF0E1
MINVOL	inspired Tidal Volume	
	Label:	
	NLS_NOM_VOL_AWAY_INSP_TIDAL	0x0002F0E0
MINVOL	Observed Value:	
	NOM_VOL_AWAY_INSP_TIDAL	0xF0E0
	Airway Minute Volum Inspiratory	
PlatTi	Label:	
	NLS_NOM_VOL_MINUTE_AWAY	0x00025148
	Observed Value (from VueLink):	
PlatTi	NOM_VOL_MINUTE_AWAY	0x5148
	Compound Observed Value:	
	NOM_VOL_MINUTE_AWAY_EXP	0x514C
PlatTi	NOM_VOL_MINUTE_AWAY_INSP	0x5150
	Units:	
	NOM_DIM_X_L_PER_MIN	0x0C00
SpMV	Plateau Time	
	Label:	
	NLS_NOM_TIME_PD_RESP_PLAT	0x0002F0FF
SpMV	Observed Value:	
	NOM_TIME_PD_RESP_PLAT	0xF0FF
	Spontaneous Minute Volume	
DeltaO2	Label:	
	NLS_NOM_VENT_VOL_MINUTE_AWAY_SPONT	0x0002F091
	Observed Value:	
DeltaO2	NOM_VENT_VOL_MINUTE_AWAY_SPONT	0xF091
	relative Dead Space	
	Label:	
PECO2	NLS_NOM_VENT_CONC_AWAY_O2_DELTA	0x00025168
	Observed Value:	
	NOM_VENT_CONC_AWAY_O2_DELTA	0x5168
PECO2	Partial O2 Venous	
	Label:	
	NLS_NOM_VENT_AWAY_CO2_EXP	0x0002517C
AWFin	Observed Value:	
	NOM_VENT_AWAY_CO2_EXP	0x517C
	Airway Flow Wave - measured in the inspiratory path	
VQI	Label:	
	NLS_NOM_VENT_FLOW_INSP	0x0002518C
	Observed Value:	
VQI	NOM_VENT_FLOW_INSP	0x518C
	Ventilation Perfusion Index	
	Label:	
Poccl	NLS_NOM_VENT_FLOW_RATIO_PERF_ALV_INDEX	0x00025190
	Observed Value:	
	NOM_VENT_FLOW_RATIO_PERF_ALV_INDEX	0x5190
Poccl	Occlusion Pressure	
	Label:	
	NLS_NOM_VENT_PRESS_OCCL	0x0002519C
Poccl	Observed Value:	
	NOM_VENT_PRESS_OCCL	0x519C

PEEP	Positive End-Expiratory Pressure PEEP	
	Label:	
	NLS_NOM_VENT_PRESS_AWAY_END_EXP_POS	0x000251A8
	Observed Value:	
	NOM_VENT_PRESS_AWAY_END_EXP_POS	0x51A8
Vd	Dead Space Volume Vd	
	Label:	
	NLS_NOM_VENT_VOL_AWAY_DEADSP	0x000251B0
	Observed Value:	
	NOM_VENT_VOL_AWAY_DEADSP	0x51B0
relVd	relative Dead Space	
	Label:	
	NLS_NOM_VENT_VOL_AWAY_DEADSP_REL	0x000251B4
	Observed Value:	
	NOM_VENT_VOL_AWAY_DEADSP_REL	0x51B4
TrpVol	Lung Volume Trapped	
	Label:	
	NLS_NOM_VENT_VOL_LUNG_TRAPD	0x000251B8
	Observed Value:	
	NOM_VENT_VOL_LUNG_TRAPD	0x51B8
Leak	Leakage	
	Label:	
	NLS_NOM_VENT_VOL_LEAK	0x00025370
	Observed Value:	
	NOM_VENT_VOL_LEAK	0x5370
ALVENT	Alveolar Ventilation ALVENT	
	Label:	
	NLS_NOM_VENT_VOL_LUNG_ALV	0x00025374
	Observed Value:	
	NOM_VENT_VOL_LUNG_ALV	0x5374
VC	Vital Lung Capacity	
	Label:	
	NLS_NOM_CAPAC_VITAL	0x00025080
	Observed Value:	
	NOM_CAPAC_VITAL	0x5080
COMP	generic label Lung Compliance	
	Label:	
	NLS_NOM_COMPL_LUNG	0x00025088
	Observed Value:	
	NOM_COMPL_LUNG	0x5088
Cdyn	Dynamic Lung Compliance	
	Label:	
	NLS_NOM_COMPL_LUNG_DYN	0x0002508C
	Observed Value:	
	NOM_COMPL_LUNG_DYN	0x508C
Cstat	Static Lung Compliance	
	Label:	
	NLS_NOM_COMPL_LUNG_STATIC	0x00025090
	Observed Value:	
	NOM_COMPL_LUNG_STATIC	0x5090
PIF	Inspiratory Peak Flow	
	Label:	
	NLS_NOM_FLOW_AWAY_INSP_MAX	0x000250DD
	Observed Value:	
	NOM_FLOW_AWAY_INSP_MAX	0x50DD
PEF	Expiratory Peak Flow	
	Label:	
	NLS_NOM_FLOW_AWAY_EXP_MAX	0x000250D9
	Observed Value:	
	NOM_FLOW_AWAY_EXP_MAX	0x50D9
BIS	Bispectral Index	
	Label:	
	NLS_NOM_EEG_BISPECTRAL_INDEX	0x0002F04E
	Observed Value:	
	NOM_EEG_BISPECTRAL_INDEX	0xF04E

SQI	Units:	
	NOM_DIM_DIMLESS	0x0200
	Signal Quality Index	
	Label:	
EMG	NLS_NOM_EEG_BIS_SIG_QUAL_INDEX	0x0002F04D
	Observed Value:	
	NOM_EEG_BIS_SIG_QUAL_INDEX	0xF04D
	Units:	
TP	NOM_DIM_PERCENT	0x0220
	Electromyography	
	Label:	
	NLS_NOM_EMG_ELEC_POTL_MUSCL	0x0002593C
TP1	Observed Value:	
	NOM_EMG_ELEC_POTL_MUSCL	0x593C
	Units:	
	NOM_DIM_DECIBEL	0x1920
TP2	Total Power	
	Label:	
	NLS_NOM_EEG_PWR_SPEC_TOT	0x000259B8
	Observed Value:	
SR	NOM_EEG_PWR_SPEC_TOT	0x59B8
	Units:	
	NOM_DIM_DECIBEL	0x1920
	Total Power channel 1	
SEF	Label:	
	NLS_EEG_NAMES_CHAN_TP1	0x800F5403
	Observed Value:	
	NOM_EEG_PWR_SPEC_TOT	0x59B8
MDF	Units:	
	NOM_DIM_NANO_WATT	0x0FD4
	Total Power channel 2	
	Label:	
PPF	NLS_EEG_NAMES_CHAN_TP2	0x800F5404
	Observed Value:	
	NOM_EEG_PWR_SPEC_TOT	0x59B8
	Units:	
SEF	NOM_DIM_NANO_WATT	0x0FD4
	Suppression Ratio	
	Label:	
	NLS_NOM_EEG_RATIO_SUPPRN	0x0002F04A
MDF	Observed Value:	
	NOM_EEG_RATIO_SUPPRN	0xF04A
	Units:	
	NOM_DIM_PERCENT	0x0220
PPF	Spectral Edge Frequency	
	Label:	
	NLS_NOM_EEG_FREQ_PWR_SPEC_CRTX_SPECTRAL_EDGE	0x00025988
	Observed Value:	
MDF	NOM_EEG_FREQ_PWR_SPEC_CRTX_SPECTRAL_EDGE	0x5988
	Units:	
	NOM_DIM_HZ	0x09C0
	Mean Dominant Frequency	
PPF	Label:	
	NLS_NOM_EEG_FREQ_PWR_SPEC_CRTX_DOM_MEAN	0x0002597C
	Observed Value:	
	NOM_EEG_FREQ_PWR_SPEC_CRTX_DOM_MEAN	0x597C
PPF	Units:	
	NOM_DIM_HZ	0x09C0
	Peak Power Frequency	
	Label:	
PPF	NLS_NOM_EEG_FREQ_PWR_SPEC_CRTX_PEAK	0x00025984
	Observed Value:	
	NOM_EEG_FREQ_PWR_SPEC_CRTX_PEAK	0x5984
	Units:	
PPF	NOM_DIM_HZ	0x09C0

Frequ1	generic label for EEG channel 1	
	Label:	
	NLS_EEG_NAMES_CHAN_FREQ1	0x800F5413
	Compound Observed Value:	
	NOM_EEG_FREQ_PWR_SPEC_CRTX_SPECTRAL_EDGE	0x5988
	NOM_EEG_FREQ_PWR_SPEC_CRTX_DOM_MEAN	0x597C
	NOM_EEG_FREQ_PWR_SPEC_CRTX_PEAK	0x5984
	Units:	
	NOM_DIM_HZ	0x09C0
Frequ2	generic label for EEG channel 2	
	Label:	
	NLS_EEG_NAMES_CHAN_FREQ2	0x800F5414
	Compound Observed Value:	
	NOM_EEG_FREQ_PWR_SPEC_CRTX_SPECTRAL_EDGE	0x5988
	NOM_EEG_FREQ_PWR_SPEC_CRTX_DOM_MEAN	0x597C
	NOM_EEG_FREQ_PWR_SPEC_CRTX_PEAK	0x5984
	Units:	
	NOM_DIM_HZ	0x09C0
Prcnt1	generic label for EEG channel 1	
	Label:	
	NLS_EEG_NAMES_CHAN_PCNT1	0x800F5415
	Compound Observed Value:	
	NOM_EEG_PWR_SPEC_ALPHA_REL	0x59D4
	NOM_EEG_PWR_SPEC_BETA_REL	0x59D8
	NOM_EEG_PWR_SPEC_DELTA_REL	0x59DC
	NOM_EEG_PWR_SPEC_THETA_REL	0x59E0
	Units:	
	NOM_DIM_PERCENT	0x0220
Prcnt2	generic label for EEG channel 2	
	Label:	
	NLS_EEG_NAMES_CHAN_PCNT2	0x800F5416
	Compound Observed Value:	
	NOM_EEG_PWR_SPEC_ALPHA_REL	0x59D4
	NOM_EEG_PWR_SPEC_BETA_REL	0x59D8
	NOM_EEG_PWR_SPEC_DELTA_REL	0x59DC
	NOM_EEG_PWR_SPEC_THETA_REL	0x59E0
	Units:	
	NOM_DIM_PERCENT	0x0220
AAI	A-Line ARX Index	
	Label:	
	NLS_NOM_ELEC_EVOK_POTL_CRTX_ACOUSTIC_AAI	0x0002F873
BSI	Observed Value:	
	NOM_ELEC_EVOK_POTL_CRTX_ACOUSTIC_AAI	0xF873
	Burst Suppression Indicator	
	Label:	
	NLS_NOM_EEG_BURST_SUPPRN_INDEX	0x0002F840
	Observed Value:	
Temp	NOM_EEG_BURST_SUPPRN_INDEX	0xF840
	Unspecific Temperature	
	Label:	
	NLS_NOM_TEMP	0x00024B48
	Observed Value:	
	NOM_TEMP	0x4B48
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
Trect	Rectal Temperature	
	Label:	
	NLS_NOM_TEMP_RECT	0x0002E004
	Observed Value:	
	NOM_TEMP_RECT	0xE004
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
Tblood	Tblood	

Tcore	Label:	
	NLS_NOM_TEMP_BLD	0x0002E014
	Observed Value:	
	NOM_TEMP_BLD	0xE014
	Units:	
Tskin	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Core (Body) Temperature	
	Label:	
	NLS_NOM_TEMP_CORE	0x00024B60
Tesoph	Observed Value:	
	NOM_TEMP_CORE	0x4B60
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
Tnasoph	Skin Temperature	
	Label:	
	NLS_NOM_TEMP_SKIN	0x00024B74
	Observed Value:	
	NOM_TEMP_SKIN	0x4B74
Tnasoph	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Esophagial Temperature	
	Label:	
Tnasoph	NLS_NOM_TEMP_ESOPH	0x00024B64
	Observed Value:	
	NOM_TEMP_ESOPH	0x4B64
	Units:	
	NOM_DIM_DEGC	0x17A0
Tnasoph	NOM_DIM_FAHR	0x1140
	Naso pharyngial Temperature	
	Label:	
	NLS_NOM_TEMP_NASOPH	0x00024B6C
	Observed Value:	
Tnasoph	NOM_TEMP_NASOPH	0x4B6C
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
Tnasoph	Areterial Temperature	
	Label:	
	NLS_NOM_TEMP_ART	0x00024B50
	Observed Value:	
	NOM_TEMP_ART	0x4B50
Tnasoph	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Venous Temperature	
	Label:	
Tnasoph	NLS_NOM_TEMP_VEN	0x00024B7C
	Observed Value:	
	NOM_TEMP_VEN	0x4B7C
	Units:	
	NOM_DIM_DEGC	0x17A0
Tnasoph	NOM_DIM_FAHR	0x1140
	Temperature of the Urine fluid	
	Label:	
	NLS_NOM_TEMP_VESICAL	0x0002F0C4
	Observed Value:	
Tnasoph	NOM_TEMP_VESICAL	0xF0C4
	Tympanic Temperature	
	Label:	
	NLS_NOM_TEMP_TYMP	0x00024B78
	Observed Value:	
Tnasoph	NOM_TEMP_TYMP	0x4B78

Tcereb	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Cerebral Temperature	
	Label:	
Tamb	NLS_NOM_TEMP_CEREBRAL	0x0002F0C5
	Observed Value:	
	NOM_TEMP_CEREBRAL	0xF0C5
	Units:	
	NOM_DIM_DEGC	0x17A0
Tairwy	NOM_DIM_FAHR	0x1140
	Ambient Temperature	
	Label:	
	NLS_NOM_TEMP_AMBIENT	0x0002F0C6
	Observed Value:	
Tinj	NOM_TEMP_AMBIENT	0xF0C6
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Airway Temperature	
T1Core	Label:	
	NLS_NOM_TEMP_AWAY	0x00024B54
	Observed Value:	
	NOM_TEMP_AWAY	0x4B54
	Units:	
T2Core	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Injectate Temperature	
	Label:	
	NLS_NOM_TEMP_INJ	0x00024B68
DeltaTemp	Observed Value:	
	NOM_TEMP_INJ	0x4B68
	Core Temperature 1 (generic)	
	Label:	
	NLS_NOM_TEMP_CORE_GEN_1	0x0002F966
Tbody	Observed Value:	
	NOM_TEMP_CORE_GEN_1	0xF966
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
pTrect	Core Temperature 2 (generic)	
	Label:	
	NLS_NOM_TEMP_CORE_GEN_2	0x0002F967
	Observed Value:	
	NOM_TEMP_CORE_GEN_2	0xF967
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Difference Temperature	
	Label:	
	NLS_NOM_TEMP_DIFF	0x0002E018
	Observed Value:	
	NOM_TEMP_DIFF	0xE018
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Patient Temperature	
	Label:	
	NLS_NOM_TEMP_BODY	0x00024B5C
	Observed Value:	
	NOM_TEMP	0x4B48
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Predictive Rectal Temperature	

	Label:	
	NLS_NOM_TEMP_RECT_PRED	0x0002F114
	Observed Value:	
	NOM_TEMP_RECT_PRED	0xF114
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Predictive Oral Temperature	
pToral	Label:	
	NLS_NOM_TEMP_ORAL_PRED	0x0002F110
	Observed Value:	
	NOM_TEMP_ORAL_PRED	0xF110
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Predictive Axillary Temperature	
pTaxil	Label:	
	NLS_NOM_TEMP_AXIL_PRED	0x0002F118
	Observed Value:	
	NOM_TEMP_AXIL_PRED	0xF118
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Generic Temperature 1 (T1)	
T1	Label:	
	NLS_NOM_TEMP_GEN_1	0x0002F0C7
	Observed Value:	
	NOM_TEMP_GEN_1	0xF0C7
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Generic Temperature 2 (T2)	
T2	Label:	
	NLS_NOM_TEMP_GEN_2	0x0002F0C8
	Observed Value:	
	NOM_TEMP_GEN_2	0xF0C8
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Generic Temperature 3 (T3)	
T3	Label:	
	NLS_NOM_TEMP_GEN_3	0x0002F0C9
	Observed Value:	
	NOM_TEMP_GEN_3	0xF0C9
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	Generic Temperature 4 (T4)	
T4	Label:	
	NLS_NOM_TEMP_GEN_4	0x0002F0CA
	Observed Value:	
	NOM_TEMP_GEN_4	0xF0CA
	Units:	
	NOM_DIM_DEGC	0x17A0
	NOM_DIM_FAHR	0x1140
	generic N2 label	
N2	Label:	
	NLS_NOM_CONC_AWAY_N2	0x0002537C
	Observed Value (from VueLink):	
	NOM_CONC_AWAY_N2	0x537C
	Compound Observed Value:	
	NOM_CONC_AWAY_N2_ET	0x5380
	NOM_CONC_AWAY_N2_INSP	0x5384
	Units:	
	NOM_DIM_MMHG	0x0F20

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Attribute IDs

The Attribute ID specifies the type of an attribute in the AttributeList. The IDs are taken from the Object Oriented Elements partition. Unknown attributes should be ignored.

Device P-Alarm List	NOM_ATTR_AL_MON_P_AL_LIST	0x0902
Device T-Alarm List	NOM_ATTR_AL_MON_T_AL_LIST	0x0904
Altitude	NOM_ATTR_ALTITUDE	0x090C
Application Area	NOM_ATTR_AREA_APPL	0x090D
Color	NOM_ATTR_COLOR	0x0911
Device Alert Condition	NOM_ATTR_DEV_AL_COND	0x0916
Display Resolution	NOM_ATTR_DISP_RES	0x0917
Visual Grid	NOM_ATTR_GRID_VIS_I16	0x091A
Association Invoke Id	NOM_ATTR_ID_ASSOC_NO	0x091D
Bed Label	NOM_ATTR_ID_BED_LABEL	0x091E
Object Handle	NOM_ATTR_ID_HANDLE	0x0921
Label	NOM_ATTR_ID_LABEL	0x0924
Label String	NOM_ATTR_ID_LABEL_STRING	0x0927
System Model	NOM_ATTR_ID_MODEL	0x0928
Product Specification	NOM_ATTR_ID_PROD_SPECN	0x092D
Object Type	NOM_ATTR_ID_TYPE	0x092F
Line Frequency	NOM_ATTR_LINE_FREQ	0x0935
System Localization	NOM_ATTR_LOCALIZN	0x0937
Metric Info Label	NOM_ATTR_METRIC_INFO_LABEL	0x093C
Metric Info Label String	NOM_ATTR_METRIC_INFO_LABEL_STR	0x093D
Metric Specification	NOM_ATTR_METRIC_SPECN	0x093F
Metric State		

	NOM_ATTR_METRIC_STAT	0x0940
Measure Mode		
	NOM_ATTR_MODE_MSMT	0x0945
Operating Mode		
	NOM_ATTR_MODE_OP	0x0946
Nomenclature Version		
	NOM_ATTR_NOM_VERS	0x0948
Compound Numeric Observed Value		
	NOM_ATTR_NU_CMPD_VAL_OBS	0x094B
Numeric Observed Value		
	NOM_ATTR_NU_VAL_OBS	0x0950
Patient BSA		
	NOM_ATTR_PT_BSA	0x0956
Pat Demo State		
	NOM_ATTR_PT_DEMOG_ST	0x0957
Patient Date of Birth		
	NOM_ATTR_PT_DOB	0x0958
Patient ID		
	NOM_ATTR_PT_ID	0x095A
Family Name		
	NOM_ATTR_PT_NAME_FAMILY	0x095C
Given Name		
	NOM_ATTR_PT_NAME_GIVEN	0x095D
Patient Sex		
	NOM_ATTR_PT_SEX	0x0961
Patient Type		
	NOM_ATTR_PT_TYPE	0x0962
Sample Array Calibration Specification		
	NOM_ATTR_SA_CALIB_I16	0x0964
Compound Sample Array Observed Value		
	NOM_ATTR_SA_CMPD_VAL_OBS	0x0967
Sample Array Physiological Range		
	NOM_ATTR_SA_RANGE_PHYS_I16	0x096A
Sample Array Specification		
	NOM_ATTR_SA_SPECN	0x096D
Sample Array Observed Value		
	NOM_ATTR_SA_VAL_OBS	0x096E
Scale and Range Specification		
	NOM_ATTR_SCALE_SPECN_I16	0x096F
Safety Standard		
	NOM_ATTR_STD_SAFETY	0x0982
System ID		
	NOM_ATTR_SYS_ID	0x0984
System Specification		
	NOM_ATTR_SYS_SPECN	0x0985
System Type		
	NOM_ATTR_SYS_TYPE	0x0986
Date and Time		
	NOM_ATTR_TIME_ABS	0x0987
Sample Period		
	NOM_ATTR_TIME_PD_SAMP	0x098D
Relative Time		
	NOM_ATTR_TIME_REL	0x098F
Absolute Time Stamp		
	NOM_ATTR_TIME_STAMP_ABS	0x0990
Relative Time Stamp		

	NOM_ATTR_TIME_STAMP_REL	0x0991
Unit Code	NOM_ATTR_UNIT_CODE	0x0996
Enumeration Observed Value	NOM_ATTR_VAL_ENUM_OBS	0x099E
MDS Status	NOM_ATTR_VMS_MDS_STAT	0x09A7
Patient Age	NOM_ATTR_PT_AGE	0x09D8
Patient Height	NOM_ATTR_PT_HEIGHT	0x09DC
Patient Weight	NOM_ATTR_PT_WEIGHT	0x09DF
Sample Array Fixed Values Specification	NOM_ATTR_SA_FIXED_VAL_SPECN	0x0A16
Patient Paced Mode	NOM_ATTR_PT_PACED_MODE	0x0A1E
Internal Patient ID	NOM_ATTR_PT_ID_INT	0xF001
Private Attribute	NOM_SAT_O2_TONE_FREQ	0xF008
Private Attribute	NOM_ATTR_CMPD_REF_LIST	0xF009
IP Address Information	NOM_ATTR_NET_ADDR_INFO	0xF100
Protocol Support	NOM_ATTR_PCOL_SUPPORT	0xF101
Notes1	NOM_ATTR_PT_NOTES1	0xF129
Notes2	NOM_ATTR_PT_NOTES2	0xF12A
Time for Periodic Polling	NOM_ATTR_TIME_PD_POLL	0xF13E
Patient BSA Formula	NOM_ATTR_PT_BSA_FORMULA	0xF1EC
Mds General System Info	NOM_ATTR_MDS_GEN_INFO	0xF1FA
no of prioritized objects for poll request	NOM_ATTR_POLL_OBJ_PRIO_NUM	0xF228
Numeric Object Priority List	NOM_ATTR_POLL_NU_PRIO_LIST	0xF239
Wave Object Priority List	NOM_ATTR_POLL_RTSA_PRIO_LIST	0xF23A
Metric Modality	NOM_ATTR_METRIC_MODALITY	0xF294

The attributes are arranged in the following attribute groups:

Alert Monitor Group	NOM_ATTR_GRP_AL_MON	0x0801
Metric Observed Value Group	NOM_ATTR_GRP_METRIC_VAL_OBS	0x0803
Patient Demographics Attribute Group	NOM_ATTR_GRP_PT_DEMOG	0x0807
System Application Attribute Group	NOM_ATTR_GRP_SYS_APPL	0x080A

System Identification Attribute Group	
NOM_ATTR_GRP_SYS_ID	0x080B
System Production Attribute Group	
NOM_ATTR_GRP_SYS_PROD	0x080C
VMO Dynamic Attribute Group	
NOM_ATTR_GRP_VMO_DYN	0x0810
VMO Static Attribute Group	
NOM_ATTR_GRP_VMO_STATIC	0x0811

Component IDs

The Component IDs specify system components such as the entries in the Production Specification attribute of the Medical Device Service object. A Component ID is a PrivateOid and is not assigned to any nomenclature partition.

for the overall product	
ID_COMP_PRODUCT	0x0008
for the specific bundle	
ID_COMP_CONFIG	0x0010
for the boot code	
ID_COMP_BOOT	0x0018
mainboard component	
ID_COMP_MAIN_BD	0x0050
application software component	
ID_COMP_APPL_SW	0x0058

Unit Codes

The Unit Codes describe the dimension of a physiological measurement. They are grouped in the Units partition.

NOS	(no dimension)	
	NOM_DIM_NOS	0
/	(/)	
	NOM_DIM_DIV	2
-	(no dimension)	
	NOM_DIM_DIMLESS	512
%	(percentage)	
	NOM_DIM_PERCENT	544
ppth	(parts per thousand)	
	NOM_DIM_PARTS_PER_THOUSAND	576
ppm	(parts per million)	
	NOM_DIM_PARTS_PER_MILLION	608
mol/mol	(mole per mole)	
	NOM_DIM_X_MOLE_PER_MOLE	864
ppb	(parts per billion)	
	NOM_DIM_PARTS_PER_BILLION	672
ppt	(parts per trillion)	
	NOM_DIM_PARTS_PER_TRILLION	704
pH	(pH)	
	NOM_DIM_PH	992
drop	(vital signs count drop)	
	NOM_DIM_DROP	1024
rbc	(vital signs count red blood cells)	
	NOM_DIM_RBC	1056
beat	(vital signs count beat)	
	NOM_DIM_BEAT	1088
breath	(vital signs count breath)	
	NOM_DIM_BREATH	1120
cell	(vital signs count cells)	
	NOM_DIM_CELL	1152
cough	(vital signs count cough)	
	NOM_DIM_COUGH	1184
sigh	(vital signs count sigh)	
	NOM_DIM_SIGH	1216
%PCV	(percent of packed cell volume)	
	NOM_DIM_PCT_PCV	1248
m	(meter)	
	NOM_DIM_X_M	1280
cm	(centimeter)	
	NOM_DIM_CENTI_M	1297
mm	(millimeter)	
	NOM_DIM_MILLI_M	1298
µm	(micro-meter)	
	NOM_DIM_MICRO_M	1299
in	(inch)	
	NOM_DIM_X_INCH	1376
ml/m2	(used e.g. for SI and ITBVI)	
	NOM_DIM_MILLI_L_PER_M_SQ	1426
/m	(per meter)	
	NOM_DIM_PER_X_M	1440
/mm	(per millimeter)	

	NOM_DIM_PER_MILLI_M	1458
m2	(used e.g. for BSA calculation) NOM_DIM_SQ_X_M	1472
in2	(used e.g. for BSA calculation) NOM_DIM_SQ_X_INCH	1504
m3	(cubic meter) NOM_DIM_CUBIC_X_M	1568
cm3	(cubic centimeter) NOM_DIM_CUBIC_CENTI_M	1585
l	(liter) NOM_DIM_X_L	1600
ml	(milli-liters used e.g. for EVLW ITBV SV) NOM_DIM_MILLI_L	1618
ml/breath	(milli-liter per breath) NOM_DIM_MILLI_L_PER_BREATH	1650
/cm3	(per cubic centimeter) NOM_DIM_PER_CUBIC_CENTI_M	1681
/l	(per liter) NOM_DIM_PER_X_L	1696
1/nl	(per nano-liter) NOM_DIM_PER_NANO_LITER	1716
g	(gram) NOM_DIM_X_G	1728
kg	(kilo-gram) NOM_DIM_KILO_G	1731
mg	(milli-gram) NOM_DIM_MILLI_G	1746
µg	(micro-gram) NOM_DIM_MICRO_G	1747
ng	(nono-gram) NOM_DIM_NANO_G	1748
lb	(pound) NOM_DIM_X_LB	1760
oz	(ounce) NOM_DIM_X_OZ	1792
/g	(per gram) NOM_DIM_PER_X_G	1824
g-m	(used e.g. for LVSW RVSW) NOM_DIM_X_G_M	1856
kg-m	(used e.g. for RCW LCW) NOM_DIM_KILO_G_M	1859
g-m/m2	(used e.g. for LVSWI and RVSWI) NOM_DIM_X_G_M_PER_M_SQ	1888
kg-m/m2	(used e.g. for LCWI and RCWI) NOM_DIM_KILO_G_M_PER_M_SQ	1891
kg-m2	(gram meter squared) NOM_DIM_KILO_G_M_SQ	1923
kg/m2	(kilo-gram per square meter) NOM_DIM_KG_PER_M_SQ	1955
kg/m3	(kilo-gram per cubic meter) NOM_DIM_KILO_G_PER_M_CUBE	1987
g/cm3	(gram per cubic meter) NOM_DIM_X_G_PER_CM_CUBE	2016
mg/cm3	(milli-gram per cubic centimeter) NOM_DIM_MILLI_G_PER_CM_CUBE	2034
µg/cm3	(micro-gram per cubic centimeter)	

	NOM_DIM_MICRO_G_PER_CM_CUBE	2035
ng/cm ³	(nano-gram per cubic centimeter)	
	NOM_DIM_NANO_G_PER_CM_CUBE	2036
g/l	(gram per liter)	
	NOM_DIM_X_G_PER_L	2048
g/dl	(used e.g. for Hb)	
	NOM_DIM_X_G_PER_DL	2112
mg/dl	(milli-gram per deciliter)	
	NOM_DIM_MILLI_G_PER_DL	2130
g/ml	(gram per milli-liter)	
	NOM_DIM_X_G_PER_ML	2144
mg/ml	(milli-gram per milli-liter)	
	NOM_DIM_MILLI_G_PER_ML	2162
µg/ml	(micro-gram per milli-liter)	
	NOM_DIM_MICRO_G_PER_ML	2163
ng/ml	(nano-gram per milli-liter)	
	NOM_DIM_NANO_G_PER_ML	2164
sec	(seconds)	
	NOM_DIM_SEC	2176
msec	(milli-seconds)	
	NOM_DIM_MILLI_SEC	2194
µsec	(micro-seconds)	
	NOM_DIM_MICRO_SEC	2195
min	(minutes)	
	NOM_DIM_MIN	2208
hrs	(hours)	
	NOM_DIM_HR	2240
days	(days)	
	NOM_DIM_DAY	2272
weeks	(weeks)	
	NOM_DIM_WEEKS	2304
months	(months)	
	NOM_DIM_MON	2336
years	(years)	
	NOM_DIM_YR	2368
TOD	(time of day)	
	NOM_DIM_TOD	2400
date	(date)	
	NOM_DIM_DATE	2432
/sec	(per second)	
	NOM_DIM_PER_X_SEC	2464
Hz	(hertz)	
	NOM_DIM_HZ	2496
/min	(per minute used e.g. for the PVC count numerical value)	
	NOM_DIM_PER_MIN	2528
/hour	(per hour)	
	NOM_DIM_PER_HR	2560
/day	(per day)	
	NOM_DIM_PER_DAY	2592
/week	(per week)	
	NOM_DIM_PER_WK	2624
/month	(per month)	
	NOM_DIM_PER_MO	2656
/year	(per year)	
	NOM_DIM_PER_YR	2688
bpm	(beats per minute used e.g. for HR/PULSE)	

	NOM_DIM_BEAT_PER_MIN	2720
puls/min	(puls per minute)	
	NOM_DIM_PULS_PER_MIN	2752
rpm	(respiration breathes per minute)	
	NOM_DIM_RESP_PER_MIN	2784
m/sec	(meter per second)	
	NOM_DIM_X_M_PER_SEC	2816
mm/sec	(speed for recordings)	
	NOM_DIM_MILLI_M_PER_SEC	2834
l/min/m2	(used for CI)	
	NOM_DIM_X_L_PER_MIN_PER_M_SQ	2848
ml/min/m2	(used for DO2I VO2I O2AVI)	
	NOM_DIM_MILLI_L_PER_MIN_PER_M_SQ	2866
m2/sec	(square meter per second)	
	NOM_DIM_SQ_X_M_PER_SEC	2880
cm2/sec	(square centimeter per second)	
	NOM_DIM_SQ_CENTI_M_PER_SEC	2897
m3/sec	(cubic meter per second)	
	NOM_DIM_CUBIC_X_M_PER_SEC	2912
cm3/sec	(cubic centimeter per second)	
	NOM_DIM_CUBIC_CENTI_M_PER_SEC	2929
l/sec	(liter per second)	
	NOM_DIM_X_L_PER_SEC	3040
l/min	(liter per minutes)	
	NOM_DIM_X_L_PER_MIN	3072
dl/min	(deciliter per second)	
	NOM_DIM_DECI_L_PER_MIN	3088
ml/min	(used for DO2 VO2 ALVENT)	
	NOM_DIM_MILLI_L_PER_MIN	3090
l/hour	(liter per hour)	
	NOM_DIM_X_L_PER_HR	3104
ml/hour	(milli-liter per hour)	
	NOM_DIM_MILLI_L_PER_HR	3122
l/day	(liter per day)	
	NOM_DIM_X_L_PER_DAY	3136
ml/day	(milli-liter per day)	
	NOM_DIM_MILLI_L_PER_DAY	3154
ml/kg	(used e.g. for EVLWI)	
	NOM_DIM_MILLI_L_PER_KG	3186
kg/sec	(kilo-gram per second)	
	NOM_DIM_KILO_G_PER_SEC	3299
g/min	(gram per minute)	
	NOM_DIM_X_G_PER_MIN	3328
kg/min	(kilo-gram per minute)	
	NOM_DIM_KILO_G_PER_MIN	3331
mg/min	(milli-gram per minute)	
	NOM_DIM_MILLI_G_PER_MIN	3346
µg/min	(micro-gram per minute)	
	NOM_DIM_MICRO_G_PER_MIN	3347
ng/min	(nano-gram per minute)	
	NOM_DIM_NANO_G_PER_MIN	3348
g/hour	(gram per hour)	
	NOM_DIM_X_G_PER_HR	3360
kg/hour	(kilo-gram per hour)	
	NOM_DIM_KILO_G_PER_HR	3363
mg/hour	(milli-gram per hour)	

	NOM_DIM_MILLI_G_PER_HR	3378
µg/hour	(micro-gram per hour)	
	NOM_DIM_MICRO_G_PER_HR	3379
ng/hr	(nano-gram per hour)	
	NOM_DIM_NANO_G_PER_HR	3380
kg/day	(kilo-gram per day)	
	NOM_DIM_KILO_G_PER_DAY	3395
g/kg/min	(gram per kilo-gram per minute)	
	NOM_DIM_X_G_PER_KG_PER_MIN	3456
mg/kg/min	(milli-gram per kilo-gram per minute)	
	NOM_DIM_MILLI_G_PER_KG_PER_MIN	3474
µg/kg/min	(micro-gram per kilo-gram per minute)	
	NOM_DIM_MICRO_G_PER_KG_PER_MIN	3475
ng/kg/min	(nano-gram per kilo-gram per minute)	
	NOM_DIM_NANO_G_PER_KG_PER_MIN	3476
g/kg/hour	(gram per kilo-gram per hour)	
	NOM_DIM_X_G_PER_KG_PER_HR	3488
mg/kg/hour	(mili-gram per kilo-gram per hour)	
	NOM_DIM_MILLI_G_PER_KG_PER_HR	3506
µg/kg/hour	(micro-gram per kilo-gram per hour)	
	NOM_DIM_MICRO_G_PER_KG_PER_HR	3507
ng/kg/hour	(nano-gram per kilo-gram per hour)	
	NOM_DIM_NANO_G_PER_KG_PER_HR	3508
kg/l/sec	(kilo-gram per liter per second)	
	NOM_DIM_KILO_G_PER_L_SEC	3555
kg/m/sec	(kilo-gram per meter per second)	
	NOM_DIM_KILO_G_PER_M_PER_SEC	3683
kg-m/sec	(kilo-gram meter per second)	
	NOM_DIM_KILO_G_M_PER_SEC	3715
N-s	(newton seconds)	
	NOM_DIM_X_NEWTON_SEC	3744
N	(newton)	
	NOM_DIM_X_NEWTON	3776
Pa	(pascal)	
	NOM_DIM_X_PASCAL	3840
hPa	(hekto-pascal)	
	NOM_DIM_HECTO_PASCAL	3842
kPa	(kilo-pascal)	
	NOM_DIM_KILO_PASCAL	3843
mmHg	(mm mercury)	
	NOM_DIM_MMHG	3872
cmH2O	(centimeter H2O)	
	NOM_DIM_CM_H2O	3904
mBar	(milli-bar)	
	NOM_DIM_MILLI_BAR	3954
J	(Joules)	
	NOM_DIM_X_JOULES	3968
eV	(electronvolts)	
	NOM_DIM_EVOLT	4000
W	(watt)	
	NOM_DIM_X_WATT	4032
mW	(milli-watt)	
	NOM_DIM_MILLI_WATT	4050
nW	(nano-watt)	
	NOM_DIM_NANO_WATT	4052
pW	(pico-watt)	

	NOM_DIM_PICO_WATT	4053
Dyn-sec/cm ⁵	(dyne second per cm ⁵)	
	NOM_DIM_X_DYNE_PER_SEC_PER_CM5	4128
A	(ampere)	
	NOM_DIM_X_AMPS	4160
mA	(milli-ampere used e.g. for the battery indications)	
	NOM_DIM_MILLI_AMPS	4178
C	(coulomb)	
	NOM_DIM_X_COULOMB	4192
μC	(micro-coulomb)	
	NOM_DIM_MICRO_COULOMB	4211
V	(volts)	
	NOM_DIM_X_VOLT	4256
mV	(milli-volt)	
	NOM_DIM_MILLI_VOLT	4274
μV	(micro-volt)	
	NOM_DIM_MICRO_VOLT	4275
Ohm	(Ohm)	
	NOM_DIM_X_OHM	4288
kOhm	(kilo-ohm)	
	NOM_DIM_OHM_K	4291
F	(farad)	
	NOM_DIM_X_FARAD	4352
°K	(kelvin)	
	NOM_DIM_KELVIN	4384
°F	(degree-fahrenheit)	
	NOM_DIM_FAHR	4416
cd	(candela)	
	NOM_DIM_X_CANDELA	4480
mOsm	(milli-osmole)	
	NOM_DIM_MILLI_OSM	4530
mol	(mole)	
	NOM_DIM_X_MOLE	4544
mmol	(milli-mole)	
	NOM_DIM_MILLI_MOLE	4562
mEq	(milli-equivalents)	
	NOM_DIM_MILLI_EQUIV	4594
mOsm/l	(milli-osmole per liter)	
	NOM_DIM_MILLI_OSM_PER_L	4626
mmol/l	(used for HB)	
	NOM_DIM_MILLI_MOLE_PER_L	4722
μmol/l	(micro-mol per liter)	
	NOM_DIM_MICRO_MOLE_PER_L	4723
mEq/l	(milli-equivalents per liter)	
	NOM_DIM_MILLI_EQUIV_PER_L	4850
mEq/day	(milli-equivalents per day)	
	NOM_DIM_MILLI_EQUIV_PER_DAY	5202
i.u.	(international unit)	
	NOM_DIM_X_INTL_UNIT	5472
mi.u.	(mili-international unit)	
	NOM_DIM_MILLI_INTL_UNIT	5490
i.u./cm ³	(international unit per cubic centimeter)	
	NOM_DIM_X_INTL_UNIT_PER_CM_CUBE	5504
mi.u./cm ³	(mili-international unit per cubic centimeter)	
	NOM_DIM_MILLI_INTL_UNIT_PER_CM_CUBE	5522
i.u./ml	(international unit per milli-liter)	

	NOM_DIM_X_INTL_UNIT_PER_ML	5600
i.u./min	(international unit per minute)	
	NOM_DIM_X_INTL_UNIT_PER_MIN	5664
mi.u./ml	(milli-international unit per milli-liter)	
	NOM_DIM_MILLI_INTL_UNIT_PER_ML	5618
mi.u./min	(milli-international unit per minute)	
	NOM_DIM_MILLI_INTL_UNIT_PER_MIN	5682
i.u./hour	(international unit per hour)	
	NOM_DIM_X_INTL_UNIT_PER_HR	5696
mi.u./hour	(milli-international unit per hour)	
	NOM_DIM_MILLI_INTL_UNIT_PER_HR	5714
i.u./kg/min	(international unit per kilo-gram per minute)	
	NOM_DIM_X_INTL_UNIT_PER_KG_PER_MIN	5792
mi.u./kg/min	(milli-international unit per kilo-gram per minute)	
	NOM_DIM_MILLI_INTL_UNIT_PER_KG_PER_MIN	5810
i.u./kg/hour	(international unit per kilo-gram per hour)	
	NOM_DIM_X_INTL_UNIT_PER_KG_PER_HR	5824
mi.u./kg/hour	(milli-international unit per kilo-gram per hour)	
	NOM_DIM_MILLI_INTL_UNIT_PER_KG_PER_HR	5842
ml/cmH2O	(milli-liter per centimeter H2O)	
	NOM_DIM_MILLI_L_PER_CM_H2O	5906
cmH2O/l/sec	(centimeter H2O per second)	
	NOM_DIM_CM_H2O_PER_L_PER_SEC	5920
ml2/sec	(milli-liter per second)	
	NOM_DIM_MILLI_L_SQ_PER_SEC	5970
cmH2O/%	(centimeter H2O per percent)	
	NOM_DIM_CM_H2O_PER_PERCENT	5984
DS*m2/cm5	(used for SVRI and PVRI)	
	NOM_DIM_DYNE_SEC_PER_M_SQ_PER_CM_5	6016
°C	(degree-celsius)	
	NOM_DIM_DEGC	6048
cmH2O/l	(centimeter H2O per liter)	
	NOM_DIM_CM_H2O_PER_L	6144
mmHg/%	(milli-meter mercury per percent)	
	NOM_DIM_MM_HG_PER_PERCENT	6176
kPa/%	(kilo-pascal per percent)	
	NOM_DIM_KILO_PA_PER_PERCENT	6211
l/mmHg	(liter per mmHg)	
	NOM_DIM_X_L_PER_MM_HG	6272
ml/mmHg	(milli-liter per milli-meter Hg)	
	NOM_DIM_MILLI_L_PER_MM_HG	6290
mAh	(milli-ampere per hour used e.g. for the battery indications)	
	NOM_DIM_MILLI_AMP_HR	6098
ml/dl	(used for CaO2 CvO2 Ca-vO2)	
	NOM_DIM_MILLI_L_PER_DL	6418
dB	(decibel)	
	NOM_DIM_DECIBEL	6432
g/mg	(gram per milli-gram)	
	NOM_DIM_X_G_PER_MILLI_G	6464
mg/mg	(milli-gram per milli-gram)	
	NOM_DIM_MILLI_G_PER_MILLI_G	6482
bpm/l	(beats per minute per liter)	
	NOM_DIM_BEAT_PER_MIN_PER_X_L	6496
bpm/ml	(beats per minute per milli-liter)	
	NOM_DIM_BEAT_PER_MIN_PER_MILLI_L	6514
1/(min*l)	(per minute per liter)	

	NOM_DIM_PER_X_L_PER_MIN	6528
m/min	(meter per minute)	
	NOM_DIM_X_M_PER_MIN	6560
cm/min	(speed for recordings)	
	NOM_DIM_CENTI_M_PER_MIN	6577
pg/ml	(pico-gram per milli-liter)	
	NOM_DIM_PICO_G_PER_ML	2165
ug/l	(micro-gram per liter)	
	NOM_DIM_MICRO_G_PER_L	2067
ng/l	(nano-gram per liter)	
	NOM_DIM_NANO_G_PER_L	2068
/mm ³	(per cubic millimeter)	
	NOM_DIM_PER_CUBIC_MILLI_M	1682
mm ³	(cubic milli-meter)	
	NOM_DIM_CUBIC_MILLI_M	1586
u/l	(intl. units per liter)	
	NOM_DIM_X_INTL_UNIT_PER_L	5568
/l	(10 ⁶ intl. units per liter)	
	NOM_DIM_MEGA_INTL_UNIT_PER_L	5573
mol/kg	(mole per kilo-gram)	
	NOM_DIM_MILLI_MOL_PER_KG	4946
mcg/dl	(micro-gram per deci-liter)	
	NOM_DIM_MICRO_G_PER_DL	2131
mg/l	(milli-gram per liter)	
	NOM_DIM_MILLI_G_PER_L	2066
/ul	(micro-liter)	
	NOM_DIM_PER_MICRO_L	1715
complx	(-)	
	NOM_DIM_COMPLEX	61440
count	(count as a dimension)	
	NOM_DIM_COUNT	61441
part	(part)	
	NOM_DIM_PART	61442
puls	(puls)	
	NOM_DIM_PULS	61443
μV p-p	(micro-volt peak to peak)	
	NOM_DIM_UV_PP	61444
μV ²	(micor-volt square)	
	NOM_DIM_UV_SQ	61445
lumen	(lumen)	
	NOM_DIM_LUMEN	61447
lb/in ²	(pound per square inch)	
	NOM_DIM_LB_PER_INCH_SQ	61448
mmHg/s	(milli-meter mercury per second)	
	NOM_DIM_MM_HG_PER_SEC	61449
ml/s	(milli-liter per second)	
	NOM_DIM_ML_PER_SEC	61450
bpm/ml	(beat per minute per milli-liter)	
	NOM_DIM_BEAT_PER_MIN_PER_ML_C	61451
J/day	(joule per day)	
	NOM_DIM_X_JOULE_PER_DAY	61536
kJ/day	(kilo joule per day)	
	NOM_DIM_KILO_JOULE_PER_DAY	61539
MJ/day	(mega joule per day)	
	NOM_DIM_MEGA_JOULE_PER_DAY	61540
cal	(calories)	

	NOM_DIM_X_CALORIE	61568
kcal	(kilo calories)	
	NOM_DIM_KILO_CALORIE	61571
10**6 cal	(million calories)	
	NOM_DIM_MEGA_CALORIE	61572
cal/day	(calories per day)	
	NOM_DIM_X_CALORIE_PER_DAY	61600
kcal/day	(kilo-calories per day)	
	NOM_DIM_KILO_CALORIE_PER_DAY	61603
Mcal/day	(mega calories per day)	
	NOM_DIM_MEGA_CALORIE_PER_DAY	61604
cal/ml	(calories per milli-liter)	
	NOM_DIM_X_CALORIE_PER_ML	61632
kcal/ml	(kilo calories per ml)	
	NOM_DIM_KILO_CALORIE_PER_ML	61635
J/ml	(Joule per milli-liter)	
	NOM_DIM_X_JOULE_PER_ML	61664
kJ/ml	(kilo-joules per milli-liter)	
	NOM_DIM_KILO_JOULE_PER_ML	61667
RPM	(revolutions per minute)	
	NOM_DIM_X_REV_PER_MIN	61696
l/(mn*l*kg)	(per minute per liter per kilo)	
	NOM_DIM_PER_L_PER_MIN_PER_KG	61728
l/mbar	(liter per milli-bar)	
	NOM_DIM_X_L_PER_MILLI_BAR	61760
ml/mbar	(milli-liter per milli-bar)	
	NOM_DIM_MILLI_L_PER_MILLI_BAR	61778
l/kg/hr	(liter per kilo-gram per hour)	
	NOM_DIM_X_L_PER_KG_PER_HR	61792
ml/kg/hr	(milli-liter per kilogram per hour)	
	NOM_DIM_MILLI_L_PER_KG_PER_HR	61810
bar/l/s	(bar per liter per sec)	
	NOM_DIM_X_BAR_PER_LITER_PER_SEC	61824
mbar/l/s	(milli-bar per liter per sec)	
	NOM_DIM_MILLI_BAR_PER_LITER_PER_SEC	61842
bar/l	(bar per liter)	
	NOM_DIM_X_BAR_PER_LITER	61856
mbar/l	(bar per liter)	
	NOM_DIM_MILLI_BAR_PER_LITER	61874
V/mV	(volt per milli-volt)	
	NOM_DIM_VOLT_PER_MILLI_VOLT	61888
cmH2O/uV	(cm H2O per micro-volt)	
	NOM_DIM_CM_H2O_PER_MICRO_VOLT	61920
J/l	(joule per liter)	
	NOM_DIM_X_JOULE_PER_LITER	61952
l/bar	(liter per bar)	
	NOM_DIM_X_L_PER_BAR	61984
m/mV	(meter per milli-volt)	
	NOM_DIM_X_M_PER_MILLI_VOLT	62016
mm/mV	(milli-meter per milli-volt)	
	NOM_DIM_MILLI_M_PER_MILLI_VOLT	62034
l/min/kg	(liter per minute per kilo-gram)	
	NOM_DIM_X_L_PER_MIN_PER_KG	62048
ml/min/kg	(milli-liter per minute per kilo-gram)	
	NOM_DIM_MILLI_L_PER_MIN_PER_KG	62066
Pa/l/s	(pascal per liter per sec)	

	NOM_DIM_X_PASCAL_PER_L_PER_SEC	62080
hPa/l/s	(hPa per liter per sec)	
	NOM_DIM_HECTO_PASCAL_PER_L_PER_SEC	62082
kPa/l/s	(kPa per liter per sec)	
	NOM_DIM_KILO_PASCAL_PER_L_PER_SEC	62083
ml/Pa	(milli-liter per pascal)	
	NOM_DIM_MILLI_L_PER_X_PASCAL	62112
ml/hPa	(milli-liter per hecto-pascal)	
	NOM_DIM_MILLI_L_PER_HECTO_PASCAL	62114
ml/kPa	(milli-liter per kilo-pascal)	
	NOM_DIM_MILLI_L_PER_KILO_PASCAL	62115
mmHg/l/s	(mm)	
	NOM_DIM_MM_HG_PER_X_L_PER_SEC	62144
mol/h	(mole per hour)	
	NOM_DIM_X_MOLE_PER_HR	62176
mmol/h	(milli-mol per hour)	
	NOM_DIM_MILLI_MOLE_PER_HR	62194
l/beat	(liter per beat)	
	NOM_DIM_X_L_PER_BEAT	62208
ml/beat	(milli-liter per beat)	
	NOM_DIM_MILLI_L_PER_BEAT	62226
l/beat/m2	(liter per beat per square meter)	
	NOM_DIM_X_L_PER_BEAT_PER_M_SQ	62240
ml/beat/m2	(milli-liter per beat per square meter)	
	NOM_DIM_MILLI_L_PER_BEAT_PER_M_SQ	62258
bar/s	(bar per second)	
	NOM_DIM_X_BAR_PER_SEC	62272
mbar/s	(milli-bar per second)	
	NOM_DIM_MILLI_BAR_PER_SEC	62290
pascal/l	(pascal per liter)	
	NOM_DIM_X_PASCAL_PER_L	62304
kpascal/l	(kilo pascal per liter)	
	NOM_DIM_KILO_PASCAL_PER_L	62307
mmHg/l	(mmHg per liter)	
	NOM_DIM_MM_HG_PER_L	62336
vol%/l	(volume percent per liter)	
	NOM_DIM_VOL_PERCENT_PER_L	62368

Alert Codes

The first column in the tables below shows the alert source, the second column shows the associated alert code and the third column contains the alert text which would be displayed by the monitor. The XXX in the alert text is a placeholder for the actual alert source. It is filled depending on the alert source. Note that the alert text depends on the localization of your monitor.

The least significant bit of the alert codes listed below is used to identify the source of an alert (refer to “Alert Monitor Object” on page 99). If the alert code is marked with a (*), the associated alert source is from the object oriented nomenclature partition and hence the least significant bit of the alert code is set to 1.

NOTE On monitors with SW Rev. G.0 or lower some alert codes will only be issued correctly, if the connected MMS or FMS has the same or a higher SW revision.

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The following code from the object oriented partition are used for the alert source:

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The following codes from the event partition are used for the alert code:

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Private Unicode Characters

The IntelliVue monitor may use the following private codes for UNICODE characters:

```
#define SUBSCRIPT_CAPITAL_E_CHAR      0xE145
/* SUBSCRIPT CAPITAL E                */
#define SUBSCRIPT_CAPITAL_L_CHAR      0xE14C
/* SUBSCRIPT CAPITAL L                */
#define LITER_PER_CHAR                 0xE400
/* LITER PER - used in 4 char unit "l/min" */
#define HYDROGEN_CHAR                  0xE401
/* HYDROGEN - Used in 4 char unit "cmH2O" */
#define ALARM_STAR_CHAR                 0xE40D
/* ALARM STAR                          */
#define CAPITAL_V_WITH_DOT_ABOVE_CHAR  0xE425
/* CAPITAL_V_WITH_DOT_ABOVE (V with dot) */
#define ZERO_WIDTH_NO_BREAK_SPACE_CHAR 0xFEFF
/* The character 0xFEFF is used as FILL character.
   For each wide asian character, a FILL character is
   appended for size calculations. */
```

List of Constants Used Within the Protocol Definition

RO Types

```
#define ROIV_APDU      1
#define RORS_APDU      2
#define ROER_APDU      3
#define ROLRS_APDU     5
```

ROLRS Identifier

```
#define RORLS_FIRST      1    /* set in the first message */
#define RORLS_NOT_FIRST_NOT_LAST 2
#define RORLS_LAST      3    /* last RORLSapdu, one RORSapdu
                               to follow */
```

ROSE Commands

```
typedef u_16             CMDType;
#define CMD_EVENT_REPORT      0
#define CMD_CONFIRMED_EVENT_REPORT 1
#define CMD_GET               3
#define CMD_SET               4
#define CMD_CONFIRMED_SET     5
#define CMD_CONFIRMED_ACTION  7
```

ROER Error Values

```
#define NO_SUCH_OBJECT_CLASS      0
#define NO_SUCH_OBJECT_INSTANCE  1
#define ACCESS_DENIED            2
#define GET_LIST_ERROR           7
#define SET_LIST_ERROR           8
#define NO_SUCH_ACTION           9
#define PROCESSING_FAILURE       10
#define INVALID_ARGUMENT_VALUE   15
#define INVALID_SCOPE            16
#define INVALID_OBJECT_INSTANCE  17
```

Action and Event Types

The Action and Event Types are defined in the Object Oriented Elements partition of the nomenclature.

```
#define NOM_NOTI_MDS_CREAT      3334
/* MDS Create Notification */
#define NOM_NOTI_CONN_INDIC    3351
/* connect indication event type */
#define NOM_ACT_POLL_MDIB_DATA 3094
/* poll data action */
#define NOM_ACT_POLL_MDIB_DATA_EXT 61755
/* extended poll data action */
```

Protocol Identification

The IDs for the protocol identification are from the Infrastructure nomenclature partition.

```
#define NOM_POLL_PROFILE_SUPPORT 1
/* id for polling profile */
#define NOM_MDIB_OBJ_SUPPORT    258
/* supported objects for the active profile */
#define NOM_ATTR_POLL_PROFILE_EXT 61441
/* id for poll profile extensions opt. package */
```

Association Control

```

#define MDDL_VERSION1                0x80000000
    /* Data Export Protocol Version */
#define NOMEN_VERSION                0x40000000
    /* Nomenclature Version */
#define SYST_CLIENT                  0x80000000
    /* System Type Client */
#define SYST_SERVER                  0x00800000
    /* System Type Server */
#define HOT_START                    0x80000000
    /* Startup Mode Hotstart */
#define WARM_START                    0x40000000
    /* Startup Mode Warmstart */
#define COLD_START                    0x20000000
    /* Startup Mode Coldstart */
#define POLL_PROFILE_REV_0           0x80000000
    /* Poll Profile Revision */
#define P_OPT_DYN_CREATE_OBJECTS     0x40000000
    /* option dynamic object creation */
#define P_OPT_DYN_DELETE_OBJECTS     0x20000000
    /* option dynamic object deletion */
#define POLL_EXT_PERIOD_NU_1SEC      0x80000000
    /* 1 sec Real-time Numerics */
#define POLL_EXT_PERIOD_NU_AVG_12SEC 0x40000000
    /* 12 sec averaged Numerics */
#define POLL_EXT_PERIOD_NU_AVG_60SEC 0x20000000
    /* 1 min. averaged Numerics */
#define POLL_EXT_PERIOD_NU_AVG_300SEC 0x10000000
    /* 5 min. averaged Numerics */
#define POLL_EXT_PERIOD_RTSA         0x08000000
    /* allow enumeration objects */
#define POLL_EXT_ENUM                0x04000000
    /* allow numeric priority list to be set */
#define POLL_EXT_NU_PRIO_LIST        0x02000000
    /* send timestamps for numerics with dynamic modalities */
#define POLL_EXT_DYN_MODALITIES      0x01000000

```

Label Mapping Table

With IntelliVue release G the nomenclature of some numeric and wave labels have been changed. The labels that previously resided in the namespace NOM_EMFC mainly used by VueLink devices have been moved into the NOM_SCADA namespace partition and the new defined NOM_SETTING namespace partition.

If your want to integrate support for the new nomenclature definitions in your existing client application you have to accept both label ids.

To guide you through the transition of the nomenclature changes introduced in release G, you may find the following table useful. The revision F label is given first followed by the new label. For further descriptions of the old labels see the Revision F of the Data Export Programmers Guide.

Label Definition	Label Id
NLS_NOM_EMFC_sAVDel	(0x040180CC)
NLS_NOM_SETT_APNEA_ALARM_DELAY	(0x0402F8D9)
NLS_NOM_EMFC_C20_PER_C	(0x04010E78)
NLS_NOM_C20_PER_C_INDEX	(0x0002F81A)
NLS_NOM_EMFC_Rf_V5	(0x0401075C)
NLS_NOM_ECG_AMPL_ST_BASELINE_V5	(0x0002F417)
NLS_NOM_EMFC_Urine	(0x04010BD8)
NLS_NOM_FLOW_URINE_PREV_24HR	(0x0002F883)
NLS_NOM_EMFC_PT	(0x040105E4)
NLS_NOM_TIME_PD_PT	(0x0002F18B)
NLS_NOM_EMFC_SerCa	(0x0401059C)
NLS_NOM_CONC_CA_SER	(0x0002F824)
NLS_NOM_EMFC_sBPAl	(0x0401A024)
NLS_NOM_SETT_PRESS_AL_ONOFF	(0x0402F8F7)
NLS_NOM_EMFC_SetTmp	(0x04010AD8)
NLS_NOM_TEMP_BODY	(0x00024B5C)
NLS_NOM_EMFC_sCO2Wm	(0x0401815C)
NLS_NOM_SETT_VENT_CO2_WARMING_MONITOR_ONOFF	(0x0402F915)
NLS_NOM_EMFC_sAPkF1	(0x04018030)
NLS_NOM_SETT_FLOW_AWAY_INSP_APNEA	(0x0402F8ED)
NLS_NOM_EMFC_SerGlc	(0x04010590)
NLS_NOM_CONC_GLU_SER	(0x0002F82A)
NLS_NOM_EMFC_RT_PCT_BE	(0x04010810)
NLS_NOM_EEG_PWR_SPEC_BETA_REL_RIGHT	(0x0002F860)
NLS_NOM_EMFC_T4	(0x04010414)
NLS_NOM_TEMP_GEN_4	(0x0002F0CA)
NLS_NOM_EMFC_GOT	(0x0401060C)
NLS_NOM_CONC_GOT	(0x0002F188)
NLS_NOM_EMFC_highO2	(0x0401A020)
NLS_NOM_SETT_VENT_CONC_AWAY_O2_LIMIT_HI	(0x0402F919)
NLS_NOM_EMFC_MCV	(0x040105D4)
NLS_NOM_VOL_CORP_MEAN	(0x0002F8C4)

NLS_NOM_EMFC_sEnTrg	(0x040180B4)
NLS_NOM_SETT_TRIG_ONOFF	(0x0402F90C)
NLS_NOM_EMFC_Plts	(0x040105D0)
NLS_NOM_PLTS_CNT	(0x0002F167)
NLS_NOM_EMFC_sLinPr	(0x04018100)
NLS_NOM_SETT_PRESS_AWAY_MIN	(0x040250F2)
NLS_NOM_EMFC_GGT	(0x04010608)
NLS_NOM_CONC_GGT	(0x0002F189)
NLS_NOM_EMFC_sAGTWm	(0x0401816C)
NLS_NOM_SETT_VENT_AGENT_WARMING_MONITOR_ONOFF	(0x0402F90D)
NLS_NOM_EMFC_sAPVhP	(0x0401807C)
NLS_NOM_SETT_VENT_PRESS_AWAY_MAX_PV_APNEA	(0x0402F931)
NLS_NOM_EMFC_sfgSEV	(0x040181AC)
NLS_NOM_SETT_CONC_AWAY_SEVOFL	(0x040251E4)
NLS_NOM_EMFC_highMV	(0x0401A02C)
NLS_NOM_SETT_VENT_VOL_MINUTE_AWAY_LIMIT_HI	(0x0402F94B)
NLS_NOM_EMFC_P6_MEAN	(0x04010407)
NLS_NOM_PRESS_GEN_6_MEAN	(0x0002F3FB)
NLS_NOM_EMFC_SpRR	(0x04010BF4)
NLS_NOM_RESP_RATE_SPONT	(0x0002F828)
NLS_NOM_EMFC_Sample	(0x04010AAC)
NLS_NOM_SETT_SAMPLE	(0x0402F956)
NLS_NOM_EMFC_CK_MM	(0x04010604)
NLS_NOM_CONC_CREA_KIN_MM	(0x0002F17F)
NLS_NOM_EMFC_sFlas	(0x040181F8)
NLS_NOM_SETT_VENT_FLOW_AWAY_ASSIST	(0x0402F91C)
NLS_NOM_EMFC_RBC	(0x040105CC)
NLS_NOM_RB_CNT	(0x0002F169)
NLS_NOM_EMFC_TOF4	(0x04010DCC)
NLS_NOM_TRAIN_OF_FOUR_4	(0x0002F8AA)
NLS_NOM_EMFC_sSens	(0x04018188)
NLS_NOM_SETT_SENS_LEVEL	(0x0402F904)
NLS_NOM_EMFC_sSIMV	(0x04018118)
NLS_NOM_SETT_VENT_MODE_SYNC_MAND_INTERMIT	(0x0402F924)
NLS_NOM_EMFC_UrCa	(0x04010624)
NLS_NOM_CONC_CA_URINE	(0x0002F19C)
NLS_NOM_EMFC_vECG	(0x0401119C)
NLS_NOM_ELEC_POTL_VECT	(0x0002F874)
NLS_NOM_EMFC_PCO2_ADJ	(0x04010A7C)
NLS_NOM_CONC_PCO2_GEN_ADJ	(0x0002F834)
NLS_NOM_EMFC_BLANK	(0x04010960)
NLS_NOM_METRIC_NOS	(0x0002EFFF)
NLS_NOM_EMFC_sPIP	(0x040180FC)
NLS_NOM_SETT_PRESS_AWAY_INSP_MAX	(0x04025109)

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Building a Computer Client

Interfacing the LAN interface with UDP/IP

When setting up a Computer Client, a network traffic analyzing tool can be useful to verify the success of each step. Widely used tools are:

- Microsoft® Network Monitor
- Tcpdump (available under the GNU Public License from <ftp://ftp.ee.lbl.gov/>)

Setting Up the BootP Server

Step 1: Connect the Computer Client to the IntelliVue monitor.

The Computer Client and the IntelliVue monitor should be connected with a crossover LAN cable. If you need a dedicated system to run the BootP server, use a hub/switch to connect the devices. It is strongly recommended that a dedicated network is used for the data export. Do not connect any additional devices.

Step 2: Start the BootP server.

Please refer to the documentation of your BootP/DHCP server for installation guidelines. If you use a DHCP server, make sure that the server supports BootP clients.

Step 3: Verify that the IntelliVue monitor receives a valid IP address.

Use a network monitor to verify that the IntelliVue monitor receives the correct IP address. If the IntelliVue monitor shows an **Unsupported LAN INOP**, it has not received a valid IP address.

If the IntelliVue monitor does not receive an IP address,

- make sure that there is no IP address conflict on the network
- try to reboot the IntelliVue monitor.

Parsing the Connect Indication Message

Step1: Verify that the Connect Indication message is sent.

Use a networked monitor and check that the IntelliVue monitor sends a subnet broadcast message to the Connect Indication port (24005). If the IntelliVue monitor does not send the message, verify that the IntelliVue monitor received a valid IP address from the BootP server (see “Setting Up the BootP Server” on page 279).

Step 2: Receive the Connect Indication message on the Computer Client.

Open a socket on the Computer Client that receives the subnet broadcast message. If the Computer Client does not receive the Connect Indication message, verify the correct network connection, use an ICMP echo (ping) to check connectivity of the IntelliVue monitor.

Step 3: Parse the Data Export Protocol Command.

The Computer Client must parse the Connect Indication message to determine the port for the Data Export Protocol. The message also contains the IP address of the IntelliVue monitor.

The Computer Client should check that all length and type fields in the message are set correctly, otherwise the message must be discarded.

Then the Computer Client should parse the appended *AttributeList* and extract the IP address and port information (refer to “Connect Indication Attributes” on page 107 for the specification of these attributes).

Interfacing the MIB/RS232 Interface with the Fixed Baudrate Protocol

Step 1: Connect the Computer Client to the IntelliVue monitor.

It may be useful to try out the Association Request/Response mechanism on the LAN interface before working with the MIB/RS232 interface. This might help to find out whether an error is related to a ill-formatted Data Export message or if it is related to a transport layer problem.

Step 2: Implement the framing algorithm.

The section “The Fixed Baudrate Protocol, RS232 Port Settings” on page 30 contains some examples which can be used to check if your framing algorithm works correctly. Remember to apply the framing algorithm to both the *Hdr* and *User Data* part of the message.

If you have tried out the Association Request message on the LAN interface, you can try to send the message within the Fixed Baudrate protocol. Just add the *Hdr* information and apply the framing algorithm.

You should keep the following points in mind when implementing the the Fixed Baudrate protocol:

- Verify that the checksum algorithm works correctly for received messages, i.e., make sure that received messages with a corrupt checksum are discarded.
- Make sure that you implement an exception handling in case the received message grows larger than your receive buffer (e.g., if an end of frame character is lost somewhere during communication).
- The Fixed Baudrate Protocol is not connection oriented. After starting your application, there may be an existing Data Export Association (either from running your own application previously or from another system which has been connected to the MIB/RS232 Interface before). This may have some unexpected consequences for your application.

Interfacing the MIB/RS232 Interface with the AutoSpeed Protocol

Step 1: Connect the Computer Client and the IntelliVue monitor.

If your operating system comes with an IrDA stack, please refer to the documentation of your operating system. The operating system will cover most of the steps below automatically.

It may be useful to try out the Association Request/Response mechanism on the LAN interface before working with the MIB/RS232 interface. This may help you to find out whether an error is related to a wrongly-formatted Data Export message or if it is related to a transport layer problem.

Step 2: Establish an IrDA connection

The IrDA protocol supports a device detection procedure. If the detection is successful, it will return information about the detected device. This information contains a device nickname and a service hints field which indicates that the device supports the IEEE 1073 standard.

After this the Computer Client can establish an IrLAP connection with the device. This involves the negotiation of the baudrate and packet size for the lower layers. Refer to the Serial Infrared Link Access Protocol (IrLAP) specification (see page 32) for more information on this topic.

Step 3: Query the IAS database

The IAS database contains the object “IEEE:1073:3:2:MDDL” with the attribute “IrDA:TinyTP:LsapSel”. This attribute contains the number of the TinyTP Service Access point for the Data Export protocol. The value type of the attribute is an *integer*. The value should be equal to 1 if the MIB/RS232 Interface is used for Data Export.

The database also contains an object named “IEEE:1073:3:2” with the attribute “NodeType”. This attribute is of type integer and specifies the type of driver which resides on the interface. A value of 1 indicates that it is a data source, i.e. it is used to export data from the monitor.

After finishing the IAS query, the Computer Client should close the IAS connection before connecting to the TinyTP Service Access Point.

Step 4: Connect to the IEEE:1073:3:2:MDDL TinyTP Service Access Point

After connecting to the TinyTP Service Access Point, the connection can be used to send Association Control and Data Export Protocol messages within TinyTP data packets.

You should check the following points for your IrDA protocol stack:

- The connection may be interrupted or reset due to communication problems (e.g., if the cable is disconnected, or the monitor is rebooted). The Computer Client should be able to recover from such problems and initiate a new connection. Note: when a disconnect occurs on the IrDA protocol layer, an Association on the Data Export protocol layer will be terminated automatically.
- The Data Export protocol is packet oriented, this means that data is exchanged as a sequence of packets. Your IrDA stack may or may not provide a packet oriented interface to the TinyTP layer. The Data Export software requires that a received IrDA packet contains only one Data Export Protocol message.

Establishing an Association

Step 1: Send an Association Request message to the IntelliVue monitor.

Format an Association Request message as described in the section “Association Request Message” on page 67. Make sure that all length fields are set correctly, the right byte order is used, and the compiler does not insert extra bytes for structure alignment.

Step 2: Parse the Association Response message sent by the IntelliVue monitor.

Verify that the IntelliVue monitor sends an Association Response message.

If the IntelliVue monitor does not send a Response message, this can have the following reasons:

- The Association Request message has been sent to the wrong port.
- The IntelliVue monitor is connected to a central station or has been connected to one (reboot the IntelliVue monitor).
- The Association Request message was not formatted correctly.

If the IntelliVue monitor sends a Refuse message, this can have the following reasons:

- The Association Request message was not formatted correctly or requested a protocol that is not supported by the IntelliVue monitor.
- The IntelliVue monitor already has an association with a different Computer Client on the same interface.
- The IntelliVue monitor already has an association with a different Computer Client on another interface and the active association uses a different source for the numeric data (only one source for numeric data may be active at a time). Please refer to “Association Request Message” on page 67 for more information on the different sources for numeric data.

If the Computer Client has an association with the IntelliVue monitor and sends a second Association Request from the same source port, the message is discarded.

Look for the byte sequence described in “Association Response Message” on page 73 to find the beginning of the User Data. Parse the User Data and make sure that the IntelliVue monitor sets the protocol versions and options as expected. Check that the requested optional packages are present.

Step 3: Parse the MDS Create Event message.

The IntelliVue monitor will send the MDS Create Event message shortly after the Association Response message. The Computer Client should parse the message and extract all necessary information. Refer to the section “Wave Objects” on page 82 for a description of the available attributes.

Step 4: Send an MDS Create Result message.

The Computer Client must send an MDS Create Result message to confirm the MDS Create Event message. Refer to “MDS CREATE EVENT RESULT” on page 55 to see how the message is formatted.

Make sure that the message uses the correct presentation context ID.

It is important that the result message has the same invoke ID as the MDS Create Event message.

If the IntelliVue monitor receives a correct MDS Create Result message, it stops re-sending MDS Create Event messages. Use a network monitor to verify this.

Step 5: Send a Release Request message.

Use the building blocks from the section “RELEASE REQUEST” on page 301 to build a Release Request message and send it to the IntelliVue monitor.

The IntelliVue monitor identifies a Computer Client based on its IP address and the source port of the messages. The Computer Client must use the same source port as in the Association Request for all communication during the association. If a message is sent from another source port, it will be treated as a message from a different Computer Client.

Step 6: Parse the Release Response message

The IntelliVue monitor sends the Release Response message to confirm that the association has been terminated. For the Computer Client it is sufficient to check the session header of the response and verify that it is indeed a Release Response message (see “Release Response” on page 73).

If the Computer Client does not receive the response message, it should try to resend the Release Request message.

To identify the IntelliVue Monitor software revision, poll the MDS objects system production attribute group and read the ProductSpecification attribute. (see “Attribute: Production Specification” on page 95)

Accessing Data

Step 1: Establish an association as described above.

Step 2: Send a Poll Data Request message to the IntelliVue monitor.

Message Frequencies

If the Computer Client sends Protocol Messages with a high frequency, the IntelliVue monitor is not able to process all the requests. Some of the messages will be discarded. The Computer Client can detect discarded Poll Data Request messages by checking the poll number in the response. The Computer Client must set the poll number so that it will be able to detect loss of messages.

Single and Extended Polling

If the Computer Client needs to access real-time numeric or wave data, it should use Poll Profile Extensions (see “EXTENDED POLL DATA REQUEST” on page 59). This avoids sending poll requests with a high frequency and reduces the communication overhead.

The Computer Client can use an Extended Poll Request only to access Numerics, Waves and Alarms. It must use Single Poll Data Requests to access data from Patient Demographics or from the Medical Device System object.

Receive the Poll Data Response message and parse it.

The IntelliVue monitor sends a Single or Extended Poll Data Result message if the Poll Request message was parsed correctly.

Availability of Data

Not all of the data is available right after a new association has been established. The time span until all data is collected depends on the internal update frequency of the data. Typical times are listed in the table below.

Object Type	Max. Time (Typical)
Numerics (real-time)	< 2 s
Numerics (12 second averaged)	< 18 s
Numerics (1 minute averaged)	< 70 s
Numerics (5 minute averaged)	< 310 s
Alarms	< 2 s
Patient Demographics	< 10 s
Medical Device System Object	< 1 s

During the startup phase, Poll Data Request messages on the object will result in Poll Data Response messages, which

- do not contain all the objects which are present in the IntelliVue monitor.
- do not contain all the available attributes of an object.

Numeric data is only available if a Measurement Server is connected to the IntelliVue monitor and if the system is not in stand-by mode. If a Measurement Server is connected to a running system, it may take several seconds until the data from the Measurement Server is available.

Parsing the Poll Result

The Poll Data Result message contains a checksum in the transport layer message. The Computer Client should verify that this checksum is correct. In the case of a corrupted checksum, the Computer Client must discard the message.

The Computer Client should check the poll number in the Poll Data Result message if it needs to detect lost messages. The Computer Client should check the *rel_time_stamp* which indicates the system time when the data was internally generated.

If the Computer Client needs to acquire a specific Numeric label (e.g., ABP), the preferred method is to use the *PhysioId* which is part of the Numeric Observed Value attribute (see “Numeric Objects” on page 75). The *physio_id* (physiological identifier) field contains a nomenclature code from the SCADA partition that identifies the represented value (typically a physiological measurement). It can be mapped to a label. However, for some numerics, the *physio_id* does not uniquely identify the measurement. E.g. all difference temperatures have the same *physio_id*, the numerics in the two channels of an EEG have the same *physio_ids*, the VueLink module may have numerics where the *physio_id* is not specified. However, if the label is derived by enumeration (e.g. the temperatures T1 and T2), the labels map to the same *PhysioId*. This ambiguity can be resolved if the user assigns other labels to the Numerics.

A Computer Client should not send Poll Requests for all attribute groups (*polled_attr_grp* = 0) when querying data with a high update frequency. Polling all attribute groups with a high frequency might lead to high system load and increased response latency. Future releases of the Data Export Protocol may support more attributes for each object.

If the IntelliVue monitor sends no response, check for the following causes:

- There is no association. Either the association was not established correctly or the IntelliVue monitor sent an Abort message (e.g., time-out) in the meantime.
- The Computer Client sent too many messages and messages were lost.
- The length of the transport layer message is corrupt.
- Length fields in the message are corrupt.

If the IntelliVue monitor sends a Remote Operation Error, this might have one of the following reasons:

- Wrong length field in the message.
- Wrong message type (*ro_type*, *command_type*, *action_type*).
- Wrong *managed_object* for the action (for Poll Requests, this must be the MDS object announced in the MDS Create Event).
- Wrong *polled_obj_type* (refer to “SINGLE POLL DATA REQUEST” on page 55 and “EXTENDED POLL DATA REQUEST” on page 59).
- Computer Client sent an Extended Poll Data Request, but the necessary optional package was not negotiated.
- Computer Client sent an Extended Poll Data Request with the wrong polled attribute group.
- Computer Client requested periodic Poll Data Result messages for too many objects. The Computer Client should at most send one request for Numerics (Metric Observed Value Attribute Group) and one for the AlertMonitor (Alert Monitor Attribute Group).

If the IntelliVue monitor sends a Poll Result message which does not contain all object/attributes check for the following problems:

- The Computer Client sent a Single Poll Data Request with the wrong polled attribute group. The Poll Result shows the objects with empty attribute lists (there are no attributes from the requested group).
- The association has been established and not all of the objects have been created. Wait until the objects are created.

Parsing AttributeLists

When parsing an *AttributeList*, the Computer Client should adhere to the following guidelines:

- Verify that the length fields in the *AttributeList* are consistent with other length fields in the message.
- Check both the count and length field of the *AttributeList* to detect the end of the list.
- Do not rely on the sequence of attributes in an *AttributeList*.
- Skip unknown attributes.
- Verify that the length field of each *AVAType* is consistent with its value.

If the Computer Client fails parsing the message, it is useful to compare the raw message (captured with a network monitor) with the Computer Client’s interpretation of the data. Common problems are:

- The Computer Client uses a different byte order. Wrong interpretation of length and count fields in particular can lead to problems.

- The Computer Client uses a different alignment for structures. The offset for members of a structure will be wrong, because the compiler for the Computer Client inserted bytes for alignment.
- Length fields denote the length of data appended, excluding the size of the length field.

Interpreting Data from Numerics

- Do not rely on the sequence of values within a Compound Numeric Observed Value attribute. The physiological identifiers must be interpreted.
- A triple valued pressure parameter can change to single valued (mean only), whenever the diastolic and systolic values are close together. This commonly happens when a pressure is being zeroed or when a transducer is left exposed to air. The parameter is still sent as a Compound Numeric Observed Value, even if only one value is available.
- The text in the label strings is localized. If you have a monitor with chinese localization, the strings will contain chinese UNICODE characters.

Interpreting Data from the Alert Monitor

- If the Computer Client wants to display Alarm messages, it should check the strings for UNICODE characters from the private use area (see “Definitions Shared by Protocols” on page 35).
- The text in the alarm strings is localized. If you have a monitor with chinese localization, the strings will contain chinese UNICODE characters.

Interpreting Wave Data

- The IntelliVue patient monitor supports the following wave types, which are defined by sample period, sample and array size (Sample Array Specification), and update period (Metric Specification) in the static context.

Wave Type	Sample Period	Sample Size	Array Size	Update Period	Bandwidth Requirement ¹
500 samples/s (ECG)	2 ms	16 bits	128 samples	256 ms	1064 bytes/s
250 samples/s (Compound ECG)	4 ms	16 bits	3*64 samples	256 ms	1640 bytes/s
125 samples/s	8 ms	16 bits	32 samples	256 ms	296 bytes/s
62.5 samples/s	16 ms	16 bits	16 samples	256 ms	168 bytes/s

1. Observed values, not including context data.

- The Computer Client can poll the dynamic context to determine the available waves. Because of the high amount of data, the client should specify the required wave objects before requesting wave observed values in a periodic data poll.
- Up to three ECG waves (500 samples/s) can be polled simultaneously by selecting the appropriate lead labels in the Wave object priority list. The object handle is the same for all ECG waves. Waves can be identified by their physiological identifier.
- It is possible to select up to three individual ECG waves with 500 sps each or the single ECGcompound wave (containing three channels, 250 sps each). Additionally up to eight 125 sps or 62.5 sps waves may be chosen. Bandwidth restrictions need to be considered (see table above for bandwidth usage of the individual wave types)

- In non-EASI mode, three ECG waves (250 samples/s, including the primary and secondary lead) can be polled by selecting the NLS_NOM_ECG_ELEC_POTL label in the Wave object priority list. The monitor sends poll results with a compound wave, containing three waves with common context. Waves can be identified by their physiological identifier.
- Up to eight non-ECG waves (125 or 62.5 samples/s) can be polled simultaneously by selecting the appropriate labels in the Wave object priority list.
- The Computer Client needs to keep track of the poll results time stamps to detect missing wave samples.
- Entries in the Wave object priority list are ignored if the label does not exist or the object is not available, or more than three ECG and/or more than eight non-ECG waves are specified.

The wave context can be polled separately or multiplexed with the wave observed values. If the *polled_attr_grp* is 0 in a periodic data poll request, the monitor reports one object's static and dynamic context per 1024 ms. Context attributes are included in the observation poll.

Troubleshooting

This chapter will help you identify and locate faults that may occur when using the Protocol. The procedure to locate faults uses a troubleshooting matrix.

When the fault has been identified, check the Possible Causes and corresponding Corrective Actions. Perform the corrective actions. Re-check the fault after each corrective action is performed until the fault has been cleared. It is assumed that you have a functioning Computer Client.

Fault	Possible Causes	Corrective Actions
Computer Client doesn't receive LAN messages	Cable connection is broken or wrong cable used.	<ul style="list-style-type: none"> Verify that the IntelliVue monitor is correctly connected to the network. Verify that the Computer Client is correctly connected to the network. Try to use an ICMP echo (ping) to check the monitor and Computer Client connections.
	IntelliVue monitor failure	Re-boot the IntelliVue monitor and try to make a new connection. Refer to the Troubleshooting section in the <i>Service Guide</i> of your device.
IntelliVue monitor shows an Unsupported LAN INOP	BootP server does not send a valid IP address.	Check the configuration of the BootP server. Check that the BootP server is correctly connected to the network.
	Cable connection is broken or wrong cable used.	Check the connection between the IntelliVue monitor and the BootP Server.
IntelliVue monitor shows a No Central Monitoring INOP	Central Monitoring Mandatory is configured to On in the monitor	Data Export must not be used with a central station. Configure Central Monitoring to Optional .
	Central Monitoring Mandatory is configured to On in the monitor and the connection to the central station is interrupted	Data Export must not be used with a central station. Reboot the IntelliVue monitor and make sure it is not connected to a central station.

Fault	Possible Causes	Corrective Actions
Computer Client doesn't receive messages with the AutoSpeed protocol	Cable connection is broken or wrong cable used.	Check the connection between the IntelliVue monitor and the Computer Client.
	Wrong configuration of MIB/RS232 Interface	Check if the MIB/RS232 interface is configured for the desired protocol
	IntelliVue monitor failure	Re-boot the IntelliVue monitor and try to make a new connection. Disconnect the MIB/RS232 cable for more than 60s, this will most likely reset the IrDA stack of the client system too. Refer to the Troubleshooting section in the <i>Service Guide</i> of your device.
Computer Client does not establish an association.	Another Computer Client Application is already associated with the IntelliVue monitor.	Make sure no other Computer Client Application is trying to connect to the IntelliVue monitor. Reboot the IntelliVue monitor or wait until the association is timed out.
Computer Client does not report data.	Measurement Server is disconnected.	Connect the measurement server to the IntelliVue monitor,
	Parameter is switched off.	If the Computer Client requires a specific measurement, the parameter must be switched on in the IntelliVue monitor.
	Wave label is not included in the Wave object priority list.	Specify the wave objects to be polled in the Set Priority List Request
Wave Samples are missing in a periodic data poll	Too many Wave objects polled.	Reduce the number of entries in the Wave object priority list.

NOTE The IntelliVue Data Export Interface cannot be accessed via the Local Area Network when the IntelliVue monitor is connected to the Philips LAN, e.g. to an Information Center (central station). Communication via the MIB/RS232 Interface is always possible.

Further Troubleshooting

Further troubleshooting can be done using the Philips Data Export Test Tool (DETT).

DETT is used to test the communication interface protocol, which transfers data from the Philips IntelliVue Patient Monitor via the Local Area Network (LAN) Interface or Serial Interface (MIB/RS232) to an external Computer.

Complete DETT functionality information is available in the DETT "Instructions for Use".

Download DETT 453564212161_DETT.zip file from InCenter at: <http://incenter.medical.philips.com>

A personal InCenter login account is required to access DETT.

Customers or users without a personal InCenter login requiring the DETT, please contact your local Philips Response Center for further support.

The DETT "Instructions for Use" (453564254321.pdf) can be downloaded from InCenter. See the link listed above and its description.

DETT "Instruction for Use" Information is also available within the DETT program via "Help".

Protocol Examples

Data Export Protocol Examples

CONNECT INDICATION EVENT

The Connect Indication message contains the *ConnectIndInfo* which is of variable length. The length fields in the message depend on the length of the *ConnectIndInfo*. This message is only available on the LAN interface.

```

Nomenclature      u_32          : 1.0
                   {0x00 0x00 0x01 0x00}
ROapdus           ro_type       : ROIV_APDU
                   length        : <xx>
                   {0x00 0x01 0xXX 0xXX}
ROIVapdu          invoke_id     : 0
                   command_type  : CMD_EVENT_REPORT
                   length        : <xx>
                   {0x00 0x00 0x00 0x00 0xXX 0xXX}

EventReportArg.   m_obj_class    : NOM_MOC_VMS_MDS_COMPOS_SINGLE_BED
ManagedObjectId  context_id     : 0
                  handle        : 0
RelativeTime      event_time    : 39424
OIDType           event_type    : NOM_NOTI_MDS_CONNECT_INDIC
u_16              length        : <xx>
                  {0x00 0x23 0x00 0x00 0x00 0x00 0x00 0x00 0x00
                    0x9A 0x00 0x0D 0x17 0xXX 0xXX}

ConnectIndInfo    [...]

```

MDS CREATE EVENT

The MDS Create Event message contains an *AttributeList* which is of variable length. The length fields in the message depend on the length of the *AttributeList*.

```

SPpdu             session_id     : 0xE100
                  p_context_id   : 2
                  {0xE1 0x00 0x00 0x02}
ROapdus           ro_type       : ROIV_APDU
                   length        : <xx>
                   {0x00 0x01 0xXX 0xXX}
ROIVapdu          invoke_id     : 1
                   command_type  : CMD_CONFIRMED_EVENT_REPORT
                   length        : <xx>
                   {0x00 0x01 0x00 0x01 0xXX 0xXX}

EventReportArg.   m_obj_class    : NOM_MOC_VMS_MDS
ManagedObjectId  context_id     : 0
                  handle        : 0
RelativeTime      event_time    : 126976
OIDType           event_type    : NOM_NOTI_MDS_CREAT
u_16              length        : <xx>
                  {0x00 0x21 0x00 0x00 0x00 0x00 0x00 0x00 0x01
                    0xf0 0x00 0x0d0x06 0xXX 0xXX}

MDSCreateInfo     ManagedObjectId m_obj_class    : NOM_MOC_VMS_MDS

```

```

                                context_id      : 0
                                handle         : 0
                                {0x00 0x21 0x00 0x00 0x00 0x00}
AttributeList [ ... ]

```

MDS CREATE EVENT RESULT

```

SPpdu      session_id      : 0xE100
           p_context_id    : 2
           {0xE1 0x00 0x00 0x02}
ROapdus    ro_type        : RORS_APDU
           length          : 20
           {0x00 0x02 0x00 0x14}
RORSapdu   invoke_id      : 1
           command_type    : CMD_CONFIRMED_EVENT_REPORT
           length          : 14
           {0x00 0x01 0x00 0x01 0x00 0x0e}

EventReportRes.
ManagedObjectId m_obj_class : NOM_MOC_VMS_MDS
                context_id   : 0
                handle       : 0
RelativeTime    event_time  : 4736768
OIDType        event_type  : NOM_NOTI_MDS_CREAT
u_16           length      : 0
                {0x00 0x21 0x00 0x00 0x00 0x00 0x00 0x00 0x48
                  0x47 0x00 0x0d 0x06 0x00 0x00}

```

SINGLE POLL DATA REQUEST

```

SPpdu      session_id      : 0xE100
           p_context_id    : 2
           {0xE1 0x00 0x00 0x02}
ROapdus    ro_type        : ROIV_APDU
           length          : 28
           {0x00 0x01 0x00 0x1c}
ROIVapdu   invoke_id      : 0
           command_type    : CMD_CONFIRMED_ACTION
           length          : 22
           {0x00 0x01 0x00 0x07 0x00 0x16}

ActionArgument
ManagedObjectId m_obj_class : NOM_MOC_VMS_MDS
                context_id   : 0
                handle       : 0
u_32           scope       : 0
OIDType        action_type : NOM_ACT_POLL_MDIB_DATA
u_16           length      : 8
                {0x00 0x21 0x00 0x00 0x00 0x00 0x00 0x00
                  0x00 0x00 0x0c 0x16 0x00 0x08}

PollMdibDataReq
u_16         poll_number   : 1
TYPE        partition    : NOM_PART_OBJ
           code           : NOM_MOC_VMO_METRIC_NU
OIDType     polled_attr_grp : all attribute groups
           {0x00 0x01 0x00 0x01 0x00 0x06 0x00 0x00}

```

SINGLE POLL DATA RESULT

The Single Poll Data Result message contains a *PollInfoList* which is of variable length. The length fields in the message depend on the length of the *PollInfoList*.

```

SPpdu      session_id      : 0xE100
           p_context_id    : 2
           {0xE1 0x00 0x00 0x02}
ROapdus    ro_type        : RORS_APDU
           length          : <xx>
           {0x00 0x02 0xXX 0xXX}
RORSapdu   invoke_id      : 0
           command_type    : CMD_CONFIRMED_ACTION
           length          : <xx>
           {0x00 0x00 0x00 0x07 0xXX 0xXX}

ActionResult
ManagedObjectId m_obj_class : NOM_MOC_VMS_MDS
                context_id   : 0
                handle       : 0
OIDType        action_type : NOM_ACT_POLL_MDIB_DATA
u_16           length      : <xx>

```

```

                                {0x00 0x21 0x00 0x00 0x00 0x00 0x0c 0x16
                                0xXX 0xXX}
PollMdbDataReply
u_16      poll_number          : 1
RelativeTime rel_time_stamp    : 4766464
AbsoluteTime abs_time_stamp    : 0xffffffff 0xffffffff
TYPE      partition           : NOM_PART_OBJ
          code                : NOM_MOC_VMO_METRIC_NU
OIDType   polled_attr_grp     : all attribute groups
          {0x00 0x01 0x00 0x48 0xbb 0x00 0xff 0xff
           0xff 0xff 0xff 0xff 0xff 0xff 0x00 0x01
           0x00 0x06 0x00 0x00}
PollInfoList [...]

```

SINGLE POLL DATA RESULT (LINKED)

It is assumed that the IntelliVue monitor needs two messages to encode all the data from a Poll Request.

The first message would have a linked result header:

```

SPpdu      session_id          : 0xE100
          p_context_id        : 2
          {0xE1 0x00 0x00 0x02}
ROapdus    ro_type            : ROLRS_APDU
          length              : <xx>
          {0x00 0x05 0xXX 0xXX}

ROLRSapdu
RorlsId     state              : ROLRS_FIRST
          count              : 1
u_16        invoke_id         : 0
CMDType     command_type      : CMD_CONFIRMED_ACTION
u_16        length            : <xx>
          {0x01 0x01 0x00 0x00 0x00 0x07 0xXX 0xXX}
ActionResult[...]

```

The second message would contain the rest of the data:

```

SPpdu      session_id          : 0xE100
          p_context_id        : 2
          {0xE1 0x00 0x00 0x02}
ROapdus    ro_type            : ROLRS_APDU
          length              : <xx>
          {0x00 0x05 0xXX 0xXX}

ROLRSapdu
RorlsId     state              : ROLRS_LAST
          count              : 2
u_16        invoke_id         : 0
CMDType     command_type      : CMD_CONFIRMED_ACTION
u_16        length            : <xx>
          {0x03 0x02 0x00 0x00 0x00 0x07 0xXX 0xXX}
ActionResult[...]

```

Finally, the monitor sends a Remote Operation Result message:

```

SPpdu      session_id          : 0xE100
          p_context_id        : 2
          {0xE1 0x00 0x00 0x02}
ROapdus    ro_type            : RORS_APDU
          length              : <xx>
          {0x00 0x02 0xXX 0xXX}
RORSapdu   invoke_id         : 0
          command_type        : CMD_CONFIRMED_ACTION
          length              : <xx>
          {0x00 0x00 0x00 0x07 0xXX 0xXX}
ActionResult [...]

```

Note that all messages contain a fully encoded *ActionResult* data structure. The last Remote Operation Result message, however, would contain a *PollInfoList* structure with the *count* and *length* field set to 0. A client system should not depend on the terminating Remote Operation Result to have an empty *PollInfoList*. The message should be parsed as any other message.

EXTENDED POLL DATA REQUEST

The next example shows a message which could be used to access averaged data. The message will only be accepted if the optional package for Poll Profile Extensions has been negotiated during the association phase.

```

SPpdu          session_id      : 0xE100
                p_context_id   : 2
                {0xE1 0x00 0x00 0x02}
ROapdus        ro_type        : ROIV_APDU
                length         : 32
                {0x00 0x01 0x00 0x20}
ROIvApdu        invoke_id      : 0
                command_type    : CMD_CONFIRMED_ACTION
                length          : 26
                {0x00 0x01 0x00 0x07 0x00 0x1a}

ActionArgument
ManagedObjectId m_obj_class    : NOM_MOC_VMS_MDS
                context_id     : 0
                handle         : 0
                scope          : 0
OIDType         action_type     : NOM_ACT_POLL_MDIB_DATA_EXT
u_16            length         : 12
                {0x00 0x21 0x00 0x00 0x00 0x00 0x00 0x00
                  0x00 0x00 0xf1 0x3b 0x00 0x0c}

PollMdibDataReqExt
u_16            poll_number     : 1
TYPE           partition       : NOM_PART_OBJ
                code           : NOM_MOC_VMO_METRIC_NU
OIDType         polled_attr_grp : all attribute groups
AttributeList
u_16            count          : 0
u_16            length         : 0
                {0x00 0x01 0x00 0x01 0x00 0x06 0x00 0x00
                  0x00 0x00 0x00 0x00}

```

EXTENDED POLL DATA RESULT

The Extended Poll Data Result message contains an additional *sequence_no*, which is used if the client requests periodic replies.

```

SPpdu          session_id      : 0xE100
                p_context_id   : 2
                {0xE1 0x00 0x00 0x02}
ROapdus        ro_type        : RORS_APDU
                length         : <xx>
                {0x00 0x02 0xXX 0xXX}
RORSapdu        invoke_id      : 0
                command_type    : CMD_CONFIRMED_ACTION
                length          : <xx>
                {0x00 0x00 0x00 0x07 0xXX 0xXX}

ActionResult
ManagedObjectId m_obj_class    : NOM_MOC_VMS_MDS
                context_id     : 0
                handle         : 0
OIDType         action_type     : NOM_ACT_POLL_MDIB_DATA_EXT
u_16            length         : <xx>
                {0x00 0x21 0x00 0x00 0x00 0x00 0xf1 0x3b
                  0xXX 0xXX}

PollMdibDataReplyExt
u_16            poll_number     : 1
u_16            sequence_no    : 0
RelativeTime    rel_time_stamp : 4766464
AbsoluteTime    abs_time_stamp : 0xffffffff 0xffffffff
TYPE           partition       : NOM_PART_OBJ
                code           : NOM_MOC_VMO_METRIC_NU
OIDType         polled_attr_grp : all attribute groups
                {0x00 0x01 0x00 0x00 0x00 0x48 0xbb 0x00
                  0xff 0xff 0xff 0xff 0xff 0xff 0xff 0xff
                  0x00 0x01 0x00 0x06 0x00 0x00}

PollInfoList    [...]

```

GET PRIORITY LIST REQUEST

```

SPpdu          session_id      : 0xE100

```

```

                                p_context_id      : 2
                                {0xE1 0x00 0x00 0x02}
ROapdus                        ro_type           : ROIV_APDU
                                length            : 22
                                {0x00 0x01 0x00 0x16}
ROIVapdu                       invoke_id         : 0
                                command_type      : CMD_GET
                                length            : 16
                                {0x00 0x00 0x00 0x03 0x00 0x10}
GetArgument
ManagedObjectId              m_obj_class        : NOM_MOC_VMS_MDS
                                context_id        : 0
                                handle            : 0
u_32                          scope             : 0
AttributeIdList              count             : 1
                                length            : 2
                                OIDType          : NOM_ATTR_POLL_RTSA_PRIO_LIST
                                {0x00 0x21 0x00 0x00 0x00 0x00 0x00 0x00
                                0x00 0x00 0x00 0x01 0x00 0x02 0xF2 0x3A}

```

GET PRIORITY LIST RESULT

```

SPpdu                          session_id        : 0xE100
                                p_context_id      : 2
                                {0xE1 0x00 0x00 0x02}
ROapdus                        ro_type           : RORS_APDU
                                length            : <xx>
                                {0x00 0x02 0xFF 0xFF}
RORSapdu                       invoke_id         : 0
                                command_type      : CMD_GET
                                length            : <xx>
                                {0x00 0x00 0x00 0x03 0xFF 0xFF}
GetResult
ManagedObjectId              m_obj_class        : NOM_MOC_VMS_MDS
                                context_id        : 0
                                handle            : 0
                                {0x00 0x21 0x00 0x00 0x00 0x00}
AttributeList                 count             : 1
                                length            : <xx>
AvaType                       attribute_id       : NOM_ATTR_POLL_RTSA_PRIO_LIST
                                length            : <xx>
                                {0x00 0x01 0xFF 0xFF 0xF2 0x3A 0xFF 0xFF}
TextIdList                    [...]

```

SET PRIORITY LIST REQUEST

```

SPpdu                          session_id        : 0xE100
                                p_context_id      : 2
                                {0xE1 0x00 0x00 0x02}
ROapdus                        ro_type           : ROIV_APDU
                                length            : <xx>
                                {0x00 0x01 0xFF 0xFF}
ROIVapdu                       invoke_id         : 0
                                command_type      : CMD_CONFIRMED_SET
                                length            : <xx>
                                {0x00 0x00 0x00 0x05 0xFF 0xFF}
SetArgument
ManagedObjectId              m_obj_class        : NOM_MOC_VMS_MDS
                                context_id        : 0
                                handle            : 0
u_32                          scope             : 0
                                {0x00 0x21 0x00 0x00 0x00 0x00 0x00 0x00
                                0x00 0x00}
ModificationList              count             : 1

```

```

                                length                : <xx>
AttributeModEntry              modifyOperator          : REPLACE
AvaType                        attribute_id            : NOM_ATTR_POLL_RTSA_PRIO_LIST
                                length                : <xx>
                                {0x00 0x01 0xXX 0xXX 0x00 0x00 0xF2 0x3A
                                0xXX 0xXX}
TextIdList                    [...]

```

SET PRIORITY LIST RESULT

```

SPpdu                         session_id              : 0xE100
                                p_context_id          : 2
                                {0xE1 0x00 0x00 0x02}
ROapdus                       ro_type                : RORS_APDU
                                length                : <xx>
                                {0x00 0x02 0xXX 0xXX}
RORSapdu                      invoke_id              : 0
                                command_type          : CMD_CONFIRMED_SET
                                length                : <xx>
                                {0x00 0x00 0x00 0x05 0xXX 0xXX}
SetResult
ManagedObjectId              m_obj_class             : NOM_MOC_VMS_MDS
                                context_id            : 0
                                handle                : 0
                                {0x00 0x21 0x00 0x00 0x00 0x00}
AttributeList                 count                  : 1
                                length                : <xx>
AvaType                       attribute_id            : NOM_ATTR_POLL_RTSA_PRIO_LIST
                                length                : <xx>
                                {0x00 0x01 0xXX 0xXX 0xF2 0x3A 0xXX 0xXX}
TextIdList                    [...]

```

AttributeList

This example shows an AttributeList which contains attributes from the Alert Monitor.

```

AttributeList                 count                  : 5
                                length                : 248
                                {0x00 0x05 0x00 0xf8}
AVAType                      attribute_id            : NOM_ATTR_ID_HANDLE
                                length                : 2
                                attribute_val         : 0x835d
                                {0x09 0x21 0x00 0x02 0x83 0x5d}
AVAType                      attribute_id            : NOM_ATTR_ID_TYPE
                                length                : 4
                                attribute_val         : 0x0001 0x0036
                                {0x09 0x2f 0x00 0x04 0x00 0x01 0x00 0x36}
AVAType                      attribute_id            : NOM_ATTR_DEV_AL_COND
                                length                : 10
                                attribute_val         : 0x1000 0x091a 0x0000 0x0002
                                                                0x0000
                                {0x09 0x16 0x00 0x0a 0x10 0x00 0x09 0x1a
                                0x00 0x00 0x00 0x02 0x00 0x00}
AVAType                      attribute_id            : NOM_ATTR_AL_MON_P_AL_LIST
                                length                : 4
                                attribute_val         : 0x0000 0x0000
                                {0x09 0x02 0x00 0x04 0x00 0x00 0x00 0x00}
AVAType                      attribute_id            : NOM_ATTR_AL_MON_T_AL_LIST
                                length                : 208
attribute_val                 : [...]
                                {0x09 0x04 0x00 0xd0 0x00 0x03 0x00 0xcc
                                0x4b 0xb8 0x01 0xba 0x00 0x02 0x10 0x00
                                0x00 0x02 0x00 0x00 0x83 0x3a 0x02 0x04
                                0x00 0x32 0x00 0x01 0x80 0x15 0x04 0x02
                                0x00 0x07 0x78 0x00 0x00 0x26 0x00 0x53
                                0x00 0x70 0x00 0x4f 0x20 0x82 0x00 0x20
                                0x00 0x4e 0x00 0x4f 0x00 0x4e 0x00 0x2d
                                0x00 0x50 0x00 0x55 0x00 0x4c 0x00 0x53
                                0x00 0x41 0x00 0x54 0x00 0x49 0x00 0x4c
                                0x00 0x45 0x00 0x00 0x50 0x00 0x01 0x12
                                0x00 0x02 0x10 0x00 0x00 0x09 0x00 0x00
                                0x02 0x91 0x02 0x04 0x00 0x32 0x00 0x01}

```

```
0x00 0x03 0x01 0x0c 0x00 0x00 0x78 0x00
0x00 0x26 0x00 0x52 0x00 0x65 0x00 0x73
0x00 0x70 0x00 0x20 0x00 0x20 0x00 0x20
0x00 0x4c 0x00 0x45 0x00 0x41 0x00 0x44
0x00 0x53 0x00 0x20 0x00 0x4f 0x00 0x46
0x00 0x46 0x00 0x20 0x00 0x20 0x00 0x00
0x4a 0x04 0x00 0xf2 0x00 0x02 0x10 0x00
0x00 0x02 0x00 0x00 0x82 0x63 0x02 0x04
0x00 0x32 0x00 0x01 0x00 0x03 0x00 0x46
0x00 0x02 0x78 0x00 0x00 0x26 0x00 0x4e
0x00 0x42 0x00 0x50 0x00 0x20 0x00 0x20
0x00 0x20 0x00 0x20 0x00 0x45 0x00 0x51
0x00 0x55 0x00 0x49 0x00 0x50 0x00 0x20
0x00 0x4d 0x00 0x41 0x00 0x4c 0x00 0x46
0x00 0x20 0x00 0x00 }
```

Association Control Protocol Examples

ASSOCIATION REQUEST

The following building blocks can be used to format an Association Request message:

AssocReqSessionHeader

```
0x0D <LI>
```

AssocReqSessionData

```
0x05 0x08 0x13 0x01 0x00 0x16 0x01 0x02
0x80 0x00 0x14 0x02 0x00 0x02
```

AssocReqPresentationHeader

```
0xC1 <LI> 0x31 0x80 0xA0 0x80 0x80 0x01
0x01 0x00 0x00 0xA2 0x80 0xA0 0x03 0x00
0x00 0x01 0xA4 0x80 0x30 0x80 0x02 0x01
0x01 0x06 0x04 0x52 0x01 0x00 0x01 0x30
0x80 0x06 0x02 0x51 0x01 0x00 0x00 0x00
0x00 0x30 0x80 0x02 0x01 0x02 0x06 0x0C
0x2A 0x86 0x48 0xCE 0x14 0x02 0x01 0x00
0x00 0x00 0x01 0x01 0x30 0x80 0x06 0x0C
0x2A 0x86 0x48 0xCE 0x14 0x02 0x01 0x00
0x00 0x00 0x02 0x01 0x00 0x00 0x00 0x00
0x00 0x00 0x61 0x80 0x30 0x80 0x02 0x01
0x01 0xA0 0x80 0x60 0x80 0xA1 0x80 0x06
0x0C 0x2A 0x86 0x48 0xCE 0x14 0x02 0x01
0x00 0x00 0x00 0x03 0x01 0x00 0x00 0xBE
0x80 0x28 0x80 0x06 0x0C 0x2A 0x86 0x48
0xCE 0x14 0x02 0x01 0x00 0x00 0x00 0x01
0x01 0x02 0x01 0x02 0x81
```

AssocReqUserData

The *AssocReqUserData* contains variable data, see see “Protocol Commands” on page 65.

AssocReqPresentationTrailer

```
0x00 0x00 0x00 0x00 0x00 0x00 0x00 0x00
0x00 0x00 0x00 0x00 0x00 0x00 0x00 0x00
```

ASSOCIATION RESPONSE

The following building blocks can be used to format an Association Response message:

AssocRespSessionHeader

```
0x0E <LI>
```

AssocRespSessionData

```
0x05 0x08 0x13 0x01 0x00 0x16 0x01 0x02
0x80 0x00 0x14 0x02 0x00 0x02
```

AssocRespPresentationHeader

```
0xC1 <LI> 0x31 0x80 0xA0 0x80 0x80 0x01
0x01 0x00 0x00 0xA2 0x80 0xA0 0x03 0x00
0x00 0x01 0xA5 0x80 0x30 0x80 0x80 0x01
0x00 0x81 0x02 0x51 0x01 0x00 0x00 0x30
0x80 0x80 0x01 0x00 0x81 0x0C 0x2A 0x86
0x48 0xCE 0x14 0x02 0x01 0x00 0x00 0x00
0x02 0x01 0x00 0x00 0x00 0x00 0x61 0x80
0x30 0x80 0x02 0x01 0x01 0xA0 0x80 0x61
0x80 0xA1 0x80 0x06 0x0C 0x2A 0x86 0x48
0xCE 0x14 0x02 0x01 0x00 0x00 0x00 0x03
0x01 0x00 0x00 0xA2 0x03 0x02 0x01 0x00
0xA3 0x05 0xA1 0x03 0x02 0x01 0x00 0xBE
0x80 0x28 0x80 0x02 0x01 0x02 0x81
```

AssocRespUserData

The *AssocRespUserData* contains variable data, see “Protocol Commands” on page 65.

AssocRespPresentationTrailer

```
0x00 0x00 0x00 0x00 0x00 0x00 0x00 0x00
0x00 0x00 0x00 0x00 0x00 0x00 0x00 0x00
```

REFUSE

The following building blocks can be used to format a Refuse message:

RefuseSessionHeader

0x0C 0x03

RefuseSessionData

0x32 0x01 0x00

RefusePresentationHeader

This block is empty in the Refuse message.

RefuseUserData

This block is empty in the Refuse message.

RefusePresentationTrailer

This block is empty in the Refuse message.

RELEASE REQUEST

The following building blocks can be used to format a Release Request message:

ReleaseReqSessionHeader

0x09 0x18

ReleaseReqSessionData

This block is empty in the Release Request message.

ReleaseReqPresentationHeader

0xC1 0x16 0x61 0x80 0x30 0x80 0x02 0x01
0x01 0xA0 0x80 0x62 0x80 0x80 0x01 0x00
0x00 0x00 0x00 0x00

ReleaseReqUserData

This block is empty in the Release Request message.

ReleaseReqPresentationTrailer

0x00 0x00 0x00 0x00

RELEASE RESPONSE

The following building blocks can be used to format a Release Response message:

ReleaseRespSessionHeader

0x0A 0x18

ReleaseRespSessionData

This block is empty in the Release Response message.

ReleaseRespPresentationHeader

0xC1 0x16 0x61 0x80 0x30 0x80 0x02 0x01
0x01 0xA0 0x80 0x63 0x80 0x80 0x01 0x00
0x00 0x00 0x00 0x00

ReleaseRespUserData

This block is empty in the Release Response message.

ReleaseRespPresentationTrailer

0x00 0x00 0x00 0x00

ASSOCIATION ABORT

The following building blocks can be used to format a Association Abort message:

AbortSessionHeader

0x19 0x2E

AbortSessionData

0x11 0x01 0x03

AbortPresentationHeader

0xC1 0x29 0xA0 0x80 0xA0 0x80 0x30 0x80
0x02 0x01 0x01 0x06 0x02 0x51 0x01 0x00
0x00 0x00 0x00 0x61 0x80 0x30 0x80 0x02
0x01 0x01 0xA0 0x80 0x64 0x80 0x80 0x01
0x01 0x00 0x00 0x00 0x00 0x00 0x00

AbortUserData

This block is empty in the Abort message.

AbortPresentationTrailer

0x00 0x00 0x00 0x00

User Data

The following section contains an example for the User Data which is contained in an Association Request message.

```

UserData
ASNLength      length      : 72
                {0x48}

MDSEUserInfoStd
ProtocolVersion protocol_version : MDDL_VERSION1
NomenclatureVers.nomenclature_version : NOMEN_VERSION
FunctionalUnits functional_units  : 0
SystemType      system_type     : SYST_CLIENT
StartupMode     startup_mode    : COLD_START
                {0x80 0x00 0x00 0x00 0x40 0x00 0x00 0x00
                 0x00 0x00 0x00 0x00 0x80 0x00 0x00 0x00
                 0x20 0x00 0x00 0x00}

Option List
AttributeList  count        : 0
                length      : 0
                {0x00 0x00 0x00 0x00}

Supported Profiles
AttributeList  count        : 1
                length      : 44
                {0x00 0x01 0x00 0x2c}

AVAType
OIDType       attribute_id   : NOM_POLL_PROFILE_SUPPORT
u_16          length        : 40
                {0x00 0x01 0x00 0x28}

PollProfileSupport (attribute_val)
PollProfileRev. poll_profile_revision : POLL_PROFILE_REV_0
RelativeTime    min_poll_period      : 800000
u_32            max_mtu_rx           : 1000
u_32            max_mtu_tx           : 1000
u_32            max_bw_tx            : 0xffff 0xffff
PollProfileOpt. options              : 0x6000 0x0000
                {0x80 0x00 0x00 0x00 0x00 0x00 0x09 0xc4
                 0x00 0x00 0x09 0xc4 0x00 0x00 0x03 0xe8
                 0xff 0xff 0xff 0xff 0x60 0x00 0x00 0x00}

Optional Packages
AttributeList  count        : 1
                length      : 12
                {0x00 0x01 0x00 0x0c}

AVAType
OIDType       attribute_id   : NOM_ATTR_POLL_PROFILE_EXT
u_16          length        : 8
                {0xf0 0x01 0x00 0x08}

PollProfileExt (attribute_val)
PollProfileExtOpt.options        : POLL_EXT_PERIOD_NU_AVG_60SEC
AttributeList  count            : 0
                length          : 0
                {0x20 0x00 0x00 0x00 0x00 0x00 0x00 0x00}

```

With this User Data, the length field of the Presentation Header must be set to 220 (0xDC) and the length field of the Session Header must be set to 236 (0xEC).